

Human Pluripotent Stem Cell-Derived Endoderm qPCR Array



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #07533
#07535

1 Plate
1 Plate

Type 3 (96 wells)
Type 5 (384 wells)

Product Description

Human Pluripotent Stem Cell (hPSC)-Derived Endoderm Quantitative Polymerase Chain Reaction (qPCR) Array is designed for characterization of hPSC-derived definitive endodermal progenitor cells and their differentiated progeny including pancreatic, hepatic, and intestinal cells. This qPCR array is designed for characterization of the gene expression profile of hPSC-derived endodermal cells following in vitro differentiation and can aid in development of protocols for directed differentiation towards endodermal lineages. Genes were selected based on their demonstrated expression in hPSCs, or in hPSC-derived ectodermal, mesodermal, and endodermal lineage cells (D'Amour et al.), including pancreatic (Kroon et al.), hepatic (Hay et al.), and intestinal progenitor cells (Grapin-Botton & Melton).

qPCR is a method for determining changes in steady-state mRNA levels of gene expression across multiple samples, generally normalized to the relative expression of internal control genes. Gene-specific primers are used in PCR to amplify target sequences within cDNA pools reverse-transcribed from mRNA. These PCR products contain hybridized sequence-specific probes that provide a fluorescent signal. Similar to TaqMan® technology, the fluorescent signal results from the 5' exonuclease activity of the Taq DNA polymerase on the probe, which is labeled with a reporter fluorophore (FAM) at the 5' end and a quencher (ZEN/IBFQ) at the 3' end. The rate of accumulation of the fluorescent signal is used to quantify the amount of cDNA present in the sample, and thereby the amount of mRNA present in the original cell lysate.

This qPCR array contains validated primers and probes for detection of 90 genes whose expression is correlated with hPSCs or hPSC-derived ectodermal, mesodermal, and endodermal cells, or pancreatic, hepatic, and intestinal progenitor cells. There are also 6 wells containing primers and probes for endogenous (housekeeping) control genes. TBP (TATA box-binding protein) qPCR Array Control Template is provided separately as a synthetic DNA positive control for use in a control well, containing primers and probes for TBP.

An annotated list of genes, as well as plate layouts and software for analysis of qPCR results, are available at www.stemcell.com/qPCRanalysis.

All kits include TBP qPCR Control Template (Component #07518). For more information on plate types and instrument compatibility, visit www.stemcell.com/DEqPCRinstruments.

Storage and Stability

Store plates at -20°C. Stable until expiry date (EXP) on boxtop label.

Store TBP qPCR Control Template at -20°C. Stable until expiry date (EXP) on label.

NOTE: Components may be shipped at room temperature (15 - 25°C) but should be stored at -20°C as indicated above.

Materials Required but Not Included

PRODUCT NAME	CATALOG #
Nuclease-Free Water (not DEPC-treated)	79001
Optical adhesive film	e.g. 38108
qPCR Master Mix Kit <ul style="list-style-type: none">• qPCR Master Mix (1 mL or 5 mL)• ROX Reference Dye (200 µL)	07516 (1 mL kit) OR 07517 (5 mL kit)
STEMscript™ cDNA Synthesis Kit with Oligo(dT) Primers OR STEMscript™ cDNA Synthesis Kit with Random Primers	79003 OR 79004
Total RNA Purification Kit	79040

Directions for Use

Please read the entire protocol before proceeding.

Use sterile technique when performing the following protocols.

Isolate RNA using the Total RNA Purification Kit. Quantify RNA by optical density at 260 nm, determine purity using $A_{260/280}$, then convert to cDNA using a STEMscript™ cDNA Synthesis Kit. Store cDNA at -20°C.

NOTE: We recommend an RNA concentration of 100 ng/μL in a 20 μL cDNA reaction. Use this amount of cDNA as is in the qPCR reaction below or dilute up to 1 in 5 (Table 2).

A. PREPARATION OF TBP qPCR CONTROL TEMPLATE AND cDNA COCKTAIL

1. Thaw qPCR Master Mix, cDNA, and ROX Reference Dye (if using) on ice.
2. If using ROX Reference Dye, add to qPCR Master Mix according to Table 1. For instruments not listed, refer to the manufacturer's instructions.

Table 1. Recommended Reference Dye Concentration Level for PCR Systems and Volume of ROX Reference Dye to Add to qPCR Master Mix

PCR SYSTEM	REFERENCE DYE CONCENTRATION LEVEL AND VOLUME OF ROX REFERENCE DYE (μL)/1 mL qPCR Master Mix		
	HIGH (40 μL)	LOW (4 μL)	NONE
Applied Biosystems • 7300 Real-Time PCR Systems	X		
Applied Biosystems • ViiA™ 7 and 7500 Real-Time PCR Systems • QuantStudio™ Systems		X	
Agilent Technologies • Mx3005P and Mx4000P		X	
Bio-Rad • CFX, iQ™, and DNA Engine Opticon® Real Time PCR Systems			X
Roche • LightCycler® Real-Time PCR System			X

3. Swirl bottle of qPCR Master Mix to mix thoroughly.
4. Prepare **TBP qPCR Control Template** as follows:
 - a. Centrifuge TBP qPCR Control Template at 3000 x *g* for 3 - 5 seconds to pellet material to the bottom of the vial.
 - b. Add 20 μL of nuclease-free water to the vial. Vortex the vial gently and thoroughly to resuspend the pellet.
 - c. Centrifuge at 3000 x *g* for 3 - 5 seconds to bring the liquid to the bottom of the vial.
5. Prepare **cDNA Cocktail** as follows:
 - a. Mix cDNA by gently pipetting up and down. Centrifuge at 3000 x *g* for 3 - 5 seconds to bring liquid to the bottom of the vial.
 - b. To a 15 mL conical tube (e.g. Catalog #38009), add components according to Table 2.

Table 2. Preparation of cDNA Cocktail

cDNA COCKTAIL COMPONENTS	VOLUME (μL)		
	1 well	96 wells*	384 wells*
cDNA	1	108	432
qPCR Master Mix	5	540	2160
Nuclease-free water	4	432	1728
Total Volume	10	1080	4320

*12.5% excess volume added to account for pipetting dead volume

- c. Cap the tube then gently vortex to mix thoroughly.
- d. Centrifuge at 3000 x *g* for 3 - 5 seconds to bring the liquid to the bottom of the tube.

B. PREPARATION OF qPCR PLATE

1. Carefully remove qPCR array plate from the box and plastic bag. Leave adhesive seal attached.
2. Centrifuge the plate at 1000 x *g* for 1 minute in a swinging bucket rotor fitted with plate holders.
3. Carefully remove and discard the adhesive seal on the plate.
4. Using a multichannel pipettor (e.g. Catalog #38064) and reagent reservoir (e.g. Catalog #38080), dispense reagents (from section A) into the plate wells as described below.
 - 5 μ L TBP qPCR Control Template + 5 μ L qPCR Master Mix in well H12 (see Figure 1)
 - 10 μ L cDNA Cocktail in all other wells

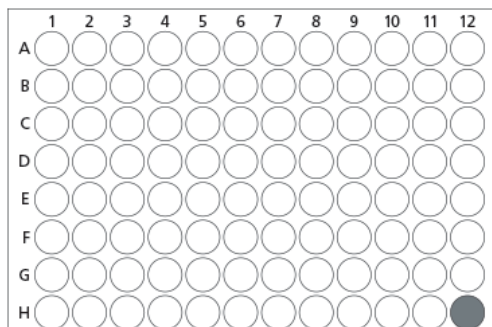


Figure 1. 96-Well Plate Diagram Indicating Well Containing TBP qPCR Control Template

- 5 μ L TBP qPCR Control Template + 5 μ L qPCR Master Mix in each of wells H12, H24, P12, and P24 (see Figure 2)
- 10 μ L cDNA Cocktail in all other wells

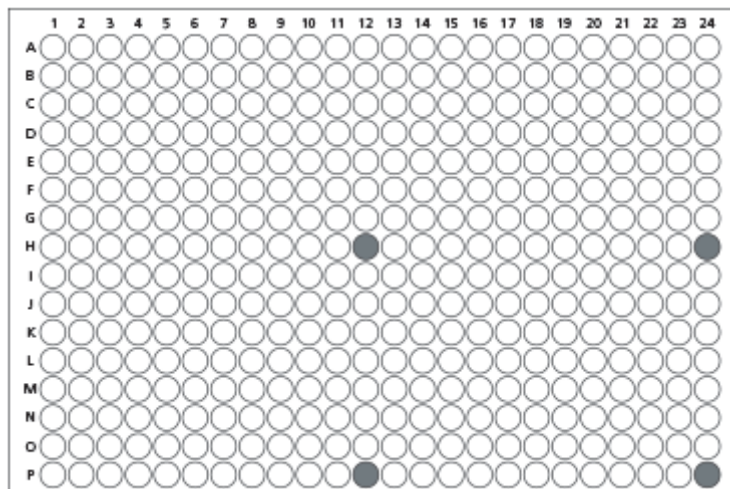


Figure 2. 384-Well Plate Diagram Indicating Wells Containing TBP qPCR Control Template

5. Carefully cover and seal the plate using optical adhesive film.
6. Centrifuge the sealed plate at 1000 x *g* for 1 minute at room temperature (15 - 25°C) to remove bubbles from the bottom of the wells.
NOTE: Bubbles in the bottom of the wells will interfere with results.
7. Place the plate on ice.

- C. qPCR
1. Program the thermocycler as indicated in Table 3.
- Table 3. Recommended qPCR Cycling Conditions**
- | STEP | TEMPERATURE | TIME | |
|---|-------------|----------------|------------------|
| | | Fast Cycling | Standard Cycling |
| Polymerase activation (1 cycle) | 95°C | 3 minutes | |
| Denaturation and annealing/extension (40 - 45 cycles) | 95°C | 5 seconds | 15 seconds |
| | 60°C* | 30 seconds* | 1 minute* |
| Hold | 4°C | Up to 24 hours | |
- *Annealing/extension temperature or time may need to be adjusted based on primer sequences.
2. If ROX Reference Dye is being used, calibrate thermocycler.
 3. Add plate and run PCR program.
 4. Save file including cycle threshold (Ct) values.
 5. Import the Ct data from the qPCR instrument to the analysis tool available at www.stemcell.com/qPCRanalysis. This analysis tool can rapidly and accurately quantitate relative gene expression, and the user can change analysis settings with ease.

Related Products

For related products, including specialized media, matrices, antibodies, cytokines, and small molecules, visit www.stemcell.com/DEworkflow, or contact us at techsupport@stemcell.com.

References

D'Amour KA et al. (2005) Efficient differentiation of human embryonic stem cells to definitive endoderm. *Nat Biotechnol* 23(12): 1534–41.

Grapin-Botton A & Melton DA. (2000) Endoderm development: from patterning to organogenesis. *Trends Genet* 16(3): 124–30.

Hay DC et al. (2008) Highly efficient differentiation of hESCs to functional hepatic endoderm requires ActivinA and Wnt3a signaling. *Proc Natl Acad Sci USA* 105(34): 12301–6.

Kroon E et al. (2008) Pancreatic endoderm derived from human embryonic stem cells generates glucose-responsive insulin-secreting cells in vivo. *Nat Biotechnol* 26(4): 443–52.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2024 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. TaqMan is a registered trademark of Roche Molecular Systems, Inc. QuantStudio and ViiA are trademarks of Thermo Fisher Scientific Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.