

Cyclic Pifithrin-Alpha


p53 inhibitor

Catalog #72062	5 mg
Catalog #72064	10 mg

Product Description

Cyclic Pifithrin-Alpha is a cell-permeable and reversible inhibitor of p53-mediated apoptosis and p53-dependent gene transcription. It is a more stable and less cytotoxic analog of the non-cyclic form of pifithrin-alpha. Cyclic Pifithrin-Alpha has also been reported to activate the aryl hydrocarbon receptor (Fernandez-Cruz et al.; Gary & Jensen; Komarov et al.). This product is supplied as a hydrobromide salt of the molecule.

Molecular Name:	Cyclic Pifithrin-Alpha (Hydrobromide)
Alternative Names:	Cyclic PFT- α , Cyclic Pifithrin- α , PFT- β , Pifithrin- β
CAS Number:	511296-88-1
Chemical Formula:	$C_{16}H_{16}N_2S \cdot HBr$
Molecular Weight:	349.3 g/mol
Purity:	$\geq 95\%$
Chemical Name:	5, 6, 7, 8-tetrahydro-2-(4-methylphenyl)-imidazo[2, 1-b]benzothiazole, monohydrobromide
Structure:	

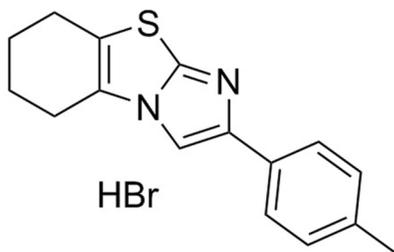

Properties

Product Format:	A crystalline solid
Stability and Storage:	Product stable at -20°C as supplied. As a precaution, STEMCELL recommends storing all small molecules away from direct light. Stable as supplied for 12 months from date of receipt.
Preparation:	<p>Solubility:</p> <ul style="list-style-type: none">• DMSO \leq 1.5 mM• Absolute ethanol \leq 1.5 mM <p>For example, to prepare a 1 mM stock solution in DMSO, resuspend 1 mg in 2.86 mL of fresh DMSO.</p> <p>Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.</p> <p>Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO or absolute ethanol concentration above 0.1% due to potential cell toxicity.</p>

Published Applications

MAINTENANCE AND SELF-RENEWAL

- Reduces UV-induced apoptosis of mouse embryonic stem cells (Qin et al.).
- Increases the numbers of mouse hematopoietic stem and progenitor cells in vivo and in vitro, also decreases the radiation-induced death of these cells (Leonova et al.).

REPROGRAMMING

- Increases efficiency of reprogramming mouse embryonic fibroblasts to induced pluripotent stem cells (Liao et al.).

References

Fernández-Cruz ML et al. (2011) Biological and chemical studies on aryl hydrocarbon receptor induction by the p53 inhibitor pifithrin- α and its condensation product pifithrin- β . *Life Sci* 88(17–18): 774–83.

Gary RK & Jensen DA. (2005) The p53 inhibitor pifithrin- α forms a sparingly soluble Derivative via intramolecular cyclization under physiological conditions. *Mol Pharm* 2(6): 462–74.

Komarov PG et al. (1999) A chemical inhibitor of p53 that protects mice from the side effects of cancer therapy. *Science* 285(5434): 1733–7.

Leonova KI et al. (2010) A small molecule inhibitor of p53 stimulates amplification of hematopoietic stem cells but does not promote tumor development in mice. *Cell Cycle* 9(7): 1434–43.

Liao J et al. (2013) Inhibition of PTEN tumor suppressor promotes the generation of induced pluripotent stem cells. *Mol Ther* 21(6): 1242–50.

Qin H et al. (2007) Regulation of apoptosis and differentiation by p53 in human embryonic stem cells. *J Biol Chem* 282(8): 5842–52.

Related Products

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2024 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.