EasySep™ Human CD19 Positive Selection Kit II

For processing 1 x 10⁹ cells

Catalog #17854

Positive Selection

Document #10000005308 | Version 01



Scientists Helping Scientists[™] | www.stemcell.com

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713
INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM
FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Description

Isolate highly purified CD19+ cells from fresh or previously frozen human peripheral blood mononuclear cells (PBMCs) or washed leukapheresis samples in as little as 18 minutes by immunomagnetic positive selection.

- · Fast and easy-to-use
- · Up to 99% purity
- · No columns required

This kit targets CD19+ cells for positive selection with antibodies recognizing the CD19 surface marker. Desired cells are labeled with antibodies and magnetic particles, and separated without columns using an EasySep™ magnet. Unwanted cells are simply poured off, while desired cells remain in the tube. Isolated cells are immediately available for downstream applications such as flow cytometry, culture, or DNA/RNA extraction.

Component Descriptions

COMPONENT NAME	COMPONENT #	QUANTITY	STORAGE	SHELF LIFE	FORMAT
EasySep™ Human CD19 Positive Selection Cocktail II	17854C	1 x 1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A combination of monoclonal antibodies in PBS. Includes an Fc receptor blocking antibody.
EasySep™ Dextran RapidSpheres™ 50100	50100	1 x 1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A suspension of magnetic particles in water.

PBS - phosphate-buffered saline

Components may be shipped at room temperature (15 - 25°C) but should be stored as indicated above.

Sample Preparation

For available fresh and frozen samples, see www.stemcell.com/primarycells.

PERIPHERAL BLOOD

Prepare a PBMC suspension from whole blood by centrifugation over a density gradient medium (e.g. Lymphoprep[™], Catalog #07801). For more rapid PBMC preparation, use the SepMate[™] RUO (Catalog #86450/86415) or SepMate[™] IVD* (Catalog #85450/85415) cell isolation tube.

If using previously frozen PBMCs, incubate the cells with DNase I Solution (Catalog #07900) at a concentration of 100 μ g/mL at room temperature (15 - 25°C) for at least 15 minutes prior to labeling and separation. Filter aggregated suspensions through a 37 μ m cell strainer (e.g. Catalog #27250) for optimal results.

After preparation, resuspend cells at 1 x 10⁸ cells/mL in recommended medium.

* SepMate™ IVD is only available in select regions where it is registered as an In Vitro Diagnostic (IVD) device for the isolation of mononuclear cells (MNCs) from whole blood or bone marrow by density gradient centrifugation. In all other regions SepMate™ is available for research use only (RUO).

LEUKAPHERESIS

Wash the peripheral blood leukapheresis sample by adding an equivalent volume of recommended medium or PBS containing 2% fetal bovine serum (FBS). Centrifuge at 500 x g for 10 minutes at room temperature (15 - 25°C). If red blood cell (RBC) lysis is necessary, lyse with Ammonium Chloride Solution (Catalog #07800). If platelet removal is necessary, centrifuge at 120 x g for 10 minutes with the brake off. Remove the supernatant and resuspend the cells at 1 x 10^8 cells/mL in recommended medium.

Recommended Medium

EasySep™ Buffer (Catalog #20144), RoboSep™ Buffer (Catalog #20104), or PBS containing 2% FBS and 1 mM EDTA. Medium should be free of Ca++ and Mg++.



Directions for Use - Manual EasySep™ Protocols

See page 1 for Sample Preparation and Recommended Medium. Refer to Tables 1 and 2 for detailed instructions regarding the EasySep™ procedure for each magnet.

Table 1. EasySep™ Human CD19 Positive Selection Kit II Protocol

		EASYSEP™ MAGNETS			
STEP	INSTRUCTIONS	EasySep™ (Catalog #18000)	"The Big Easy" (Catalog #18001)		
1	Prepare sample at the indicated cell concentration within the volume range.	1 x 10^8 cells/mL 0.1 - 2 mL	1 x 10^8 cells/mL 0.25 - 8 mL		
	Add sample to required tube.	5 mL (12 x 75 mm) polystyrene round-bottom tube (e.g. Catalog #38007)	14 mL (17 x 95 mm) polystyrene round-bottom tube (e.g. Catalog #38008)		
2	Add Selection Cocktail to sample. NOTE: Do not vortex cocktail.	100 μL/mL of sample	100 μL/mL of sample		
	Mix and incubate.	RT for 3 minutes	RT for 3 minutes		
3	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds	30 seconds		
4	Add RapidSpheres™ to sample.	100 μL/mL of sample	100 μL/mL of sample		
	Mix and incubate.	RT for 3 minutes	RT for 3 minutes		
5	Add recommended medium to top up the sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 2.5 mL	 Top up to 5 mL for samples < 1 mL Top up to 10 mL for samples ≥ 1 mL 		
	Place the tube (without lid) into the magnet and incubate.	RT for 3 minutes	RT for 3 minutes		
6	Pick up the magnet, and in one continuous motion invert the magnet and tube, pouring off the supernatant.* Remove the tube from the magnet; this tube contains the isolated cells.	Discard supernatant Discard supernatant			
7	Repeat steps as indicated.	Steps 5 and 6, three more times (total of 4 x 3-minute separations) Steps 5 and 6, two more times (total of 3 x 3-minute separations)			
8	Resuspend cells in desired medium. Be sure to collect cells from the sides of the tube.	Isolated cells are ready for use Isolated cells are ready for use			

RT - room temperature (15 - 25°C)
* Leave the magnet and tube inverted for 2 - 3 seconds, then return upright. Do not shake or blot off any drops that may remain hanging from the mouth of the tube.



Table 2. EasySep™ Human CD19 Positive Selection Kit II Protocol

		EASYSEP™ MAGNETS			
			EasyEights™ (Catalog #18103)	
STEP	INSTRUCTIONS	EasyPlate™ (Catalog #18102)	5 mL tube	14 mL tube	Easy 50 (Catalog #18002)
1	Prepare sample at the indicated cell concentration within the volume range.	1 x 10^8 cells/mL 0.05 - 0.2 mL	1 x 10^8 cells/mL 0.25 - 2 mL	1 x 10^8 cells/mL 1 - 8 mL	1 x 10^8 cells/mL > 5 - 40 mL
	Add sample to required tube (or plate if using the EasyPlate™ EasySep™ Magnet).	Round-bottom, non-tissue culture-treated 96-well plate (e.g. Catalog #38018)	5 mL (12 x 75 mm) polystyrene round-bottom tube (e.g. Catalog #38007)	14 mL (17 x 95 mm) polystyrene round-bottom tube (e.g. Catalog #38008)	50 mL (30 x 115 mm) conical tube (e.g. Catalog #38010)
2	Add Selection Cocktail to sample. NOTE: Do not vortex cocktail.	100 μL/mL of sample	100 μL/mL of sample	100 μL/mL of sample	100 μL/mL of sample
	Mix and incubate.	RT for 3 minutes	RT for 3 minutes	RT for 3 minutes	RT for 3 minutes
3	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds	30 seconds	30 seconds	30 seconds
4	Add RapidSpheres™ to sample.	100 μL/mL of sample	100 μL/mL of sample	100 μL/mL of sample	100 μL/mL of sample
4	Mix and incubate.	RT for 3 minutes	RT for 3 minutes	RT for 3 minutes	RT for 3 minutes
5	Add recommended medium to top up sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 0.25 mL	Top up to 2.5 mL	 Top up to 5 mL for samples ≤ 3 mL Top up to 10 mL for samples > 3 mL 	Top up to: • 10 mL for samples ≤ 5 mL • 20 mL for samples > 5 - 10 mL • 30 mL for samples > 10 - 15 mL • 40 mL for samples > 15 - 20 mL • 50 mL for samples > 20 - 40 mL
	Place the tube or plate (without lid) into the magnet and incubate.	RT for 5 minutes	RT for 10 minutes	RT for 10 minutes	RT for 10 minutes
6	Carefully pipette** (do not pour) off the supernatant. Remove the tube or plate, containing the isolated cells, from the magnet.	Discard supernatant	Discard supernatant	Discard supernatant	Discard supernatant
7	Add recommended medium to top up sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 0.25 mL	Top up to 2.5 mL	 Top up to 5 mL for samples ≤ 3 mL Top up to 10 mL for samples > 3 mL 	Top up to: • 10 mL for samples ≤ 5 mL • 20 mL for samples > 5 - 10 mL • 30 mL for samples > 10 - 15 mL • 40 mL for samples > 15 - 20 mL • 50 mL for samples > 20 - 40 mL
	Place the tube or plate (without lid) into the magnet and incubate.	RT for 5 minutes	RT for 5 minutes	RT for 5 minutes	RT for 5 minutes
8	Carefully pipette** (do not pour) off the supernatant. Remove the tube or plate, containing the isolated cells, from the magnet.	Discard supernatant	Discard supernatant	Discard supernatant	Discard supernatant
9	Repeat steps as indicated.	Steps 7 and 8 (total of 3 x 5-minute separations)	Steps 7 and 8 (total of 1 x 10 minute and 2 x 5-minute separations)	Steps 7 and 8 (total of 1 x 10 minute and 2 x 5-minute separations)	Steps 7 and 8 (total of 1 x 10 minute and 2 x 5-minute separations)
10	Resuspend cells in desired medium. Be sure to collect cells from the sides of the tube or plate.	Isolated cells are ready for use	Isolated cells are ready for use	Isolated cells are ready for use	Isolated cells are ready for use

RT - room temperature (15 - 25°C)

^{**} Collect the entire supernatant, all at once, into a single pipette (e.g. for EasyEights™ 5 mL tube use a 2 mL serological pipette [Catalog #38002]; for EasyEights™ 14 mL tube use a 10 mL serological pipette [Catalog #38004]).



Directions for Use – Fully Automated RoboSep™ Protocol

See page 1 for Sample Preparation and Recommended Medium. Refer to Table 3 for detailed instructions regarding the RoboSep™ procedure.

Table 3. RoboSep™ Human CD19 Positive Selection Kit II Protocol

STEP	INSTRUCTIONS	RoboSep™ (Catalog #21000)		
	Prepare sample at the indicated cell concentration within the volume range.	1 x 10^8 cells/mL 0.25 - 8 mL		
	Add sample to required tube.	14 mL (17 x 95 mm) polystyrene round-bottom tube (e.g. Catalog #38008)		
2	Select protocol.	Human CD19 Positive Selection II 17854		
3	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds		
4	Load the carousel.	Follow on-screen prompts		
4	Start the protocol.	Press the green "Run" button		
5	Unload the carousel when the run is complete. Remove the tube containing the isolated cells and resuspend in desired medium. Be sure to collect cells from the sides of the tube.	Isolated cells are ready for use		

Notes and Tips

ASSESSING PURITY


For purity assessment of CD19+ cells by flow cytometry, use one of the following fluorochrome-conjugated antibodies clones:


- · Anti-Human CD19 Antibody, Clone HIB19 (Catalog #60005; partially blocked), or
- · Anti-human CD19 antibody, clone 4G7 (partially blocked), or
- · Anti-human CD19 antibody, clone FMC63 (partially blocked)

One of the following methods can also be used:

- Use an alternative marker such as fluorochrome-conjugated Anti-Human CD20 Antibody, Clone 2H7 (Catalog #60008); this may underestimate CD19-positive
 purity by up to 15%.
- Use a fluorochrome-conjugated secondary antibody, such as Goat Anti-Mouse IgG (H+L) Antibody, Polyclonal (Catalog #60138).

Data

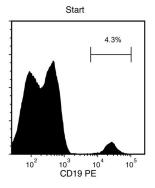

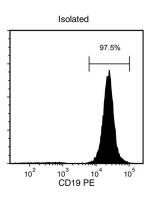

Starting with human PBMCs, the CD19+ cell content of the isolated fraction is typically 98 ± 1% (mean ± SD using the purple EasySep™ Magnet). In the above example, the purities of the start and final isolated fractions are 4.3% and 97.5%, respectively.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED. FOR ADDITIONAL INFORMATION ON QUALITY AT STEMCELL, REFER TO WWW.STEMCELL.COM/COMPLIANCE.

Copyright © 2021 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, EasyEights, EasyPlate, EasySep, RapidSpheres, RoboSep, and SepMate are trademarks of STEMCELL Technologies Canada Inc. Lymphoprep is a trademark of Alere Technologies. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.