

STEMdiff™ Endothelial Differentiation Kit

STEMdiff™ Endothelial Expansion Medium Kit



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #08005 1 Kit
Catalog #08007 1 Kit

Product Description

STEMdiff™ Endothelial Differentiation Kit (Catalog #08005) includes attachment substrate, animal component-free (ACF) endothelial induction medium, and endothelial expansion medium. It is optimized for the differentiation of human pluripotent stem cells (hPSCs) maintained in mTeSR™1 (Catalog #85850), mTeSR™ Plus (Catalog #05825), or TeSR™-E8™ (Catalog #05990) on Corning® Matrigel® to endothelial-like cells. The kit is designed to be used immediately after early mesoderm induction with STEMdiff™ Mesoderm Induction Medium (Catalog #05220), available for purchase separately.

STEMdiff™ Endothelial Expansion Medium Kit (Catalog #08007) includes the xeno-free endothelial expansion basal medium and supplement, and must be used in conjunction with Animal Component-Free Cell Attachment Substrate. STEMdiff™ Endothelial Expansion Medium is optimized for the expansion of hPSC-derived endothelial cells and for the maintenance of endothelial marker expression. For passaging, Animal Component-Free Cell Dissociation Kit (Catalog #05426) is required.

NOTE: Heparin Solution (Catalog #07980) is required for preparation of complete STEMdiff™ Endothelial Expansion Medium.

Product Information

PRODUCT NAME	CATALOG #	COMPONENTS
STEMdiff™ Endothelial Differentiation Kit	08005	<ul style="list-style-type: none"> STEMdiff™ Endothelial Induction Medium STEMdiff™ Endothelial Expansion Basal Medium STEMdiff™ Endothelial Expansion 5X Supplement Animal Component-Free Cell Attachment Substrate
STEMdiff™ Endothelial Expansion Medium Kit	08007	<ul style="list-style-type: none"> STEMdiff™ Endothelial Expansion Basal Medium STEMdiff™ Endothelial Expansion 5X Supplement

Component Storage and Stability

The following components are sold as part of complete kits (Catalog #08005 and 08007) and are not available for individual sale.

COMPONENT NAME	COMPONENT #	SIZE	STORAGE	SHELF LIFE
STEMdiff™ Endothelial Induction Medium	08006	100 mL	Store at -20°C.	Stable for 12 months from date of manufacture (MFG) on label.
STEMdiff™ Endothelial Expansion Basal Medium	08008	120 mL	Store at 2 - 8°C.	Stable until expiry date (EXP) on label.
STEMdiff™ Endothelial Expansion 5X Supplement*	08009	30 mL	Store at -20°C.	Stable until expiry date (EXP) on label.
Animal Component-Free Cell Attachment Substrate	07130	1 mL	Store at 2 - 8°C.	Stable until expiry date (EXP) on label.

*This component contains material derived from human plasma. Donors have been tested and found negative for HIV-1 and -2, hepatitis B, and hepatitis C prior to donation. However, this product should be considered potentially infectious and treated in accordance with universal handling precautions.

Materials Required But Not Included

PRODUCT NAME	CATALOG #
STEMdiff™ Mesoderm Induction Medium	05220
Heparin Solution	07980
Polyethersulfone (PES) filter unit (0.2 - 0.22 µm)	e.g. Fisher 09-741-04 (0.2 µm, 250 mL) OR Fisher SCGP00525 (0.22 µm, 50 mL)
D-PBS (Without Ca ⁺⁺ and Mg ⁺⁺)	37350
Corning® Matrigel® hESC-Qualified Matrix	Corning 354277
mTeSR™1	85850
Y-27632	72302
Tissue culture-treated 6-well plates	e.g. 38016
Conical tubes	e.g. 38009 (15 mL) or 38010 (50 mL)
ACCUASE™	07920
Trypan Blue	07050
Animal Component-Free Cell Dissociation Kit* • ACF Enzymatic Dissociation Solution • ACF Enzyme Inhibition Solution	05426

*For passaging and expansion of endothelial cells (section B of Directions for Use).

Preparation of Reagents and Materials

A. STEMdiff™ ENDOTHELIAL INDUCTION MEDIUM

Thaw STEMdiff™ Endothelial Induction Medium at 2 - 8°C overnight or at room temperature (15 - 25°C). If not used immediately, store at 2 - 8°C for up to 2 weeks.

B. STEMdiff™ ENDOTHELIAL EXPANSION MEDIUM

Use sterile technique to prepare STEMdiff™ Endothelial Expansion Medium (STEMdiff™ Endothelial Expansion Basal Medium + STEMdiff™ Endothelial Expansion 5X Supplement + Heparin Solution). The following example is for preparing 100 mL of complete medium. If preparing other volumes, adjust accordingly.

- Thaw 5X Supplement either at 2 - 8°C overnight or at 37°C until fully thawed. Mix thoroughly but do not vortex.
NOTE: Some precipitate may form. This will not affect product performance and will be removed when the complete medium is filtered (step 4).
NOTE: Once thawed, use immediately or aliquot and store at -20°C. Do not exceed the shelf life of the supplement. After thawing the aliquots, use immediately. Do not re-freeze.
- Warm STEMdiff™ Endothelial Expansion Basal Medium and Heparin Solution to room temperature (15 - 25°C).
- Add 20 mL of 5X Supplement to 79.4 mL of Basal Medium. Add 625 µL of Heparin Solution (final concentration 12.5 µg/mL). Mix thoroughly.
NOTE: Since Heparin Solution contains non-human animal-derived components, the complete medium will not be xeno-free.
- Filter the complete medium through a 0.2 - 0.22 µm PES filter unit.
NOTE: If not used immediately, store STEMdiff™ Endothelial Expansion Medium at 2 - 8°C for up to 2 weeks. If a precipitate forms, filter again as described above. This will not affect performance of the medium.

C. COATING CULTUREWARE WITH CORNING® MATRIGEL®

For complete instructions for coating cultureware with Corning® Matrigel®, refer to the Technical Manual: Maintenance of Human Pluripotent Stem Cells in mTeSR™1, available at www.stemcell.com or contact us to request a copy.

D. COATING CULTUREWARE WITH ANIMAL COMPONENT-FREE (ACF) CELL ATTACHMENT SUBSTRATE

For passaging and expansion of endothelial cells, coat cultureware with ACF Cell Attachment Substrate as described below.

NOTE: Use sterile technique; only use tissue culture-treated cultureware.

1. Dilute ACF Cell Attachment Substrate 1 in 100 in D-PBS (Without Ca⁺⁺ and Mg⁺⁺) (PBS).
For example, add 100 µL of ACF Cell Attachment Substrate to 9.9 mL of PBS.
2. Gently mix diluted ACF Cell Attachment Substrate. Do not vortex.
3. Immediately use diluted ACF Cell Attachment Substrate to coat cultureware. Refer to Table 1 for recommended coating volumes.

Table 1. Recommended Volumes for Coating Cultureware with Diluted ACF Cell Attachment Substrate


CULTUREWARE	VOLUME OF DILUTED ACF CELL ATTACHMENT SUBSTRATE
24-well plate	0.4 mL/well
6-well plate	1 mL/well
T-25 cm ² flask	2 - 3 mL/flask
T-75 cm ² flask	5 - 6 mL/flask

4. Gently rock cultureware back and forth to spread ACF Cell Attachment Substrate evenly across the surface. Seal with Parafilm®.
5. Incubate sealed cultureware at room temperature (15 - 25°C) for at least 2 hours before use. Do not let ACF Cell Attachment Substrate evaporate.

NOTE: Sealed coated cultureware can be stored at 2 - 8°C for up to 3 days after coating. Allow stored coated cultureware to warm to room temperature (15 - 25°C) for 30 minutes before proceeding to step 6.

6. Gently tilt cultureware onto one side and allow excess ACF Cell Attachment Substrate to collect at the edge. Remove excess solution using a serological pipette or by aspiration. Ensure that the coated surface is not scratched.
7. Wash cultureware twice with PBS (e.g. use 2 x 2 mL/well if using a 6-well plate).
8. Aspirate PBS. The coated cultureware is now ready for use.

Protocol Diagram

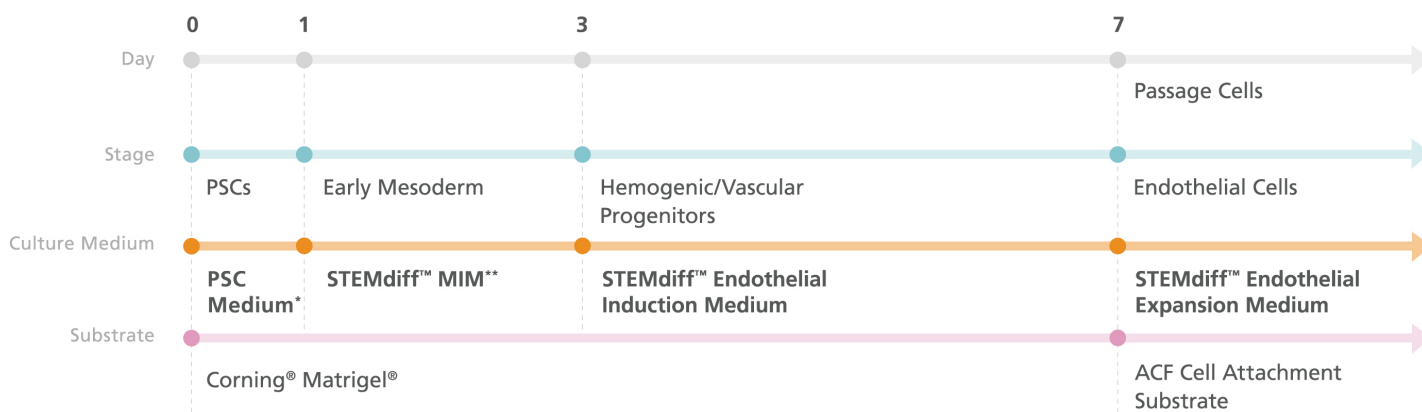

*mTeSR™1, mTeSR™ Plus, or TeSR™-E8™
 **MIM: Mesoderm Induction Medium

Directions for Use

For differentiation of hPSCs to endothelial cells, proceed to section A; for passaging and expansion of endothelial cells, proceed to section B.

A. DIFFERENTIATION OF hPSCs TO ENDOTHELIAL CELLS (MONOLAYER PROTOCOL)

The following protocol is for differentiation of hPSCs that have been maintained in mTeSR™1 on Corning® Matrigel®. mTeSR™ Plus or TeSR™-E8™ may also be used for hPSC maintenance; in the protocol below, use the same TeSR™ medium that was used for hPSC maintenance. When hPSC colonies are confluent and ready for passaging, proceed to Day 0: hPSC seeding.

Day 0: hPSC seeding

1. Prepare a single-cell suspension as described in section 7.2 of the Technical Manual: Maintenance of Human Pluripotent Stem Cells in mTeSR™1, available at www.stemcell.com or contact us to request a copy.
2. Seed hPSCs at a density of 5,000 to 10,000 cells/cm² in mTeSR™1 + 10 μM Y-27632 on a Matrigel®-coated 6-well plate.

NOTE: The seeding density may need to be adjusted depending on the cell line and hPSC maintenance medium used. Refer to Table 2 for recommended plating densities.

Table 2. Recommended Cell Plating Densities for Setting up Endothelial Induction

hPSC MAINTENANCE MEDIUM	PLATING DENSITY (cells/cm ²)	EXAMPLE OF NUMBER OF CELLS PER WELL OF A 6-WELL PLATE
mTeSR™1	5,000 to 10,000 cells/cm ²	5.0 x 10 ⁴ 7.5 x 10 ⁴ 10.0 x 10 ⁴
mTeSR™ Plus	3,500 to 7,500 cells/cm ²	3.5 x 10 ⁴ 5.0 x 10 ⁴ 7.5 x 10 ⁴
TeSR™-E8™	5,000 to 10,000 cells/cm ²	5.0 x 10 ⁴ 7.5 x 10 ⁴ 10.0 x 10 ⁴

Day 1: Mesoderm induction

For complete instructions on storage, handling, and thawing of STEMdiff™ Mesoderm Induction Medium, refer to the corresponding Product Information Sheet. The mesoderm induction stage generates early mesoderm cells from hPSCs.

The following instructions are for one well of a 6-well plate. For other cultureware, adjust volumes accordingly.

1. Warm the bottle of STEMdiff™ Mesoderm Induction Medium to room temperature (15 - 25°C).
2. Aspirate medium from the well. Slowly add 3 mL of STEMdiff™ Mesoderm Induction Medium down the side of the well.
3. Incubate at 37°C for 24 hours.

NOTE: Some degree of cell loss is expected at this stage.

Day 2

4. Repeat steps 1 - 3 above.

Day 3: Endothelial induction

5. Thaw STEMdiff™ Endothelial Induction Medium (see Preparation of Reagents and Materials).
6. Warm STEMdiff™ Endothelial Induction Medium to room temperature.
7. Aspirate medium from the well. Slowly add 4 mL of STEMdiff™ Endothelial Induction Medium down the side of the well.
8. Incubate at 37°C for 48 hours. Observe the cultures daily for changes in morphology.

Day 5

9. Repeat steps 6 - 8.

Day 7: Harvest & subculture of hPSC-derived endothelial cells

10. Prepare complete STEMdiff™ Endothelial Expansion Medium (see Preparation section).
11. Coat the desired culture vessel with ACF Cell Attachment Substrate (see Preparation section).
12. Harvest cells using ACCUTASE™ as follows:
 - a. Warm ACCUTASE™ to room temperature. Do not incubate at 37°C.
 - b. Gently wash the well twice with D-PBS (Without Ca⁺⁺ and Mg⁺⁺) (PBS).
 - c. Add 1 mL ACCUTASE™ per well.
 - d. Incubate at 37°C for 4 - 5 minutes.

- e. Gently tap the cultureware to detach cells. If less than 90% of the cells have detached, incubate at 37°C for an additional 1 - 2 minutes and tap again. Do not exceed 7 minutes of incubation.
 - f. Add 1 mL PBS to each well.
 - g. Gently pipette up and down to create a single-cell suspension and to lift any remaining cells. Add the cell suspension to a conical tube. Duplicate wells may be combined at this time.
 - h. Wash the well or flask with PBS and add the wash to the corresponding tube from step g.
 - i. Centrifuge the cell suspension at 300 x *g* for 5 minutes with the brake on.
 - j. Remove and discard the supernatant.
 - k. Resuspend each cell pellet in STEMdiff™ Endothelial Expansion Medium.
13. Perform a viable cell count using Trypan Blue and a hemocytometer.
 OPTIONAL: Perform flow cytometry on the harvested cells to assess surface marker expression.
14. Add cells to coated plates or flasks at the desired density in STEMdiff™ Endothelial Expansion Medium. Refer to Table 3 for recommended volumes. The optimal density is 10,000 cells/cm² for cells to reach confluency in 3 - 5 days.
 NOTE: Appropriate seeding density range is 5,000 - 15,000 cells/cm²; this may require optimization depending on the proliferative potential of the hPSC-derived endothelial cells.

Table 3. Recommended Volume of Medium for Various Cultureware

CULTUREWARE	VOLUME OF STEMdiff™ ENDOTHELIAL EXPANSION MEDIUM
6-well plate	2 mL/well
T-25 cm ² flask	4 mL/flask
T-75 cm ² flask	10 mL/flask

15. Incubate at 37°C and 5% CO₂. If the cells have not grown to confluence after 3 days, perform a full-medium change with STEMdiff™ Endothelial Expansion Medium and continue incubation. When cells have reached 90 - 100% confluency, proceed to section B for passaging/expansion. Refer to Figure 1 for a representative image.

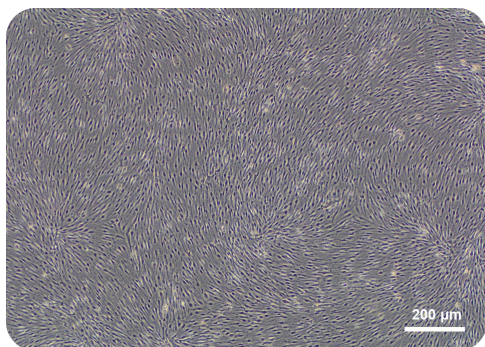


Figure 1. 90 - 100% Confluent Monolayer of H9-Derived Endothelial Cells at Passage 2 (4X Magnification)

B. PASSAGING AND EXPANSION OF ENDOTHELIAL CELLS

1. Coat the desired culture vessel(s) for expansion (e.g. T-25 cm² or T-75 cm² flask) with ACF Cell Attachment Substrate (see Preparation section).
2. Detach endothelial cells using Animal Component-Free Cell Dissociation Kit as follows:
 - a. Warm ACF Enzymatic Dissociation Solution and ACF Enzyme Inhibition Solution to room temperature (15 - 25°C). Do not incubate at 37°C.
 - b. Wash the well or flask twice with PBS.
 - c. Add ACF Enzymatic Dissociation Solution to the well or flask according to Table 4.

Table 4. Recommended Volume of ACF Enzymatic Dissociation Solution or ACF Enzyme Inhibition Solution for Various Cultureware

CULTUREWARE	VOLUME OF ACF ENZYMATIC DISSOCIATION SOLUTION OR ACF ENZYME INHIBITION SOLUTION
6-well plate	0.5 mL/well
T-25 cm ² flask	1 mL/flask
T-75 cm ² flask	2 mL/flask

- d. Incubate at 37°C for 4 - 5 minutes.
 - e. Gently tap the cultureware to detach cells. If < 90% of cells have detached, incubate at 37°C for an additional 1 - 2 minutes and tap again.
 - f. Add ACF Enzyme Inhibition Solution to the well or flask as indicated in Table 4 and collect cells in a polypropylene tube (e.g. 15 mL conical tube [Catalog #38009]).
 - g. Wash the well or flask with PBS and pipette up and down to lift any remaining cells. Add the wash to the corresponding tube from step f.
 - h. Centrifuge the cell suspension at 300 x *g* for 5 minutes with the **brake on**.
 - i. Remove and discard the supernatants. Resuspend each cell pellet in STEMdiff™ Endothelial Expansion Medium as follows:
 - For one well of a 6-well plate, resuspend in 0.3 - 0.5 mL medium
 - For a T-25 cm² flask, resuspend in 1 mL medium
 - For a T-75 cm² flask, resuspend in 2 mL medium
3. Perform a viable cell count using Trypan Blue and a hemocytometer.
 OPTIONAL: Perform flow cytometry on the harvested cells to assess surface marker expression.
4. Add cells to coated plates or flasks at the desired density in STEMdiff™ Endothelial Expansion Medium. The optimal density is 10,000 cells/cm² to reach confluency in 3 - 5 days.
 NOTE: Appropriate seeding density range is 5,000 - 15,000 cells/cm²; this may require optimization, depending on the proliferative potential of the hPSC-derived endothelial cells.
5. Incubate at 37°C and 5% CO₂. If the cells have not grown to confluency after 3 days, perform a full-medium change.

Assessment of hPSC-Derived Endothelial Cells

Assessment of hPSC-derived endothelial cells can be verified on **day 7** by flow cytometry after labeling with fluorochrome-conjugated antibodies (see list below for examples). On day 7, > 50% of cells express CD34, CD31, and CD144 and do not express lymphocyte marker CD45. The absence of undifferentiated cells can be confirmed by flow cytometry after labeling with fluorochrome-conjugated anti-OCT4 and anti-TRA-1-60. Results may vary depending on cell line used.

- Anti-Human CD34 Antibody, Clone 581 (Catalog #60013)
- Anti-Human CD31 antibody, clone WM59
- Anti-Human CD144 antibody, clone 55-7H1
- Anti-Human CD45 Antibody, Clone HI30 (Catalog #60018)
- Anti-Human CD309 antibody, clone 7D4-6
- Anti-Human CD105 Antibody, Clone 43A3 (Catalog #43A3)

Related Products

For related products, including specialized culture and storage media, supplements, antibodies, cytokines, and small molecules, visit www.stemcell.com or contact us at techsupport@stemcell.com.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED. FOR ADDITIONAL INFORMATION ON QUALITY AT STEMCELL, REFER TO WWW.STEMCELL.COM/COMPLIANCE.

Copyright © 2020 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, and STEMdiff are trademarks of STEMCELL Technologies Canada Inc. Corning and Matrigel are registered trademarks of Corning Inc. mTeSR, TeSR, and E8 are trademarks of WARF. ACCUTASE is a trademark of Innovative Cell Technologies, Inc. Parafilm is a registered trademark of Bemis Company, Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.