

# STEMdiff™ Trilineage Differentiation Kit



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

**Functional assay kit to assess pluripotency by directed differentiation of human ES and iPS cells to all three germ layers**

Catalog #05230

1 Kit

## Product Description

The STEMdiff™ Trilineage Differentiation Kit provides a simple culture assay to functionally validate the ability of new or established human embryonic stem (ES) and induced pluripotent stem (iPS) cell lines to differentiate to the three germ layers: ectoderm, mesoderm, and endoderm. This kit includes specialized, complete media and monolayer-based protocols to perform parallel in vitro directed differentiation experiments for each germ layer, clearly and reproducibly establishing trilineage differentiation potential within one week.

STEMdiff™ Trilineage Differentiation Kit is intended to be an endpoint assay and is not optimized for the generation of cells for downstream differentiation or other applications. STEMdiff™ Trilineage Differentiation Kit has been optimized to assess cells maintained in mTeSR™1.

## Product Information

The following components are sold as part of the STEMdiff™ Trilineage Differentiation Kit and are not available for individual sale.

COMPONENT NAME	COMPONENT #	SIZE	STORAGE	SHELF LIFE
STEMdiff™ Trilineage Ectoderm Medium	05231	175 mL	Store at -20°C.	Stable until expiry date (EXP) on label.
STEMdiff™ Trilineage Mesoderm Medium	05232	100 mL	Store at -20°C.	Stable until expiry date (EXP) on label.
STEMdiff™ Trilineage Endoderm Medium	05233	100 mL	Store at -20°C.	Stable until expiry date (EXP) on label.

## Materials Required But Not Included

PRODUCT NAME	CATALOG #
mTeSR™1	05850
DMEM/F-12 with 15 mM HEPES	36254
D-PBS (Without Ca++ and Mg++)	37350
Gentle Cell Dissociation Reagent	07174
Corning® Matrigel® hESC-Qualified Matrix	Corning 354277
Y-27632	72302
Trypan Blue	07050
ACCUTASE™	07920


## Preparation of Media

### STEMdiff™ Trilineage Ectoderm, Mesoderm, and Endoderm Media

Thaw at room temperature (15 - 25°C) or overnight at 2 - 8°C. Mix thoroughly.

NOTE: Do not filter medium. Once thawed, use immediately or store at 2 - 8°C for up to 2 weeks for Ectoderm and Mesoderm Medium or 1 week for Endoderm Medium. Alternatively, aliquot into polypropylene or PET-E tubes or bottles and store at -20°C. Do not exceed the expiry date as indicated on the label. After thawing the aliquots, use immediately or store at 2 - 8°C for up to 2 weeks for Ectoderm and Mesoderm Medium or 1 week for Endoderm Medium. Do not re-aliquot into additional tubes or bottles or re-freeze.

## Protocol Diagram

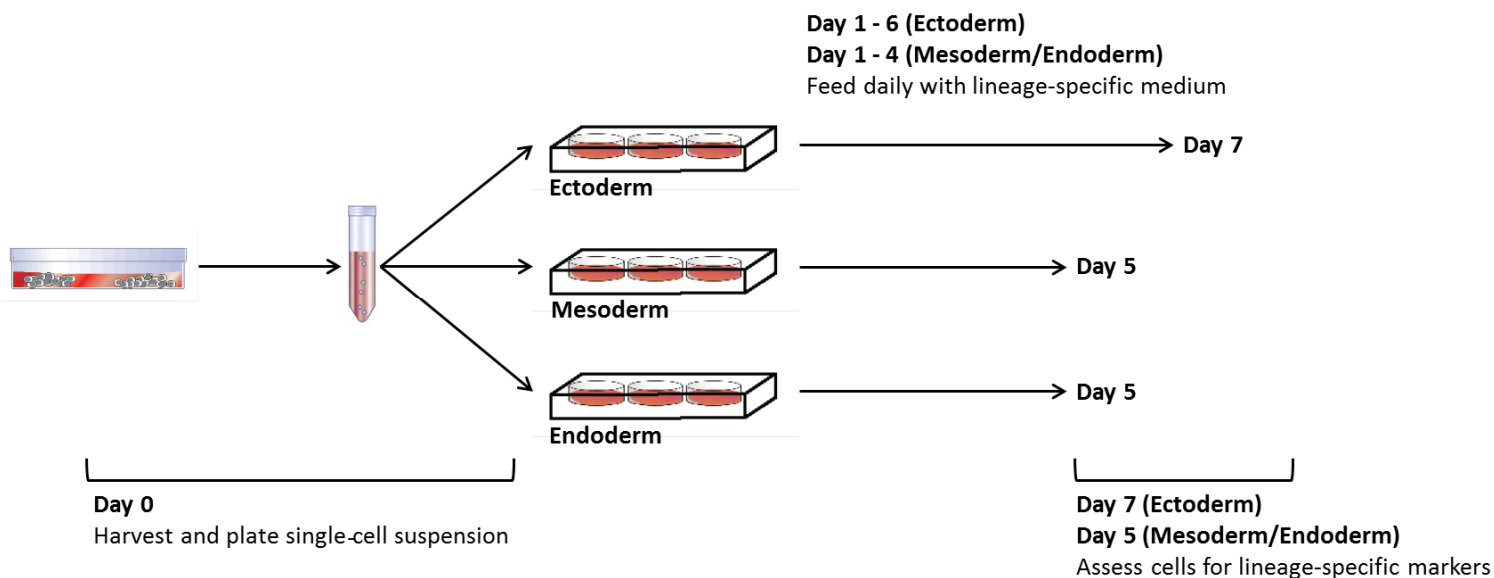

## Directions for Use

Please read the entire protocol before proceeding.

The following protocol is for generating cells of ectoderm, mesoderm, and endoderm lineages from human ES or iPS cells previously cultured in mTeSR™1. When ES/iPS cells are approximately 70% confluent, they are ready for passaging. The majority of colonies should be large, compact, and exhibit multi-layering at the center.

NOTE: For complete instructions on coating plates with Corning® Matrigel® and maintaining high-quality human ES and iPS cells for use in differentiation, refer to the Technical Manual: Maintenance of Human Pluripotent Stem Cells in mTeSR™1 (Document #29106) available on our website at [www.stemcell.com](http://www.stemcell.com) or contact us to request a copy.

### A. PLATING HUMAN ES/iPS CELLS FOR TRILINEAGE DIFFERENTIATION

The following are instructions for passaging human ES/iPS cells from wells of a 6-well plate to set up the Trilineage Differentiation assay in 24-well plates. If using alternative cultureware, adjust volumes accordingly.

Coat cultureware with Corning® Matrigel® hESC-Qualified Matrix and bring to room temperature (15 - 25°C) for at least 30 minutes prior to use.

- On **day 0**, warm (15 - 25°C) sufficient volumes of mTeSR™1, DMEM/F-12, STEMdiff™ Trilineage Ectoderm Medium, and Gentle Cell Dissociation Reagent for passaging.
- Prepare Single-Cell Plating Media as follows:  
**Ectoderm Single-Cell Plating Medium:** Add Y-27632 to STEMdiff™ Trilineage Ectoderm Medium to reach a final concentration of 10 µM.  
**Mesoderm/Endoderm Single-Cell Plating Medium:** Add Y-27632 to mTeSR™1 to reach a final concentration of 10 µM.
- Use a microscope to visually identify regions of differentiation in the wells to be passaged. Mark these using a felt tip or lens marker on the bottom of the plate. Remove regions of differentiation by scraping with a pipette tip or by aspiration.
- Wash each well to be passaged with 1 mL of D-PBS (Without Ca<sup>++</sup> and Mg<sup>++</sup>).
- Aspirate the wash medium and add 1 mL/well of Gentle Cell Dissociation Reagent.
- Incubate at 37°C for 8 - 10 minutes.
- In each well, dislodge cells by pipetting up and down 1 - 3 times using a pipette with a 1000 µL tip. Ensure all remaining cell aggregates are broken up into single cells.
- Immediately transfer cells to a tube containing 1 mL of DMEM/F-12 per well harvested. Wash each well once with 1 mL of DMEM/F-12 to collect any remaining cells and transfer to the tube. Centrifuge the tube at 300 x *g* for 5 minutes.
- Remove supernatant. Resuspend cells in 1 mL of Mesoderm/Endoderm Single-Cell Plating Medium (prepared in step 2).

NOTE: Cells from this suspension will be added into wells containing the appropriate lineage-specific plating medium in step 12.

10. Count viable cells using Trypan Blue and a hemocytometer.
11. Aspirate Matrigel® from pre-coated 24-well plates. Add 0.5 mL of the appropriate Single-Cell Plating Medium per well. Recommended volumes for plating in various types of cultureware are shown in Table 1.

**Table 1. Recommended Volumes of STEMdiff™ Trilineage Medium (Ectoderm, Mesoderm, or Endoderm) for Various Cultureware**

CULTUREWARE	VOLUME OF SINGLE-CELL PLATING MEDIUM	FEED VOLUME OF STEMdiff™ TRILINEAGE MEDIUM
24-well plate	0.5 mL/well	1 mL/well
12-well plate	1 mL/well	1.5 mL/well
6-well plate	2 mL/well	3 mL/well

12. Add the appropriate number of cells to the medium-containing wells. The required densities for each lineage are indicated in Table 2.  
NOTE: The volume of cells added should not exceed 20% of the total final volume in the well.

**Table 2. Plating Densities for Ectoderm, Mesoderm, and Endoderm Lineages**

LINEAGE	CELL DENSITY (cells/cm <sup>2</sup> )	TOTAL NUMBER OF CELLS PER WELL			SINGLE-CELL PLATING MEDIUM
		24-well plate	12-well plate	6-well plate	
Ectoderm	200,000	400,000	800,000	2,000,000	STEMdiff™ Trilineage Ectoderm Medium + 10 μM Y-27632
Mesoderm	50,000	100,000	200,000	500,000	mTeSR™1 + 10 μM Y-27632
Endoderm	200,000	400,000	800,000	2,000,000	mTeSR™1 + 10 μM Y-27632

13. Place the plate in a 37°C incubator. Move the plate in several quick, short, back-and-forth and side-to-side motions to evenly distribute the cells. Do not disturb the plate for 24 hours.
14. Continue to section B for differentiation.

## B. DIFFERENTIATING MONOLAYER CULTURES TO THREE GERM LINEAGES

1. On **day 1**, warm the three types of STEMdiff™ Trilineage media (see Preparation of Media) to room temperature (15 - 25°C).
2. Aspirate media from cell cultures.
3. Add the appropriate STEMdiff™ Trilineage medium to each well. Recommended feed volumes for various types of cultureware are shown in Table 1.
4. Incubate at 37°C for 24 hours.
5. Repeat steps 1 - 4 until **day 5** (mesoderm and endoderm lineages) or **day 7** (ectoderm lineage).
6. Proceed to section C for assessing differentiation.

## C. ASSESSING DIFFERENTIATION TO THREE GERM LINEAGES

Cells should be harvested and/or fixed for analysis of lineage-specific markers on **day 5** for mesoderm and endoderm lineages and **day 7** for the ectoderm lineage. The analysis method is determined by the user and may include flow cytometry, immunocytochemistry, or transcriptome analysis. If assessing for differentiation by flow cytometry or immunocytochemistry, label with combinations of fluorochrome-conjugated antibodies specific for early differentiation for each germ layer, for example:

- Ectoderm: PAX6 combined with Nestin (e.g. Anti-Human Nestin Antibody, Clone 10C2, Catalog #60091)
- Mesoderm: Brachyury (T) combined with either NCAM or CXCR4 (e.g. Anti-Human CD56 [NCAM] Antibody, Clone HCD56, Catalog #60021; or Anti-Human CD184 [CXCR4] Antibody, Clone 12G5, Catalog #60089)
- Endoderm: SOX17 combined with either CXCR4 or FOXA2

## Related Products

For related products, including specialized media, matrices, antibodies, cytokines, and small molecules, visit [www.stemcell.com/hPSCworkflow](http://www.stemcell.com/hPSCworkflow) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, and STEMdiff are trademarks of STEMCELL Technologies Canada Inc. Corning and Matrigel are registered trademarks of Corning Incorporated. mTeSR is a trademark of WARF. ACCUTASE is a trademark of Innovative Cell Technologies, Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.