

Human Pluripotent Stem Cell Trilineage Differentiation qPCR Array



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #07515

1 Kit

Product Description

The Human Pluripotent Stem Cell (hPSC) Trilineage Differentiation Quantitative Polymerase Chain Reaction (qPCR) Array is designed for characterization of hPSCs and their trilineage differentiation capacity. hPSCs, including embryonic stem (ES) and induced pluripotent stem (iPS) cells, are self-renewing and have the ability to differentiate into cells of the three embryonic germ layers: ectoderm, mesoderm, and endoderm. This qPCR array is designed for characterization of the gene expression profile of undifferentiated ES and iPS cells and their trilineage derivatives following in vitro directed or spontaneous differentiation, thereby validating the ability of a cell line to differentiate to the three germ layers. Genes were selected based on their demonstrated differential expression in ES and iPS cells compared with hPSC-derived ectodermal, mesodermal, and endodermal lineage cells (Adewumi et al.; Bock et al.).

qPCR is a method for determining changes in steady-state mRNA levels of gene expression across multiple samples, generally normalized to the relative expression of internal control genes. Gene-specific primers are used in PCR to amplify target sequences within cDNA pools reverse-transcribed from mRNA. These PCR products contain hybridized sequence-specific probes that provide a fluorescent signal. Similar to TaqMan® technology, the fluorescent signal results from the 5' exonuclease activity of the Taq DNA polymerase on the probe, which is labeled with a reporter fluorophore at the 5' end and a quencher fluorophore at the 3' end. The rate of accumulation of the fluorescent signal is used to quantify the amount of cDNA present in the sample, and thereby the amount of mRNA present in the original cell lysate.

This 384-well qPCR array consists of 4 x 96-well block replicates. Each 96-well quadrant contains validated primers and probes for detection of 90 genes whose expression is correlated with undifferentiated hPSCs or their derivatives undergoing the early stages of differentiation. There are also 24 wells (6 per quadrant) containing primers and probes for endogenous (housekeeping) control genes. TATA box-binding protein (TBP) qPCR Array Control Template is provided separately as a synthetic DNA positive control, for use in a control well in each quadrant that contains primers and probes for TBP.

An annotated list of genes, as well as plate layouts and software for analysis of qPCR results, are available at www.stemcell.com/qPCRanalysis.

Product Information

The following components are sold as part of a kit (Catalog #07515) and are not available for individual sale.

NAME	COMPONENT #	SIZE	STORAGE	SHELF LIFE
TBP qPCR Control Template	07518	10 ⁶ Copies	Store at -20°C.	Stable until expiry date (EXP) on label.
hPSC Trilineage Differentiation 384-Well qPCR Array Plate*	07502	1 Plate	Store at -20°C.	Stable until expiry date (EXP) on boxtop label.

*For instrument compatibility, visit www.stemcell.com/trilineageqPCRinstruments.

NOTE: Components may be shipped at room temperature (15 - 25°C) but should be stored at -20°C as indicated above.

Materials Required But Not Included

PRODUCT NAME	CATALOG #
qPCR Master Mix Kit <ul style="list-style-type: none">• qPCR Master Mix (1 mL or 5 mL)• ROX Reference Dye (200 µL)	07516 (1 mL kit) OR 07517 (5 mL kit)
Nuclease-free water (not DEPC-treated)	e.g. Thermo Fisher AM9932
Optical adhesive film	e.g. Thermo Fisher 4311971

Directions for Use

Please read the entire protocol before proceeding.

Use sterile techniques when performing the following protocols.

Isolate RNA using standard protocols. Quantify RNA by optical density at 260 nm, determine purity using $A_{260/280}$, then convert to cDNA. Store cDNA at -20°C .

NOTE: Optimal concentration of cDNA for qPCR amplification is 20 - 100 ng/ μL .

A. PREPARATION OF TBP qPCR CONTROL COCKTAIL AND cDNA COCKTAIL

1. Thaw qPCR Master Mix, cDNA, and ROX Reference Dye (if using) on ice.
2. If using ROX Reference Dye, add to qPCR Master Mix according to Table 1. For instruments not listed, refer to the manufacturer's instructions.

Table 1. Volume of ROX Reference Dye to add to qPCR Master Mix

PCR SYSTEM	VOLUME OF ROX REFERENCE DYE (μL)	
	1 mL qPCR Master Mix	5 mL qPCR Master Mix
7900HT Fast (Applied Biosystems)	40	200
ViiA™ 7 (Applied Biosystems)	4	20
QuantStudio™ Flex (Applied Biosystems)	4	20

3. Swirl bottle of qPCR Master Mix to mix thoroughly.
4. Prepare **TBP qPCR Control Cocktail** as follows:
 - a. Centrifuge TBP qPCR Control Template at $3000 \times g$ for 3 - 5 seconds to pellet material to the bottom of the vial.
 - b. Add 20 μL of nuclease-free water and 20 μL of qPCR Master Mix to the vial. Vortex the vial gently and thoroughly to resuspend the pellet.
 - c. Centrifuge at $3000 \times g$ for 3 - 5 seconds to bring the liquid to the bottom of the vial.
5. Prepare **cDNA Cocktail** as follows:
 - a. Mix cDNA by gently pipetting up and down. Centrifuge at $3000 \times g$ for 3 - 5 seconds to bring liquid to the bottom of the vial.
 - b. To a 15 mL conical tube (e.g. Catalog #38009), add components according to Table 2.

Table 2. Preparation of cDNA Cocktail

cDNA Cocktail Components	VOLUME (μL)	
	1 well	384 wells*
cDNA	1	432
qPCR Master Mix	5	2160
Nuclease-free water	4	1728
Total Volume	10	4320

*12.5% excess volume added to account for pipetting dead volume

- c. Cap the tube then gently vortex to mix thoroughly.
- d. Centrifuge at $3000 \times g$ for 3 - 5 seconds to bring the liquid to the bottom of the tube.

B. PREPARATION OF qPCR PLATE

1. Carefully remove qPCR array plate from the box and plastic bag. Leave adhesive seal attached.
2. Centrifuge the plate at $1000 \times g$ for 1 minute in a swinging bucket rotor fitted with plate holders. Ensure the plate is well balanced.
3. Carefully remove and discard the adhesive seal on the plate.

4. Dispense cocktails (from section A) into the plate wells as described below.
 - 10 μ L **TBP qPCR Control Cocktail** in each of wells H12, H24, P12, and P24 (see Figure 1)
 - 10 μ L **cDNA Cocktail** in all other wells

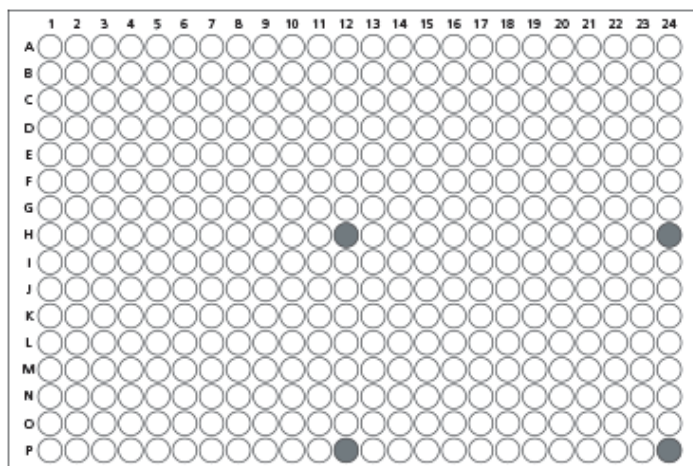


Figure 1. 384-Well Plate Diagram Indicating Wells Containing TBP qPCR Control Cocktail

5. Carefully cover and seal the plate using optical adhesive film.
 6. Centrifuge the sealed plate at 1000 x *g* for 1 minute at room temperature (15 - 25°C) to remove bubbles from the bottom of the wells.
NOTE: Bubbles in the bottom of the wells will interfere with results.
 7. Place the plate on ice.
- C. qPCR
1. Program the thermocycler according to manufacturer's instructions.
 2. If ROX Reference Dye is being used, calibrate thermocycler.
 3. Add plate and run PCR program.
 4. Save file including Ct (cycle threshold) values.
 5. Import the Ct data from the qPCR instrument to the analysis tool available at www.stemcell.com/qPCRanalysis. This analysis tool can rapidly and accurately quantitate relative gene expression, and the user can change analysis settings with ease.

Related Products

For related products, including specialized cell culture and storage media, matrices, antibodies, cytokines, and small molecules, visit www.stemcell.com/hPSCworkflow or contact us at techsupport@stemcell.com.

References

- Adewumi O et al. (2007) Characterization of human embryonic stem cell lines by the International Stem Cell Initiative. *Nat Biotechnol* 25(7): 803–16.
- Bock C et al. (2011) Reference maps of human ES and iPS cell variation enable high-throughput characterization of pluripotent cell lines. *Cell* 144(3): 439–52.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. TaqMan is a registered trademark of Roche Molecular Systems, Inc. QuantStudio and ViiA are trademarks of Thermo Fisher Scientific Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.