

# Small Molecules

PD173074

Tyrosine kinase inhibitor; Inhibits FGFR

Catalog #72164

10 mg



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

PD173074 is a selective and potent, ATP-competitive inhibitor of fibroblast growth factor receptor (FGFR). It acts on both FGFR3 and FGFR1 ( $IC_{50}$  = 5 and 21.5 nM, respectively), and also inhibits FGFR2, FGFR4, and vascular endothelial growth factor receptor 2 (VEGFR2). It is approximately 1000 times more potent than another common FGFR inhibitor SU5402 (Catalog #73912). PD173074 shows little to no activity against PDGFR, EGFR, MEK, or PKC (Koziczak et al.; Mohammadi et al.; Trudel et al.).

Molecular Name:	PD173074
Alternative Names:	Not applicable
CAS Number:	219580-11-7
Chemical Formula:	$C_{28}H_{41}N_7O_3$
Molecular Weight:	523.7 g/mol
Purity:	≥ 95%
Chemical Name:	N-[2-[[4-(diethylamino)butyl]amino]-6-(3,5-dimethoxyphenyl)pyrido[2,3-d]pyrimidin-7-yl]-N'-(1,1-dimethylethyl)-urea

Structure:

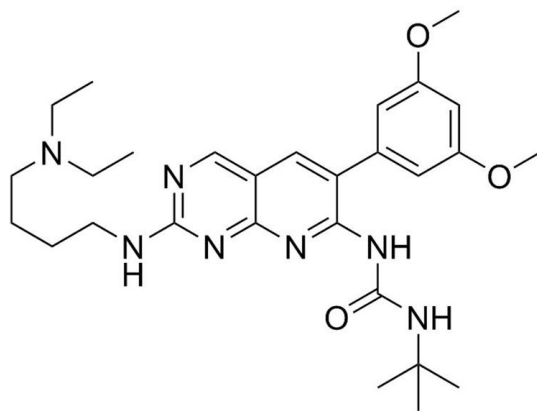

## Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt.
Solubility:	• DMSO ≤ 20 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 190 μL of fresh DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### MAINTENANCE AND SELF-RENEWAL

- Suppresses the differentiation of mouse embryonic stem (ES) cells and maintains the undifferentiated state (Kunath et al.; Ying et al.).

### REPROGRAMMING

- Prevents excision-mediated differentiation of mouse induced pluripotent stem cells generated using piggyBac transposons (Kaji et al.).
- Promotes reprogramming of human ES cells to naïve cells, or their maintenance in a naïve state, in combination with OCT4, KLF4, KLF2, LIF (Catalog #78055), CHIR99021 (Catalog #72052), and PD0325901 (Catalog #72182) (Hanna et al.).

### DIFFERENTIATION

- Blocks neural differentiation of mouse ES cells (Stavridis et al.).
- Promotes differentiation of human ES cells, but not when they are in a naïve or ground state (Hanna et al.).

## References

- Hanna J et al. (2010) Human embryonic stem cells with biological and epigenetic characteristics similar to those of mouse ESCs. *Proc Natl Acad Sci USA* 107(20): 9222–7.
- Kaji K et al. (2009) Virus-free induction of pluripotency and subsequent excision of reprogramming factors. *Nature* 458(7239): 771–5.
- Koziczak M et al. (2004) Blocking of FGFR signaling inhibits breast cancer cell proliferation through downregulation of D-type cyclins. *Oncogene* 23(20): 3501–8.
- Kunath T et al. (2007) FGF stimulation of the Erk1/2 signalling cascade triggers transition of pluripotent embryonic stem cells from self-renewal to lineage commitment. *Development* 134(16): 2895–902.
- Mohammadi M et al. (1998) Crystal structure of an angiogenesis inhibitor bound to the FGF receptor tyrosine kinase domain. *EMBO J* 17(20): 5896–904.
- Stavridis MP et al. (2007) A discrete period of FGF-induced Erk1/2 signalling is required for vertebrate neural specification. *Development* 134(16): 2889–94.
- Trudel S et al. (2004) Inhibition of fibroblast growth factor receptor 3 induces differentiation and apoptosis in t(4;14) myeloma. *Blood* 103(9): 3521–8.
- Ying Q-L et al. (2008) The ground state of embryonic stem cell self-renewal. *Nature* 453(7194): 519–23.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

**This product is hazardous. Please refer to the Safety Data Sheet (SDS).**

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2022 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.