

Small Molecules

Thiazovivin

RHO/ROCK pathway inhibitor; Inhibits ROCK

Catalog #	72252	1 mg
	72254	5 mg
	100-0247	25 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Thiazovivin is a selective inhibitor of Rho-associated coiled-coil containing protein kinase (ROCK), a serine/threonine kinase that plays a role in cell polarity, contraction, and actin cytoskeleton reorganization. Thiazovivin is effective at 5-fold-lower concentrations than another common ROCK inhibitor Y-27632 (Catalog #72302; Xu et al.).

Molecular Name:	Thiazovivin
Alternative Names:	Tzv
CAS Number:	1226056-71-8
Chemical Formula:	C ₁₅ H ₁₃ N ₅ OS
Molecular Weight:	311.4 g/mol
Purity:	≥ 98%
Chemical Name:	N-(phenylmethyl)-2-(4-pyrimidinylamino)-4-thiazolecarboxamide
Structure:	

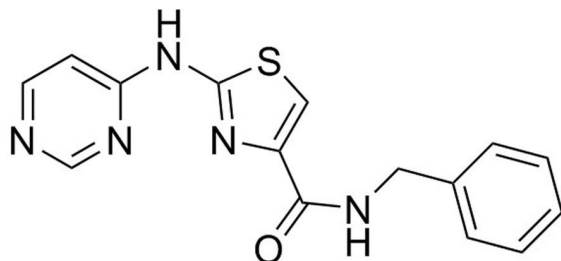

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt.
Solubility:	· DMSO ≤ 40 mM · Absolute ethanol ≤ 3.2 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 321 µL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

MAINTENANCE AND SELF-RENEWAL

- Promotes survival of human embryonic stem (ES) cells during dissociation by stabilizing E-cadherin and improves cell attachment (Xu et al.).
- Promotes survival of single human induced pluripotent stem (iPS) cells during transfection for TALEN-mediated genome editing (Sun and Zhao).

REPROGRAMMING

- Increases the efficiency of reprogramming human somatic cells to iPS cells, in combination with PD0325091 (Catalog #72182) and SB431542 (Catalog #72232) (Lin et al.).
- Increases the efficiency of reprogramming human cord blood mononuclear cells to iPS cells (Hu et al.). Promotes survival of human embryonic stem (ES) cells during dissociation by stabilizing E-cadherin and improves cell attachment (Xu et al.).

References

- Hu K et al. (2011) Efficient generation of transgene-free induced pluripotent stem cells from normal and neoplastic bone marrow and cord blood mononuclear cells. *Blood* 117(14): e109–19.
- Lin T et al. (2009) A chemical platform for improved induction of human iPSCs. *Nat Methods* 6(11): 805–8.
- Sun N & Zhao H. (2014) Seamless correction of the sickle cell disease mutation of the HBB gene in human induced pluripotent stem cells using TALENs. *Biotechnol Bioeng* 111(5): 1048–53.
- Xu Y et al. (2010) Revealing a core signaling regulatory mechanism for pluripotent stem cell survival and self-renewal by small molecules. *Proc Natl Acad Sci USA* 107(18): 8129–34.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2020 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.