

Small Molecules

StemRegenin 1 (Hydrochloride)

Aryl hydrocarbon receptor (AHR) antagonist

Catalog # 72352
72354

1 mg
5 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

StemRegenin1 (SR1) is an antagonist of the aryl hydrocarbon receptor (AhR). It promotes ex vivo expansion of CD34+ human hematopoietic stem cells (Boitano et al.; Csaszar et al.) and the generation of CD34+ hematopoietic progenitor cells from non-human primate induced pluripotent stem (iPS) cells (Gori et al.). SR1 has been shown to collaborate with UM729 (Catalog #72332) in preventing differentiation of acute myeloid leukemia (AML) cells in culture (Pabst et al.). SR1 also stimulates the proliferation and differentiation of CD34+ hematopoietic progenitor cells into dendritic cells (Thordardottir et al.).

Molecular Name:	StemRegenin 1 (Hydrochloride)
Alternative Names:	SR1 (Hydrochloride)
CAS Number:	N/A
Chemical Formula:	C ₂₄ H ₂₃ N ₅ OS · HCl
Molecular Weight:	466.0 g/mol
Purity:	≥ 98%
Chemical Name:	4-(2-((2-(benzo[b]thiophen-3-yl)-9-isopropyl-9H-purin-6-yl)amino)ethyl) phenol, monohydrochloride
Structure:	

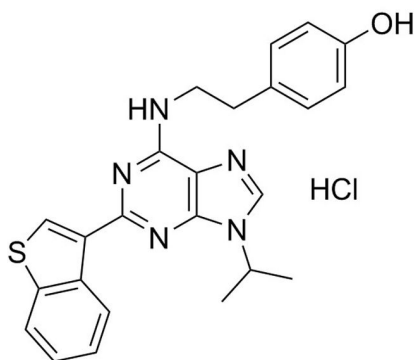

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt.
Solubility:	· DMSO ≤ 10 mM For example, to prepare a 5 mM stock solution in DMSO, resuspend 1 mg in 429 µL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

This compound has low solubility in aqueous solution. However, aqueous concentrations up to 1 µM can be achieved by dilution from a DMSO stock solution; ensure complete solubilization prior to use as a cell culture supplement. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

MAINTENANCE AND SELF-RENEWAL

- Promotes maintenance and expansion of human hematopoietic stem cells in culture (Boitano et al.; Csaszar et al.).

DIFFERENTIATION

- Stimulates differentiation of CD34+ hematopoietic progenitor cells into functional human dendritic cells (Thordardottir et al.).
- Promotes hematopoietic differentiation of iPS cells (Gori et al.).

CANCER RESEARCH

- Collaborates with UM729 in preventing differentiation of AML cells in culture (Pabst et al.).

References

Boitano AE et al. (2010) Aryl hydrocarbon receptor antagonists promote the expansion of human hematopoietic stem cells. *Science* 329(5997): 1345–8.

Csaszar E et al. (2012) Rapid expansion of human hematopoietic stem cells by automated control of inhibitory feedback signaling. *Cell Stem Cell* 10(2): 218–29.

Gori JL et al. (2012) Efficient generation, purification, and expansion of CD34(+) hematopoietic progenitor cells from nonhuman primate-induced pluripotent stem cells. *Blood* 120(13): e35–44.

Pabst C et al. (2014) Identification of small molecules that support human leukemia stem cell activity ex vivo. *Nat Methods* 11(4): 436–42.

Thordardottir S et al. (2014) The aryl hydrocarbon receptor antagonist StemRegenin 1 promotes human plasmacytoid and myeloid dendritic cell development from CD34+ hematopoietic progenitor cells. *Stem Cells Dev* 23(9): 955–67.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2021 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.