

Small Molecules

Quercetin

mTOR, PI3K/AKT, NF- κ B, and tyrosine kinase pathway inhibitor; Inhibits PI3K and SRC kinases

Catalog # 73932
73934

5 g
10 g



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Quercetin inhibits phosphoinositide 3-kinase isoforms (beta: $IC_{50} = 5.4 \mu\text{M}$; gamma: $IC_{50} = 2.4 \mu\text{M}$; delta: $IC_{50} = 3.0 \mu\text{M}$) and SRC kinases but only shows slight inhibitory activity against serine/threonine kinases ERK1/2 and p38 alpha, as well as protein kinase C (Navarro-Núñez et al.). This product is supplied as the hydrate form of the molecule.

Molecular Name:	Quercetin (Hydrate)
Alternative Names:	Not applicable
CAS Number:	849061-97-8
Chemical Formula:	$C_{15}H_{10}O_7 \cdot XH_2O$
Molecular Weight:	302.2 g/mol
Purity:	$\geq 95\%$
Chemical Name:	2-(3,4-dihydroxyphenyl)-3,5,7-trihydroxy-4H-1-benzopyran-4-one
Structure:	

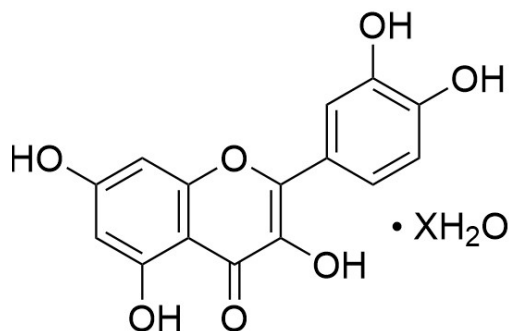

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at room temperature (15 - 25°C) as supplied. Protect from prolonged exposure to light. For long-term storage store with a desiccant. Stable as supplied for 12 months from date of receipt.
Solubility:	<ul style="list-style-type: none">· DMSO ≤ 95 mM· Absolute ethanol ≤ 6 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 10 mg in 3.31 mL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

REPROGRAMMING

- Stimulates glycolytic metabolism and promotes reprogramming of somatic cells (Zhang et al.).
- Enhances reprogramming of human primary somatic cells to induced pluripotent stem cells (Zhu et al.).

CANCER RESEARCH

- Induces hypoxia-inducible factor-1 α (HIF-1 α) and inhibits cell proliferation (Triantafyllou et al.).
- Promotes apoptosis of glioblastoma cells when used together with temozolomide (Sang et al.).
- Reduces tumorigenicity of pancreatic cancer cells (Fan et al.).

References

- Fan P et al. (2016) Continuous exposure of pancreatic cancer cells to dietary bioactive agents does not induce drug resistance unlike chemotherapy. *Cell Death Dis* 7(6): e2246.
- Navarro-Núñez L et al. (2010) Effect of quercetin on platelet spreading on collagen and fibrinogen and on multiple platelet kinases. *Fitoterapia* 81(2): 75–80.
- Sang D et al. (2014) Quercetin sensitizes human glioblastoma cells to temozolomide in vitro via inhibition of Hsp27. *Acta Pharmacol Sin* 35(6): 832–8.
- Triantafyllou A et al. (2007) The flavonoid quercetin induces hypoxia-inducible factor-1 α (HIF-1 α) and inhibits cell proliferation by depleting intracellular iron. *Free Radic Res* 41(3): 342–56.
- Zhang Y et al. (2012) Small molecules, big roles -- the chemical manipulation of stem cell fate and somatic cell reprogramming. *J Cell Sci* 125(Pt 23): 5609–20.
- Zhu S et al. (2010) Reprogramming of human primary somatic cells by OCT4 and chemical compounds. *Cell Stem Cell* 7(6): 651–5.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2020 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.