

Small Molecules

NU7441

NHEJ pathway inhibitor; Inhibits DNA-dependent protein kinase (DNA-PK)

Catalog # 74082
74084

5 mg
10 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

NU7441 is an inhibitor of DNA-dependent protein kinase (DNA-PK), an enzyme involved in the non-homologous end joining (NHEJ) DNA repair pathway. It is highly selective for DNA-PK, with an IC_{50} of 14 nM (Leahy et al.).

Molecular Name:	NU7441
Alternative Names:	KU57788
CAS Number:	503468-95-9
Chemical Formula:	$C_{25}H_{19}NO_3S$
Molecular Weight:	413.5 g/mol
Purity:	≥ 98%
Chemical Name:	8-(4-dibenzothiienyl)-2-(4-morpholinyl)-4H-1-benzopyran-4-one
Structure:	

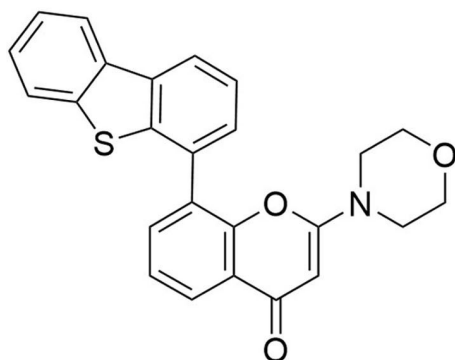

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. For long-term storage, store with a desiccant. Stable as supplied for 12 months from date of receipt.
Solubility:	· DMSO ≤ 480 μM · DMF ≤ 2.4 mM For example, to prepare a 1 mM stock solution in DMF, resuspend 1 mg in 2.42 mL of DMF.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMF at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMF concentration above 0.1% due to potential cell toxicity.

Published Applications

GENOME EDITING

- Reduces the frequency of NHEJ and increases the efficiency of homology-directed repair (HDR) in CRISPR-Cas9 genome editing (Robert et al.).

CANCER RESEARCH

- Sensitizes human cancer cell lines to DNA double-strand-break-inducing therapies (chemo- or radio-therapy) by inhibiting DNA-PK activity and delaying the repair of double-strand breaks (Ciszewski et al.; Shaheen et al.; Yang et al.; Zhao et al.).

References

Ciszewski WM et al. (2014) DNA-PK inhibition by NU7441 sensitizes breast cancer cells to ionizing radiation and doxorubicin. *Breast Cancer Res Treat* 143(1): 47–55.

Leahy JJ et al. (2004) Identification of a highly potent and selective DNA-dependent protein kinase (DNA-PK) inhibitor (NU7441) by screening of chromenone libraries. *Bioorg Med Chem Lett* 14(24): 6083–7.

Robert F et al. (2015) Pharmacological inhibition of DNA-PK stimulates Cas9-mediated genome editing. *Genome Med* 7(1): 93.

Shaheen FS et al. (2011) Targeting the DNA double strand break repair machinery in prostate cancer. *PLoS One* 6(5): e20311.

Yang C et al. (2016) NU7441 enhances the radiosensitivity of liver cancer cells. *Cell Physiol Biochem* 38(5): 1897–905.

Zhao Y et al. (2006) Preclinical evaluation of a potent novel DNA-dependent protein kinase inhibitor NU7441. *Cancer Res* 66(10): 5354–62.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2021 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.