

Cytokines

Human/Mouse Recombinant NT-3



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Neurotrophin-3

| | | |
|-----------|---------|---------|
| Catalog # | 78074 | 10 µg |
| | 78074.1 | 100 µg |
| | 78074.2 | 1000 µg |

Product Description

Neurotrophin-3 (NT-3) is a neurotrophic factor and a member of the nerve growth factor (NGF) family of proteins that includes neuron growth factor (NGF), brain-derived neurotrophic factor (BDNF), and neurotrophin-4/5. NT-3 signals a number of trophic effects through its transducing receptor tyrosine kinase TrkC. NT-3 is known to promote survival, development, and differentiation of neurons, and modulates transmitter release at several types of synapses in the peripheral and central nervous systems (Chalazonitis 1996). NT-3 has been shown to have an important role in the overall development of enteric neurons, which are crucial for gut peristalsis (Chalazonitis 2004). Studies in rats have shown the potential of NT-3 in dorsal column axonal regeneration (Bradbury et al.). NT-3 was shown to protect neurons against amyloid- β toxicity (Lesne et al.). NT-3 has applications in neuronal differentiation protocols to generate β -tubulin III+ peripheral neurons from neural crest stem cells (Menendez et al.) and oligodendrocyte precursor cells from human embryonic stem (ES) and induced pluripotent stem (iPS) cells (Douvaras et al.).

Product Information

| | |
|---------------------------|--|
| Alternative Names: | HGNF, NGF-2, Nerve growth factor 2, Neurotrophin 3, NT3 |
| Accession Number: | P20783 |
| Amino Acid Sequence: | MYAEHKSHRG EYSVCDSESL WVTDKSSAID IRGHQVTVLG EIKTGNSPVK QYFYETRCKE ARPVKNGCRG IDDKHWNSQC KTSQTYVRAL TSENNKLVGWRW IRIDTSCV CALSRKIGRT |
| Predicted Molecular Mass: | 27.5 kDa |
| Species: | Human, Mouse |
| Formulation: | Lyophilized from a sterile-filtered aqueous solution containing 0.1% trifluoroacetic acid. |
| Source: | <i>E. coli</i> |

Specifications


| | |
|------------------|--|
| Activity: | Lot #1000050644 or higher: The specific activity is $\geq 5.0 \times 10^4$ units/mg ($EC_{50} \leq 20$ ng/mL) as determined by a cell proliferation assay using a neuroblastoma cell line stably expressing TrkB (BR6). All other lots: The specific activity is ≥ 50 units/mg ($EC_{50} \leq 20$ µg/mL) as determined by a cell proliferation assay using C6 cells. |
| Purity: | $\geq 95\%$ |
| Endotoxin Level: | Measured by kinetic Limulus amoebocyte lysate (LAL) analysis and is ≤ 1 EU/µg protein. |

Preparation and Storage

| | |
|--------------|--|
| Storage: | Store at -20°C to -80°C. |
| Stability: | Stable as supplied for 12 months from date of receipt. |
| Preparation: | Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. |

OPTIONAL: After reconstitution, if product will not be used immediately, dilute with concentrated bovine serum albumin (BSA) to a final BSA concentration of 0.1%. The effect of storage of stock solution on product performance should be tested for each application. As a general guide, do not store at 2 - 8°C for more than 1 month or at -80°C for more than 3 months. Avoid repeated freeze-thaw cycles.

Data

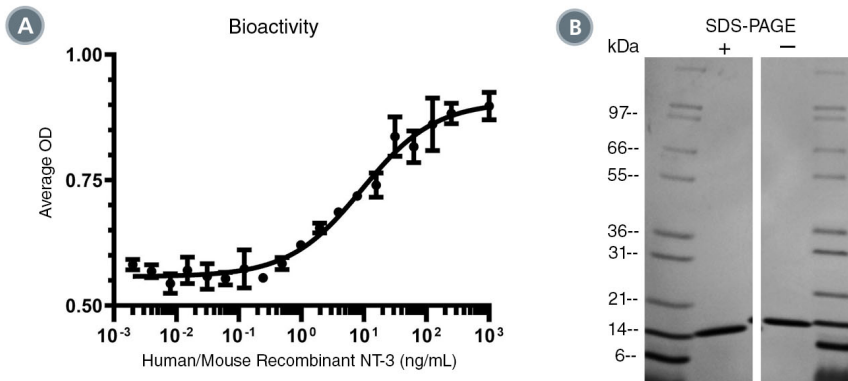

(A) The biological activity of Human/Mouse Recombinant NT-3 was tested by its ability to promote the proliferation of a neuroblastoma cell line stably expressing TrkB (BR6). Cell proliferation was measured using a colorimetric assay method. The EC₅₀ is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC₅₀ in the above example is 9.3 ng/mL. For a representative bioactivity plot of C6 cell proliferation, contact us at techsupport@stemcell.com.

(B) 1 µg of Human/Mouse Recombinant NT-3 was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human/Mouse Recombinant NT-3 is predicted to be a non-covalently linked homodimer, with each dimer having a molecular mass of 13.8 kDa for a predicted total molecular weight of 27.5 kDa.

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, visit www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Bradbury EJ et al. (1999) NT-3 promotes growth of lesioned adult rat sensory axons ascending in the dorsal columns of the spinal cord. *Eur J Neurosci* 11(11): 3873–83.
- Chalazonitis A. (1996) Neurotrophin-3 as an essential signal for the developing nervous system. *Mol Neurobiol* 12(1): 39–53.
- Chalazonitis A. (2004) Neurotrophin-3 in the development of the enteric nervous system. *Prog Brain Res* 146: 243–63.
- Douvaras P et al. (2014) Efficient generation of myelinating oligodendrocytes from primary progressive multiple sclerosis patients by induced pluripotent stem cells. *Stem cell reports* 3(2): 250–9.
- Lesne S et al. (2005) Akt-dependent expression of NAIP-1 protects neurons against amyloid-beta toxicity. *J Biol Chem* 280(26): 24941–7.
- Menendez L et al. (2013) Directed differentiation of human pluripotent cells to neural crest stem cells. *Nat Protoc* 8(1): 203–12.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2021 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.