

Human Recombinant IL-7, ACF

Interleukin 7, animal component-free

Catalog #78196	10 µg
Catalog #78196.1	100 µg
Catalog #78196.2	1000 µg

Product Description

Interleukin 7 (IL-7) is a member of the type I cytokine family that is critical for T and B cell development and survival. It is produced by non-hematopoietic cells in the thymus, lymphoid organs, and by bone marrow stromal cells (Lundström et al.). IL-7 binds to a receptor (IL-7R) composed of common gamma chain and IL-7Ra (CD127) and signals through the JAK/STAT and PI3K pathways. IL-7 regulates the survival of naïve and memory CD4+ and CD8+ T cells, $\gamma\delta$ T cells, NK T cells, innate lymphoid cells, and regulatory T cells (Carrette & Surh). Although a deficiency in IL-7R still permits the generation of normal numbers of peripheral B cells in humans, stimulation of human B cell precursors with IL-7 could promote STAT5-dependent proliferation and survival in vitro (Clark et al.; Corfe & Paige). This product is animal component-free.

Product Information

Alternative Names:	Interleukin-7, LP-1, Lymphopoietin 1, Pre-B cell factor
Accession Number:	P13232
Amino Acid Sequence:	MDCDIEGKDG KQYESVLMVS IDQLLDSMKE IGSNCLNNEF NFFKRHICDA NKEGMFLFRA ARKLRQFLKM NSTGDFDLHL LKVSEGTTL LNCTGQVKGR KPAALGEAQP TKSLLENKSL KEQKLLNDLC FLKRLLEIK TCWNKILMGT KEH
Predicted Molecular Mass:	17.5 kDa
Species:	Human
Product Formulation:	Lyophilized from a sterile-filtered aqueous solution containing sodium phosphate, pH 7.5.
Source:	E. coli
Purity:	≥ 97%


Specifications

Activity:	The specific activity is $\geq 2 \times 10^6$ units/mg ($EC_{50} \leq 0.5$ ng/mL), as determined by a cell proliferation assay using phytohemagglutinin (PHA)-stimulated human peripheral blood mononuclear cells (PBMCs).
Endotoxin Level:	Measured by kinetic Limulus amoebocyte lysate (LAL) analysis and is ≤ 1 EU/ μ g protein.

Preparation and Storage

Stability and Storage:	Store at -20 to -80°C. Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. OPTIONAL: After reconstitution, if product will not be used immediately, dilute with concentrated bovine serum albumin (BSA) to a final BSA concentration of 0.1%. The effect of storage of stock solution on product performance should be tested for each application. As a general guide, do not store at 2 - 8°C for more than 1 month or at -80°C for more than 3 months. Avoid repeated freeze-thaw cycles.

Data

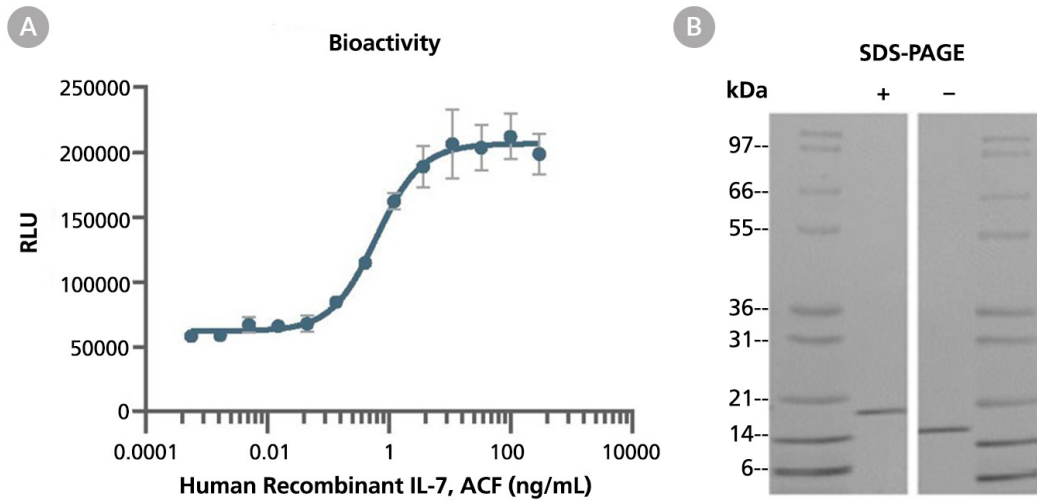

(A) The biological activity of Human Recombinant IL-7, ACF was tested by its ability to promote the proliferation of PHA-stimulated PBMCs. The EC50 is defined as the effective concentration of the cytokine at which cell proliferation is at 50% of maximum. The EC50 in the above example is ≤ 0.5 ng/mL. (B) 1 μ g of Human Recombinant IL-7, ACF was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human Recombinant IL-7, ACF has a predicted molecular mass of 17.5 kDa.

Related Products

For a complete list of cytokines or peptide pools, as well as related products available from STEMCELL Technologies, visit www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Carrette F & Surh CD. (2012) IL-7 signaling and CD127 receptor regulation in the control of T cell homeostasis. *Semin Immunol* 24(3): 209–17.
- Clark MR et al. (2014) Orchestrating B cell lymphopoiesis through interplay of IL-7 receptor and pre-B cell receptor signalling. *Nat Rev Immunol* 14(2): 69–80.
- Corfe SA & Paige CJ. (2012) The many roles of IL-7 in B cell development; mediator of survival, proliferation and differentiation. *Semin Immunol* 24(3): 198–208.
- Lundström W et al. (2012) IL-7 in human health and disease. *Semin Immunol* 24(3): 218–24.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2024 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.