

RosetteSep™ DM-M Density Medium



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713
INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM
FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #15725 100 mL

Document #10000003344 | Version 02

Intended Use

RosetteSep™ DM-M Density Medium is a density separation medium designed specifically for use as part of the RosetteSep™ HLA Myeloid Cell Enrichment Kit (Catalog #15272HLA) for the enrichment of human myeloid (CD33+) cells from whole blood.

Product Specifications

RosetteSep™ DM-M Density Medium (Catalog #15725)

Density: 1.085 g/mL

Quality Control

RosetteSep™ DM-M Density Medium is manufactured using aseptic technique and tightly controlled processes. Each lot of RosetteSep™ DM-M Density Medium is sterility tested according to USP methods.

Storage and Stability

Store at 15 - 25°C. Storage at 2 - 8°C is acceptable, but ensure that the medium equilibrates to 15 - 25°C and invert bottle to mix contents before use. Keep protected from direct light. Product stable at 15 - 25°C until expiry date as indicated on label.

Warnings and Precautions

1. For professional users only.
2. This product is for in vitro diagnostic use.
3. Do not use bottle if contents have leaked. Unused bottle may be disposed of according to standard laboratory procedures for non-hazardous liquids.
4. This product should be handled by trained personnel observing good laboratory practices. Once this product is added to human cells, treat the suspension as potentially biohazardous. Handling of reagents and disposal of wastes should observe all local, state, or national regulations.
5. This product is a potential irritant to eyes, respiratory system, and skin. This product may also be harmful if ingested. Avoid exposure through skin, eye contact, inhalation, and ingestion.

Directions for Use

For detailed instructions, refer to the Product Information Sheet for RosetteSep™ HLA Myeloid Cell Enrichment Kit (Document #10000003559).

NOTE: Do not use RosetteSep™ DM-M Density Medium with any other RosetteSep™ enrichment kits. The density of RosetteSep™ DM-M Density Medium is not optimized for the isolation of other cell types.



STEMCELL Technologies Canada Inc. | 1618 Station Street, Vancouver, BC | V6A 1B6 | Canada | www.stemcell.com

For Technical Assistance

Tel: +1.604.877.0713

European toll-free-number: 00800 7836 2355

E-mail: techsupport@stemcell.com



MDSS GmbH

Schiffgraben 41

30175 Hannover • Germany



Document #10000003344









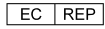
Version 02

2021

Technical Assistance

For technical support, contact us by email at techsupport@stemcell.com or call either +1.604.877.0713 or the European toll-free number 00800 7836 2355. For more information, visit www.stemcell.com.

If you require a printed copy or a translated version of this document in a certain language, contact technical support.

 Catalog or reference number	 Batch code	 Use by: YYYY-MM
 Caution, consult accompanying documents	 In Vitro Diagnostic Medical Device	 For storage within temperature limits
 CE Mark	 Manufacturer's identification (name & address)	 Authorized EC representative in the European Community

Copyright © 2021 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, and RosetteSep are trademarks of STEMCELL Technologies Canada Inc. Lymphoprep is a trademark of Alere Technologies. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.