

EasySep™ Mouse/Human Chimera Isolation Kit

For processing 1 x 10⁹ cells

Catalog #19849

Negative Selection

Document #1000003743 | Version 02



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Description

Isolate untouched and highly purified human cells from mouse splenocytes, bone marrow, or peripheral blood by immunomagnetic negative selection. When using single-cell suspensions from other tissue types, this kit may require optimization.

- Fast, easy-to-use, and column-free
- Isolated cells are untouched

This kit targets hematopoietic mouse cells for removal with biotinylated antibodies recognizing specific cell surface markers. Unwanted cells are labeled with biotinylated antibodies and streptavidin-coated magnetic particles and separated without columns using an EasySep™ magnet. Desired cells are simply poured off into a new tube. Isolated cells are immediately available for downstream applications, such as flow cytometry, culture, or DNA/RNA extraction.

Component Descriptions

COMPONENT NAME	COMPONENT #	QUANTITY	STORAGE	SHELF LIFE	FORMAT
EasySep™ Mouse/Human Chimera Isolation Cocktail	19849C	1 x 0.5 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A combination of monoclonal antibodies in PBS and 0.1% BSA.
EasySep™ Streptavidin RapiSep™ 50001	50001	1 x 1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A suspension of magnetic particles in PBS.
EasySep™ Mouse FcR Blocker	18720	1 x 0.1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A combination of monoclonal antibodies in PBS, 0.1% BSA, and < 0.1% sodium azide.

BSA - bovine serum albumin; PBS - phosphate-buffered saline

Components may be shipped at room temperature (15 - 25°C) but should be stored as indicated above.

Sample Preparation

BONE MARROW

Flush bone marrow cells from femur and tibia into recommended medium using a syringe equipped with a 23 gauge needle. Disperse aggregates by gently passing the cell suspension through the syringe several times. Alternatively, crush bones using a mortar and pestle. Remove remaining aggregates and debris by passing cell suspension through a 70 µm mesh nylon strainer. Centrifuge at 300 x g for 10 minutes and resuspend cells at 1 x 10⁸ nucleated cells/mL in recommended medium.

Ammonium chloride treatment is not recommended when preparing the cells for separation.

SPLEEN

Disrupt spleen in PBS or Hanks' Balanced Salt Solution (HBSS) containing 2% fetal bovine serum (FBS). Remove aggregates and debris by passing cell suspension through a 70 µm mesh nylon strainer (e.g. Catalog #27216). Centrifuge at 300 x g for 10 minutes and resuspend at 1 x 10⁸ nucleated cells/mL in recommended medium.

Ammonium chloride treatment is not recommended when preparing the cells for separation.

PERIPHERAL BLOOD

Blood should be lysed prior to use. Mix 1 part blood with 9 parts Ammonium Chloride Solution (Catalog #07800) and incubate on ice for 15 minutes. Centrifuge at 300 x g for 6 minutes. Discard supernatant and wash cell pellet once with recommended medium. Discard supernatant and resuspend cell pellet at 1 x 10⁸ nucleated cells/mL in recommended medium.



Recommended Medium

EasySep™ Buffer (Catalog #20144), RoboSep™ Buffer (Catalog #20104), or PBS (Catalog #37350) containing 2% FBS and 1 mM EDTA. HBSS (Catalog #37250) can be used in place of PBS. Medium should be free of Ca++, Mg++, and biotin.

Directions for Use – Manual EasySep™ Protocols

See page 1 for Sample Preparation and Recommended Medium. Refer to Tables 1 and 2 for detailed instructions regarding the EasySep™ procedure for each magnet.




Table 1. EasySep™ Mouse/Human Chimera Isolation Kit Protocol

		EASYSEP™ MAGNETS	
STEP	INSTRUCTIONS	 EasySep™ (Catalog #18000)	“The Big Easy” (Catalog #18001) 
1	Prepare sample at the indicated cell concentration within the volume range.	1 x 10 ⁸ cells/mL 0.25 - 2 mL	1 x 10 ⁸ cells/mL 0.5 - 8 mL
2	Add FcR blocker to sample.	10 µL/mL of sample	10 µL/mL of sample
3	Add sample to required tube.	5 mL (12 x 75 mm) polystyrene round-bottom tube (e.g. Catalog #38007)	14 mL (17 x 95 mm) polystyrene round-bottom tube (e.g. Catalog #38008)
4	Add Isolation Cocktail to sample. NOTE: Do not vortex cocktail.	50 µL/mL of sample	50 µL/mL of sample
	Mix and incubate.	2 - 8°C for 10 minutes	2 - 8°C for 10 minutes
5	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds	30 seconds
6	Add RapidSpheres™ to sample.	75 µL/mL of sample	75 µL/mL of sample
	Mix and incubate.	2 - 8°C for 5 minutes	2 - 8°C for 5 minutes
7	Add recommended medium to top up the sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 2.5 mL	<ul style="list-style-type: none"> • Top up to 5 mL for samples ≤ 4 mL • Top up to 10 mL for samples > 4 mL
	Place the tube (without lid) into the magnet and incubate.	<ul style="list-style-type: none"> • For bone marrow, spleen: RT for 5 minutes • For blood: RT for 10 minutes 	<ul style="list-style-type: none"> • For bone marrow, spleen: RT for 5 minutes • For blood: RT for 10 minutes
8	Pick up the magnet, and in one continuous motion invert the magnet and tube,* pouring the enriched cell suspension into a new tube.	Use new 5 mL tube	Use new 14 mL tube
9	Remove the original tube from the magnet. Place the new tube (without lid) into the magnet and incubate for a second separation.	RT for 5 minutes	RT for 5 minutes
10	Pick up the magnet, and in one continuous motion invert the magnet and tube,* pouring the enriched cell suspension into a new tube.	Isolated cells are ready for use	Isolated cells are ready for use

RT - room temperature (15 - 25°C)

* Leave the magnet and tube inverted for 2 - 3 seconds, then return upright. Do not shake or blot off any drops that may remain hanging from the mouth of the tube.

Table 2. EasySep™ Mouse/Human Chimera Isolation Kit Protocol

STEP	INSTRUCTIONS	EASYSEP™ MAGNETS		
		 EasyPlate™ (Catalog #18102)	 EasyEight™ (Catalog #18103) 5 mL tube	 14 mL tube
1	Prepare sample at the indicated cell concentration within the volume range.	1 x 10 ⁸ cells/mL 0.05 - 0.2 mL	1 x 10 ⁸ cells/mL 0.25 - 2 mL	1 x 10 ⁸ cells/mL 0.5 - 8 mL
2	Add FcR blocker to sample.	10 µL/mL of sample	10 µL/mL of sample	10 µL/mL of sample
3	Add sample to required tube (or plate when using the EasyPlate™ EasySep™ Magnet).	Round-bottom, non-tissue culture-treated 96-well plate (e.g. Catalog #38018)	5 mL (12 x 75 mm) polystyrene round-bottom tube (e.g. Catalog #38007)	14 mL (17 x 95 mm) polystyrene round-bottom tube (e.g. Catalog #38008)
4	Add Isolation Cocktail to sample. NOTE: Do not vortex cocktail.	50 µL/mL of sample	50 µL/mL of sample	50 µL/mL of sample
	Mix and incubate.	2 - 8°C for 10 minutes	2 - 8°C for 10 minutes	2 - 8°C for 10 minutes
5	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds	30 seconds	30 seconds
6	Add RapidSpheres™ to sample.	75 µL/mL of sample	75 µL/mL of sample	75 µL/mL of sample
	Mix and incubate.	2 - 8°C for 5 minutes	2 - 8°C for 5 minutes	2 - 8°C for 5 minutes
7	Add recommended medium to top up sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 0.25 mL	Top up to 2.5 mL	<ul style="list-style-type: none"> • Top up to 5 mL for samples ≤ 4 mL • Top up to 10 mL for samples > 4 mL
	Place the tube or plate (without lid) into the magnet and incubate.	RT for 10 minutes	<ul style="list-style-type: none"> • For bone marrow, spleen: RT for 5 minutes • For blood: RT for 10 minutes 	<ul style="list-style-type: none"> • For bone marrow, spleen: RT for 5 minutes • For blood: RT for 10 minutes
8	Carefully pipette** (do not pour) the enriched cell suspension into a new tube or plate.	Use new 96-well plate	Use new 5 mL tube	Use new 14 mL tube
9	Remove the tube or plate from the magnet. Place the new tube or plate (without lid) into the magnet and incubate for a second separation.	RT for 5 minutes	RT for 5 minutes	RT for 5 minutes
10	Carefully pipette** (do not pour) the enriched cell suspension into a new tube or plate.	Isolated cells are ready for use	Isolated cells are ready for use	Isolated cells are ready for use

RT - room temperature (15 - 25°C)

** Collect the entire supernatant, all at once, into a single pipette (e.g. for EasyEight™ 5 mL tube, use a 2 mL serological pipette [Catalog #38002]; for EasyEight™ 14 mL tube, use a 10 mL serological pipette [Catalog #38004]).

Notes and Tips

ASSESSING PURITY


For purity assessment of human cells by flow cytometry, use the following fluorochrome-conjugated antibody clone:

- Anti-Human CD45 Antibody, Clone HI30 (Catalog #60018)

For purity assessment of mouse cells by flow cytometry, use the following fluorochrome-conjugated antibody clones:

- Anti-Mouse CD45 Antibody, Clone 30-F11 (Catalog #60030), and
- Anti-Mouse TER119 Antibody, Clone TER-119 (Catalog #60033)

Data

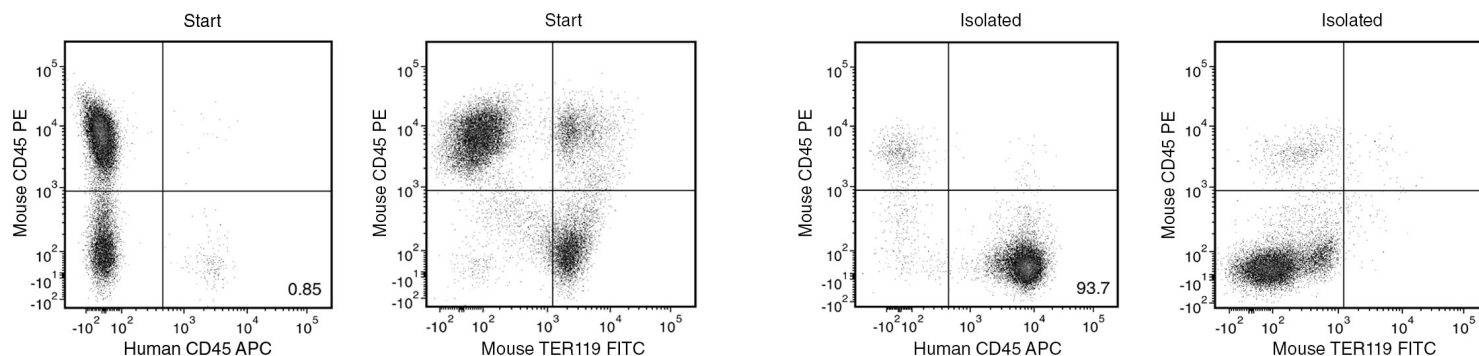

In the above example, human peripheral blood mononuclear cells (PBMCs) were seeded into mouse bone marrow at a frequency of 0.85%. The human PBMC content of the isolated fraction was 93.7%. The depletion of CD45+ mouse cells is typically 3 - 4.5 log.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED. FOR ADDITIONAL INFORMATION ON QUALITY AT STEMCELL, REFER TO WWW.STEMCELL.COM/COMPLIANCE.

Copyright © 2023 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, EasyEights, EasyPlate, EasySep, RapidSpheres, and RoboSep are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.