EasySep™ Human B Cell Isolation Kit

For processing 1 x 10⁹ cells

Catalog #17954

Negative Selection

Document #10000005312 | Version 02



Scientists Helping Scientists™ | www.stemcell.com

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713 INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Description

Isolate untouched and highly purified B cells from fresh or previously frozen human peripheral blood mononuclear cells (PBMCs) or washed leukapheresis samples in as little as 9 minutes by immunomagnetic negative selection.

- · Fast, easy-to-use, and column-free
- · Up to 96% purity with high recovery
- · Isolated cells are untouched

This kit targets non-B cells for removal with antibodies recognizing specific surface markers. Unwanted cells are labeled with antibodies and magnetic particles, and separated without columns using an EasySep™ magnet. Desired cells are simply poured off into a new tube. Isolated cells are immediately available for downstream applications such as flow cytometry, culture, or DNA/RNA extraction.

Component Descriptions

COMPONENT NAME	COMPONENT #	QUANTITY	STORAGE	SHELF LIFE	FORMAT
EasySep™ Human B Cell Isolation Cocktail	17954C	1 x 1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A combination of monoclonal antibodies in PBS.
EasySep™ Isolation Cocktail Enhancer	17900	1 x 1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A solution that enhances the performance of the isolation cocktail.
EasySep™ Dextran RapidSpheres™ 50102	50102	1 x 1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A suspension of magnetic particles in water.

PBS - phosphate-buffered saline

Components may be shipped at room temperature (15 - 25°C) but should be stored as indicated above.

Sample Preparation

For available fresh and frozen samples, see www.stemcell.com/primarycells.

PERIPHERAL BLOOD

Prepare a PBMC suspension from whole blood by centrifugation over a density gradient medium (e.g. Lymphoprep[™], Catalog #07801). For more rapid PBMC preparation, use the SepMate[™] RUO (Catalog #86450/86415) or SepMate[™] IVD* (Catalog #85450/85415) cell isolation tube.

If using previously frozen PBMCs, incubate the cells with DNase I Solution (1 mg/mL; Catalog #07900) at a concentration of 100 μg/mL at room temperature (15 - 25°C) for at least 15 minutes prior to labeling and separation. Filter agreggated suspensions through a 37 μm cell strainer (Catalog #27250) for optimal results.

After preparation, resuspend cells at 5×10^7 cells/mL in recommended medium.

* SepMate™ IVD is only available in select regions where it is registered as an In Vitro Diagnostic (IVD) device for the isolation of mononuclear cells (MNCs) from whole blood or bone marrow by density gradient centrifugation. In all other regions SepMate™ is available for research use only (RUO).

LEUKAPHERESIS

Wash the peripheral blood leukapheresis sample by adding an equivalent volume of recommended medium or PBS containing 2% fetal bovine serum (FBS). Centrifuge at 500 x g for 10 minutes at room temperature (15 - 25°C). Remove the supernatant and resuspend the cells at 5 x 10^7 cells/mL in recommended medium.

Recommended Medium

EasySep™ Buffer (Catalog #20144), RoboSep™ Buffer (Catalog #20104), or PBS containing 2% FBS and 1 mM EDTA. Medium should be free of Ca++ and Mg++.



Directions for Use – Manual EasySep™ Protocols

See page 1 for Sample Preparation and Recommended Medium. Refer to Tables 1 and 2 for detailed instructions regarding the EasySep™ procedure for each magnet.

Table 1. EasySep™ Human B Cell Isolation Kit Protocol

2		EASYSEP™ MAGNETS				
STEP	INSTRUCTIONS	EasySep™ (Catalog #18000)	"The Big Easy" (Catalog #18001)			
	Prepare sample at the indicated cell concentration within the volume range.	5 x 10^7 cells/mL 0.25 - 2 mL	5 x 10^7 cells/mL 0.5 - 8 mL			
1	Add sample to required tube.	5 mL (12 x 75 mm) polystyrene round-bottom tube (e.g. Catalog #38007)	14 mL (17 x 95 mm) polystyrene round-bottom tube (e.g. Catalog #38008)			
2	Add Cocktail Enhancer to sample.	50 μL/mL of sample	50 μL/mL of sample			
3	Add Isolation Cocktail to sample.	50 μL/mL of sample	50 μL/mL of sample			
	Mix and incubate.	RT for 5 minutes	RT for 5 minutes			
4	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds	30 seconds			
5	Add RapidSpheres™ to sample and mix.	50 μL/mL of sample NO INCUBATION, proceed immediately to next step	50 μL/mL of sample NO INCUBATION, proceed immediately to next step			
6	Add recommended medium to top up the sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 2.5 mL	 Top up to 5 mL for samples ≤ 4 mL Top up to 10 mL for samples > 4 mL 			
	Place the tube (without lid) into the magnet and incubate.	RT for 3 minutes	RT for 3 minutes			
7	Pick up the magnet, and in one continuous motion invert the magnet and tube,* pouring the enriched cell suspension into a new tube.	Use a new 5 mL tube Use a new 14 mL tube				
8	Remove the tube from the magnet and place the new tube (without lid) into the magnet and incubate for a second separation.	RT for 1 minute RT for 1 minute				
9	Pick up the magnet, and in one continuous motion invert the magnet and tube,* pouring the enriched cell suspension into a new tube.	Isolated cells are ready for use Isolated cells are ready for use				

RT - room temperature (15 - 25°C)
* Leave the magnet and tube inverted for 2 - 3 seconds, then return upright. Do not shake or blot off any drops that may remain hanging from the mouth of the tube.



Table 2. EasySep™ Human B Cell Isolation Kit Protocol

	ysep numan b cen isolation kit Protocol			E,	ASYSEP™ MAGN	ETS	
STEP	INCTRUCTIONS	EasyEights™ (Catalog #18103)				Easy 50	
	INSTRUCTIONS	The state of the s	5 mL tube	5 mL tube 14 mL tube			Easy 50 (Catalog #18002)
1	Prepare sample at the indicated cell concentration within the volume range.		5 x 10^7 cells/mL 0.5 - 2 mL		5 x 10^7 cells/mL 1 - 8.5 mL		5 x 10^7 cells/mL 10 - 45 mL
	Add sample to required tube.	pol	5 mL (12 x 75 mm) ystyrene round-bottom tube (e.g. Catalog #38007)	ро	14 mL (17 x 95 mm lystyrene round-botton (e.g. Catalog #38008	n tube	50 mL (30 x 115 mm) conical tube (e.g. Catalog #38010)
2	Add Cocktail Enhancer to sample.		50 μL/mL of sample		50 μL/mL of sample		50 μL/mL of sample
Add Isolation Cocktail to sample.			50 μL/mL of sample		50 μL/mL of sample)	50 μL/mL of sample
3	Mix and incubate.		RT for 5 minutes		RT for 5 minutes		RT for 5 minutes
4	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds			30 seconds		30 seconds
_	Add RapidSpheres™ to sample.	50 μL/mL of sample			50 μL/mL of sample		50 μL/mL of sample
5	Mix and incubate.	RT for 3 minutes			RT for 3 minutes		RT for 3 minutes
6	Add recommended medium to top up sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.		Top up to 2.5 mL		 Top up to 5 mL for samples ≤ 4 mL Top up to 10 mL for samples > 4 mL 		 Top up to 25 mL for samples ≤ 20 mL Top up to 50 mL for samples > 20 mL
	Place the tube (without lid) into the magnet and incubate.	RT for 5 minutes			RT for 5 minutes		RT for 10 minutes
7	Carefully pipette** (do not pour) the enriched cell suspension into a new tube.	Use a new 5 mL tube			Use a new 14 mL tube		Use a new 50 mL tube
8	Remove the tube from the magnet and place the new tube (without lid) into the magnet and incubate for a second separation.		RT for 5 minutes		RT for 5 minutes		RT for 10 minutes
9	Carefully pipette** (do not pour) the enriched cell suspension into a new tube.	Isolated cells are ready for use		Iso	Isolated cells are ready for use		Isolated cells are ready for use

RT - room temperature (15 - 25°C)

** Collect the entire supernatant, all at once, into a single pipette (e.g. for EasyEightsTM 5 mL tube use a 2 mL serological pipette [Catalog #38002]; for EasyEightsTM 14 mL tube use a 10 mL serological pipette [Catalog #38004]).



Directions for Use - Fully Automated RoboSep™ Protocol

See page 1 for Sample Preparation and Recommended Medium. Refer to Table 3 for detailed instructions regarding the RoboSep™ procedure.

Table 3. RoboSep™ Human B Cell Isolation Kit Protocol

STEP	INSTRUCTIONS	RoboSep™ (Catalog #20000 and #21000)		
4	Prepare sample at the indicated cell concentration within the volume range.	5 x 10^7 cells/mL 0.5 - 8 mL		
1	Add sample to required tube.	14 mL (17 x 95 mm) polystyrene round-bottom tube (e.g. Catalog #38008)		
2	Select protocol.	Human B Cell Negative Isolation Kit - 17954		
3	Vortex RapidSpheres™. NOTE: Particles should appear evenly dispersed.	30 seconds		
4	Load the carousel.	Follow on-screen prompts		
4	Start the protocol.	Press the green "Run" button		
5	Unload the carousel when the run is complete.	Isolated cells are ready for use		


Notes and Tips


ASSESSING PURITY

For purity assessment of B cells (CD3-CD19+) by flow cytometry, use the following fluorochrome-conjugated antibody clones:

- Anti-Human CD19 Antibody, Clone HIB19 (Catalog #60005)
- Anti-Human CD3 Antibody, Clone UCHT1 (Catalog #60011; optional)
- Anti-Human CD45 Antibody, Clone HI30 (Catalog #60018; optional)

Data

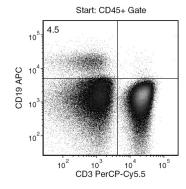

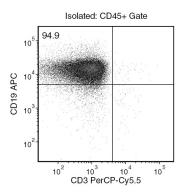

Starting with human PBMCs, the B cell content (CD3-CD19+) of the isolated fraction is typically 95.1 \pm 1.4% (mean \pm SD). In the above example, the purities of the start and final isolated fractions are 4.5% and 94.9%, respectively.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2020 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, EasyEights, EasyPlate, EasySep, RapidSpheres, RoboSep, and SepMate are trademarks of STEMCELL Technologies Canada Inc. Lymphoprep is a trademark of Alere Technologies. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.