

Small Molecules

Bleomycin

Antibiotic with antitumor properties



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #100-0550
100-0551

5 mg
10 mg

Product Description

Bleomycin is a glycopeptide-derived antibiotic isolated from the bacterium *Streptomyces verticillus* that creates DNA strand scission and exhibits antitumor activity against squamous cell carcinomas and malignant lymphomas (Galm et al.). Bleomycin induces lung injury leading to acute inflammatory response and pulmonary fibrosis (Kulkarni et al.). This product is supplied as the sulfate salt of the molecule and is a mixture of bleomycin A2 and bleomycin B2.

Alternative Names: Blenoxane; Bleomycin sulfate (salt)

CAS Number: 9041-93-4 (for A2), 9060-10-0 (for B2)


Chemical Formula: $C_{55}H_{84}N_{17}O_{21}S_3 \cdot XHSO_4$ (for A2), $C_{55}H_{84}N_{20}O_{21}S_2$ (for B2)

Molecular Weight: 1415.6 g/mol (for A2), 1425.5 g/mol (for B2)

Purity: $\geq 95\%$ (mixture of A2 and B2)

Chemical Name: 3-[[2-[2-[2-[[[(2S,3R)-2-[[[(2S,3S,4R)-4-[[[(2S,3R)-2-[[6-amino-2-[(1S)-3-amino-1-[[[(2S)-2,3-diamino-3-oxopropyl]amino]-3-oxopropyl]-5-methylpyrimidine-4-carbonyl]amino]-3-[(2R,3S,4S,5S,6S)-3-[(2R,3S,4S,5R,6R)-4-carbamoyloxy-3,5-dihydroxy-6-(hydroxymethyl)oxan-2-yl]oxy-4,5-dihydroxy-6-(hydroxymethyl)oxan-2-yl]oxy-3-(1H-imidazol-5-yl)propanoyl]amino]-3-hydroxy-2-methylpentanoyl]amino]-3-hydroxybutanoyl]amino]ethyl]-1,3-thiazol-4-yl]-1,3-thiazole-4-carbonyl]amino]propyl-dimethylsulfonium, sulfate

Structure:

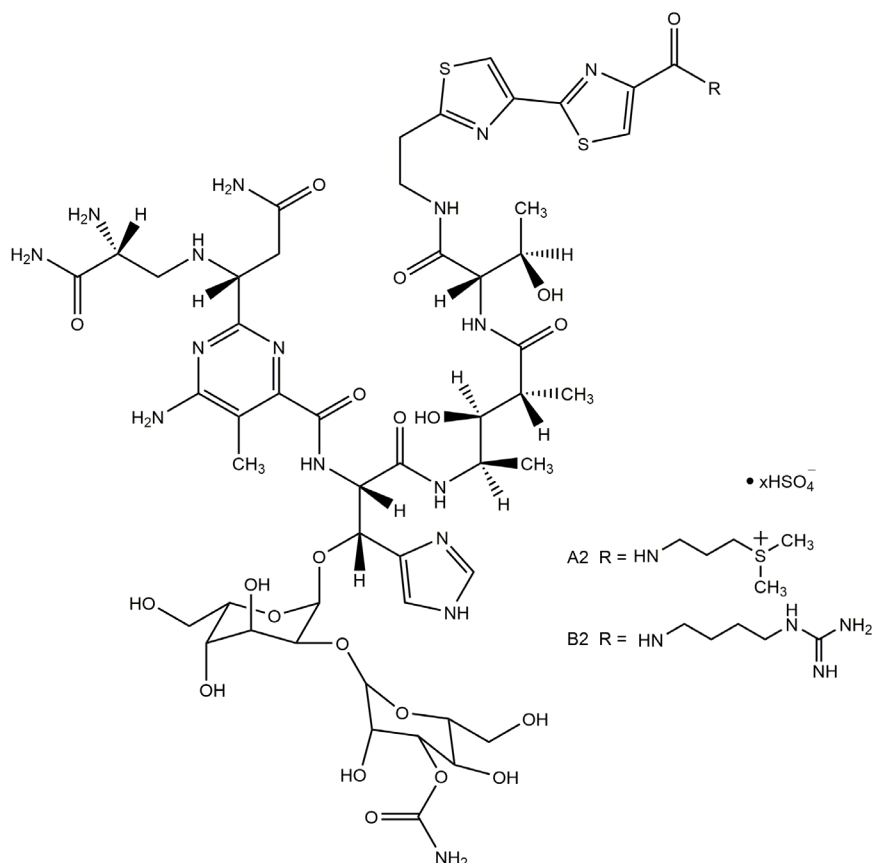

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect product from prolonged exposure to light. For long-term storage, store with a desiccant. Stable as supplied for 12 months from date of receipt.
Solubility:	<ul style="list-style-type: none">• PBS (pH 7.2) \leq 7.0 mM• DMSO \leq 9.1 mM <p>For example, to prepare a 5 mM stock solution in DMSO, resuspend 1 mg in 141 μL of DMSO.</p> <p>Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported; however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.</p> <p>For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.</p>

Published Applications

CANCER RESEARCH

- Induces single- and double-strand DNA breaks in tumor cells (Hecht).

DISEASE MODELING

- Induces lung fibrosis in mice (Kulkarni et al.).

References

Galm U et al. (2005) Antitumor antibiotics: bleomycin, enediynes, and mitomycin. Chem Rev 105(2): 739–58.

Hecht SM. (2000) Bleomycin: new perspectives on the mechanism of action. J Nat Prod 63(1): 158–68.

Kulkarni AA et al. (2013) The triterpenoid CDDO-Me inhibits bleomycin-induced lung inflammation and fibrosis. PLoS One 8(5): e63798.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules, or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2023 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.