

# Cytokines

## Human Recombinant IFN alpha 1

Interferon alpha 1, His tag



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #100-1336

100 µg

## Product Description

Interferon alpha 1 (IFNA1), also known as interferon alphaD (IFN alphaD), belongs to the type 1 interferon family of cytokines. Interferons display alpha-helical structures with four helices forming an antiparallel alpha-helix bundle (Walter). IFNA1 contains a His-residue tag at the carboxyl end of the polypeptide chain. Type 1 interferons bind to the interferon alpha receptor (IFNAR) composed of IFNAR1 and IFNAR2 subunits which are involved in interferon-induced JAK-STAT signaling (Shemesh et al.). Interferons have inhibitory effects on viral replication and pathogenesis, and studies have found that IFN alpha can prevent the spread of herpes simplex virus (HSV) (Mikloska & Cunningham) and human immunodeficiency virus (HIV) (George & Mattapallil). Pegylated forms of IFN alpha are currently approved for the treatment of chronic hepatitis B (Woo et al.). IFN alpha also displays anti-tumor effects by regulating cell growth and proliferation, and it is used in the treatment of different cancers (Tagliaferri et al.).

## Product Information

**Alternative Names:** IFL, IFN, IFNA, IFNA13, IFN alphaD

**Accession Number:** NP\_076918.1 (Cys24-Glu189) was expressed with a polyhistidine tag at the C-terminus.

**Amino Acid Sequence:** CDLPETHSLD NRRTLMLLAQ MSRISPSSCL MDRHDFGFPQ EEFDGNQFQK APAISVLHEL IQQIFNLFTT KDSSAAWDED LLDKFCTELY QQLNDLEACV MQEERVGETP LMNADSILAV KKYFRRLITLY LTEKKYSPCA WEVVR AEIMR SLSLSTNLQE RLRRKEHHHH HHHHHH

**Predicted Molecular Mass:** 20.8 kDa

**Species:** Human

**Formulation:** Lyophilized from sterile PBS, pH 7.4. Trehalose (5%), mannitol (5%), and 0.01% TWEEN® 80 are normally added as protectants before lyophilization.

**Source:** *Pichia pastoris*

## Specifications

**Activity:** The EC50 is ≤ 200 pg/mL as determined by a cell viability assay using WISH cells infected with vesicular stomatitis virus (VSV).

**Purity:** ≥ 92%

**Endotoxin Level:** Not available


## Preparation and Storage

**Storage:** Store at -20 to -80°C.

**Stability:** Stable as supplied for 12 months from date of receipt.

**Preparation:** Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.25 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. As a general guide, do not store at 2 - 8°C for more than 1 month or at -80°C for more than 3 months. Avoid repeated freeze-thaw cycles.

## Data

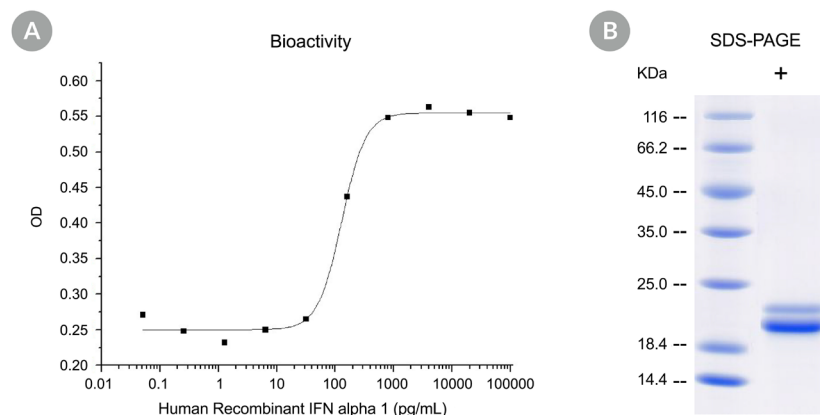

(A) The biological activity of Human Recombinant IFN alpha 1 was tested by its ability to promote the viability of WISH cells infected with VSV in a cytopathic effect (CPE) assay. Cell viability was measured using a fluorometric assay method. The EC50 is defined as the effective concentration of the cytokine at which cell survival is at 50% of maximum. The EC50 in the above example is  $\leq 200$  pg/mL.

(B) Human Recombinant IFN alpha 1 was resolved with SDS-PAGE under reducing (+) conditions and visualized by Coomassie Blue staining. Human Recombinant IFN alpha 1 has a predicted molecular mass of 20.8 kDa.

## Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, visit [www.stemcell.com/cytokines](http://www.stemcell.com/cytokines), or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

## References

- George J & Mattapallil JJ. (2018) Interferon- $\alpha$  subtypes as an adjunct therapeutic approach for human immunodeficiency virus functional cure. *Front Immunol* 9(2): 299.
- Mikloska Z & Cunningham AL. (2001) Alpha and gamma interferons inhibit herpes simplex virus type 1 infection and spread in epidermal cells after axonal transmission. *J Virol* 75(23): 11821–6.
- Shemesh M et al. (2021) IFNAR1 and IFNAR2 play distinct roles in initiating type I interferon-induced JAK-STAT signaling and activating STATs. *Sci Signal* 14(710): 4627.
- Tagliaferri P et al. (2005) New pharmacokinetic and pharmacodynamic tools for interferon-alpha (IFN- $\alpha$ ) treatment of human cancer. *Cancer Immunol Immunother* 54(1): 1–10.
- Walter MR. (2020) The role of structure in the biology of interferon signaling. *Front Immunol* 11. <https://doi.org/10.3389/fimmu.2020.606489>.
- Woo ASJ et al. (2017) Alpha-interferon treatment in hepatitis B. *Ann Transl Med* 5(7): 159.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2023 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. TWEEN is a registered trademark of Croda Americas LLC. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.