

Antibodies

Anti-Mouse CD25 Antibody, Clone 3C7

Rat monoclonal IgG2b antibody
against mouse CD25, unconjugated

Catalog #60010

500 µg 0.5 mg/mL



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

The 3C7 antibody reacts with murine CD25 (interleukin-2 receptor α chain; IL-2R α), an ~55 kDa type 1 transmembrane glycoprotein expressed on the surface of T and B cell progenitors and activated mature B and T cells. CD25 has been used as a marker to identify CD4+/FoxP3+ regulatory T cells in mice. Mutations in CD25 result in severe immunodeficiency, underscoring the importance of the roles the protein plays in the differentiation, activation and proliferation of lymphocytes, and in the maintenance of self-tolerance. CD25 per se has low affinity for its IL-2 ligand but forms dimers with CD122 (IL-2R β) and CD132 (IL-2R γ) that together associate to form the IL-2R receptor, which binds IL-2 with high affinity. CD25 acts to increase the specificity and affinity of IL-2 binding and is necessary for receptor clustering and signal transduction by the complex. Binding of the 3C7 antibody inhibits binding of IL-2 to free CD25 and the IL-2R receptor. The epitope is distinct from those recognized by the anti-CD25 antibody clones 7D4 and PC61.

Target Antigen Name:	CD25
Alternative Names:	IL-2Ralpha, Ly-43, p55, Tac
Gene ID:	16184
Species Reactivity:	Mouse
Host Species:	Rat (LEW)
Clonality:	Monoclonal
Clone:	3C7
Isotype:	IgG2b, kappa
Immunogen:	IL-2-dependent BALB/c mouse helper T cell clone HT-2
Conjugate:	Unconjugated

Applications


Verified:	FC
Reported:	CyTOF®, FA, FC, IHC
Special Applications:	This antibody clone has been verified for purity assessments of cells isolated with EasySep™ kits, including EasySep™ Mouse T Cell Enrichment Kit (Catalog #19751).

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FACS: Fluorescence activated cell sorting; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; RIA: Radioimmunoassay; WB: Western blotting

Properties

Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Purification:	The antibody was purified by affinity chromatography.
Stability and Storage:	Product stable at 2 - 8°C when stored undiluted. Do not freeze. For product expiry date, please contact techsupport@stemcell.com .
Directions for Use:	For flow cytometry the suggested use of this antibody is $\leq 1 \mu\text{g}$ per 1×10^6 cells in 100 μL volume. It is recommended that the antibody be titrated for optimal performance for each application.

Data

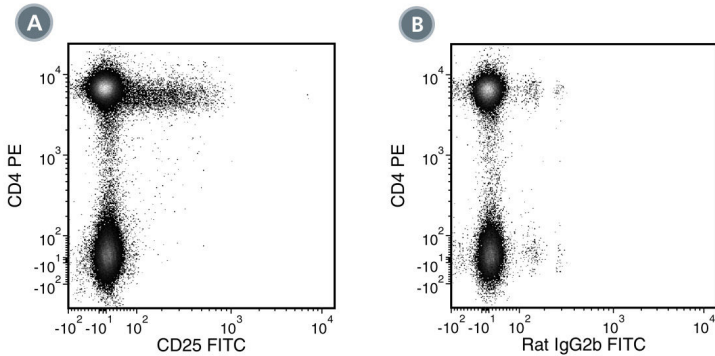

(A) Flow cytometry analysis of C57BL/6 mouse splenocytes (gated on CD3e+ cells) labeled with Anti-Mouse CD25 Antibody, Clone 3C7, followed by a mouse anti-rat IgG2b antibody, FITC and Anti-Mouse CD4, Clone RM4-5, PE (Catalog #60017PE).

(B) Flow cytometry analysis of C57BL/6 mouse splenocytes (gated on CD3e+ cells) labeled with Rat IgG2b, kappa Isotype Control Antibody, Clone RTK4530 (Catalog #60077), followed by a mouse anti-rat IgG2b antibody, FITC and Anti-Mouse CD4, Clone RM4-5, PE.

Related Products

For a complete list of antibodies, including other conjugates, sizes and clones, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/antibodies or contact us at techsupport@stemcell.com.

References

1. Read S et al. (2000) Cytotoxic T lymphocyte-associated antigen 4 plays an essential role in the function of CD25(+)CD4(+) regulatory cells that control intestinal inflammation. *J Exp Med* 192(2): 295–302.
2. Klug DB et al. (1998) Interdependence of cortical thymic epithelial cell differentiation and T-lineage commitment. *Proc Natl Acad Sci USA* 95(20): 11822–7. (IHC)
3. Lorenzo F et al. (1991) Structure-function study of the p55 subunit of murine IL-2 receptor by epitope mapping. *J Immunol* 147(9): 2970–7.
4. Moreau JL et al. (1987) Monoclonal antibodies identify three epitope clusters on the mouse p55 subunit of the interleukin 2 receptor: relationship to the interleukin 2-binding site. *Eur J Immunol* 17(7): 929–35. (Blocking, FA)
5. Malek TR et al. (1984) The murine IL 2 receptor. II. Monoclonal anti-IL 2 receptor antibodies as specific inhibitors of T cell function in vitro. *J Immunol* 133(4): 1976–82.
6. Ortega G et al. (1984) The murine IL 2 receptor. I. Monoclonal antibodies that define distinct functional epitopes on activated T cells and react with activated B cells. *J Immunol* 133(4): 1970–5. (Blocking, FA)

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists and EasySep are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.