

Small Molecules

Cyclopamine

Hedgehog pathway inhibitor; Inhibits Smoothened (SMO)

Catalog # 72072
72074

1 mg
5 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Cyclopamine is a steroid alkaloid that inhibits the Hedgehog pathway at the point of the pathway activator Smoothened. Cyclopamine binds to the heptahelical bundle of Smoothened, a G protein-coupled receptor, and prevents it from signaling further downstream (Chen et al.).

Molecular Name:	Cyclopamine
Alternative Names:	11-Deoxojervine; Jervine
CAS Number:	4449-51-8
Chemical Formula:	C ₂₇ H ₄₁ NO ₂
Molecular Weight:	411.6 g/mol
Purity:	≥ 95%
Chemical Name:	(2'R,3S,3'R,3'aS,6'S,6aS,6bS,7'aR,11aS,11bR)-1,2,3,3'a,4,4',5',6,6',6a,6b,7,7',7'a,8,11,11a,11b-octadecahydro-3',6',10,11b-tetramethyl-spiro[9H-benzo[a]fluorene-9,2'(3'H)-furo[3,2-b]pyridin]-3-ol

Structure:

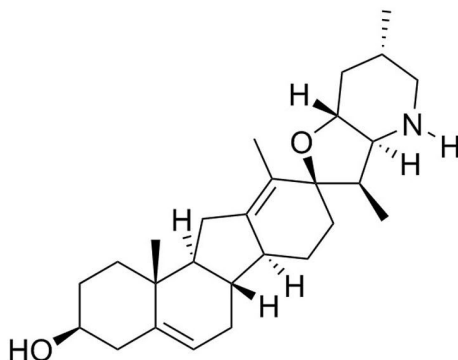

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt.
Solubility:	· Absolute ethanol ≤ 20 mM For example, to prepare a 10 mM stock solution in absolute ethanol, resuspend 1 mg in 243 µL of absolute ethanol. Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in absolute ethanol at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application. Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final ethanol concentration above 0.1% due to potential cell toxicity.

Published Applications

MAINTENANCE AND SELF-RENEWAL

- Reduces proliferation of rat neural progenitor cells and mouse neurospheres (Lai et al.; Palma & Ruiz i Altaba).
- Reduces proliferation of mouse mammospheres (Liu et al.).
- Inhibits the growth of human and mouse medulloblastoma cells, and human glioblastoma cells (Bar et al.; Berman et al.).

DIFFERENTIATION

- Promotes differentiation of pancreatic cells from human embryonic stem cells (D'Amour et al.).

CANCER RESEARCH

- Inhibits the growth of human and mouse medulloblastoma cells, and human glioblastoma cells (Bar et al.; Berman et al.).

References

- Bar EE et al. (2007) Cyclopamine-mediated hedgehog pathway inhibition depletes stem-like cancer cells in glioblastoma. *Stem Cells* 25(10): 2524–33.
- Berman DM et al. (2002) Medulloblastoma growth inhibition by hedgehog pathway blockade. *Science* 297(5586): 1559–61.
- Chen JK et al. (2002) Inhibition of Hedgehog signaling by direct binding of cyclopamine to Smoothened. *Genes & Development* 16(21): 2743–8.
- D'Amour KA et al. (2006) Production of pancreatic hormone-expressing endocrine cells from human embryonic stem cells. *Nature Biotechnology* 24(11): 1392–401.
- Lai K et al. (2003) Sonic hedgehog regulates adult neural progenitor proliferation in vitro and in vivo. *Nature Neuroscience* 6(1): 21–7.
- Liu S et al. (2006) Hedgehog signaling and Bmi-1 regulate self-renewal of normal and malignant human mammary stem cells. *Cancer Research* 66(12): 6063–71.
- Palma V & Ruiz i Altaba A. (2004) Hedgehog-GLI signaling regulates the behavior of cells with stem cell properties in the developing neocortex. *Development* 131(2): 337–45.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2018 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.