

## Small Molecules

### L-Ascorbic Acid

Antioxidant; Reducing agent

Catalog # 72132

500 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

L-Ascorbic Acid is a naturally occurring lactone that is produced by plants and many animals, but not humans or other primates. It acts as an electron donor (i.e. reducing agent), and shows antioxidant activity, particularly against reactive oxygen species. L-Ascorbic Acid is a cofactor for monooxygenase and dioxygenase as well as other enzymes (Arrighoni and De Tullio; Du et al.).

Molecular Name:	L-Ascorbic Acid
Alternative Names:	Vitamin C; (+)-Ascorbic Acid; NSC 218455; NSC 33832
CAS Number:	50-81-7
Chemical Formula:	C <sub>6</sub> H <sub>8</sub> O <sub>6</sub>
Molecular Weight:	176.1 g/mol
Purity:	≥ 95%
Chemical Name:	(R)-5-((S)-1,2-dihydroxyethyl)-3,4-dihydroxyfuran-2(5H)-one
Structure:	

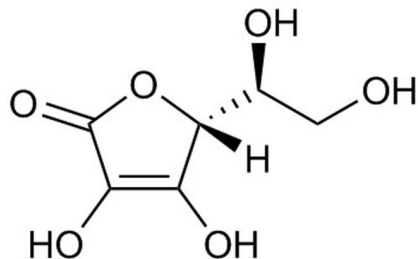

## Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com.
Solubility:	· PBS (pH 7.2) ≤ 18 mM · DMSO ≤ 55 mM For example, to prepare a 10 mM stock solution in PBS, resuspend 100 mg in 56.8 mL of PBS (pH 7.2).

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### MAINTENANCE AND SELF-RENEWAL

- Supports proliferation of mesenchymal stem cells (Choi et al.).

### REPROGRAMMING

- Increases the efficiency of reprogramming mouse and human fibroblasts to induced pluripotent stem (iPS) cells (Esteban et al.) partly through JHDM1 histone demethylase activity (Wang et al.).
- Prevents aberrant DNA methylation of the Dlk1-Dio3 locus during reprogramming of mouse somatic cells to iPS cells (Stadtfeld et al.).

### DIFFERENTIATION

- Promotes differentiation of osteoblasts from human and mouse mesenchymal cells (Pittenger et al., Tropel et al.).
- Promotes differentiation of osteoblasts from mouse embryonic stem (ES) cells (zur Nieden et al.).
- Enhances differentiation of cardiomyocytes from mouse ES cells (Takahashi et al.).

## References

- Arrigoni O & De Tullio MC. (2002) Ascorbic acid: much more than just an antioxidant. *Biochim Biophys Acta* 1569(1-3): 1–9.
- Choi K-M et al. (2008) Effect of ascorbic acid on bone marrow-derived mesenchymal stem cell proliferation and differentiation. *J Biosci Bioeng* 105(6): 586–94.
- Du J et al. (2012) Ascorbic acid: chemistry, biology and the treatment of cancer. *Biochim Biophys Acta* 1826(2): 443–57.
- Esteban MA et al. (2010) Vitamin C enhances the generation of mouse and human induced pluripotent stem cells. *Cell Stem Cell* 6(1): 71–9.
- Pittenger MF et al. (1999) Multilineage potential of adult human mesenchymal stem cells. *Science* 284(5411): 143–7.
- Stadtfeld M et al. (2012) Ascorbic acid prevents loss of Dlk1-Dio3 imprinting and facilitates generation of all-iPS cell mice from terminally differentiated B cells. *Nat Genet* 44(4): 398–405, S1–2.
- Takahashi T et al. (2003) Ascorbic acid enhances differentiation of embryonic stem cells into cardiac myocytes. *Circulation* 107(14): 1912–6.
- Tropel P et al. (2004) Isolation and characterisation of mesenchymal stem cells from adult mouse bone marrow. *Exp Cell Res* 295(2): 395–406.
- Wang T et al. (2011) The histone demethylases Jhdm1a/1b enhance somatic cell reprogramming in a vitamin-C-dependent manner. *Cell Stem Cell* 9(6): 575–87.
- zur Nieden NI et al. (2003) In vitro differentiation of embryonic stem cells into mineralized osteoblasts. *Differentiation* 71(1): 18–27.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2015 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.