

G418

Antibiotic for selecting transfected mammalian cells

Catalog #03812 250 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE


Product Description

Antibiotic G418 is an aminoglycoside related to gentamicin and is toxic to bacteria, yeast, higher plants, mammalian cells, protozoans, and helminths. G418 binds to 80S ribosomes in prokaryotic and eukaryotic cells and inhibits protein synthesis, resulting in cell death. Resistance to the antibiotic is conferred by aminoglycoside-3-phosphotransferase. Vectors containing the G418 resistance gene (neo) are commonly used for transfecting eukaryotic cells, with selection of transfected cells achieved by culturing them in the presence of G418. Cells may exhibit variable resistance to the antibiotic and can take up to a week to die (adherent cells may be more sensitive).

The G418 antibiotic is suitable for use in cell culture for the selection and maintenance of G418-resistant clones when using ClonaCell™-TCS Medium (Catalog #03814), ClonaCell™-CHO CD Medium (Catalog #03815), ClonaCell™-CHO ACF Medium (Catalog #03816), ClonaCell™-CHO CD Liquid Medium (Catalog #03817), or other cell culture media.

Molecule Name:	G418
Alternative Names:	Geneticin; (2R,3S,4R,5R,6S)-5-Amino-6-[(1R,2S,3S,4R,6S)-4,6-diamino-3-[(2R,3R,4R,5R)-3,5-dihydroxy-5-methyl-4-methylaminoxan-2-yl]oxy-2-hydroxycyclohexyl]oxy- 2-(1-hydroxyethyl)oxane-3,4-diol
CAS Number:	108321-42-2
Chemical Formula:	$C_{20}H_{40}N_4O_{10} \bullet 2H_2SO_4$
Molecular Weight:	692.71 g/mol
Chemical Name:	O-2-Amino-2,7-dideoxy-D-glycero- α -D-gluco-heptopyranosyl-(1 \rightarrow 4)-O-(3-deoxy-4-C-methyl-3-(methylamino)- β -L-arabinopyranosyl- (1 \rightarrow 6))-D-streptamin

Structure:

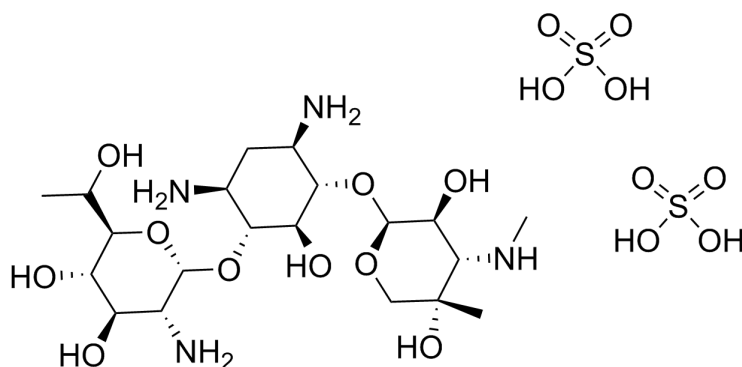

Physical Appearance:	White to off-white powder.
Storage:	Store at 2 - 8°C.
Solubility:	Soluble in water and aqueous buffers.

Please refer to the Safety Data Sheet (SDS) for hazard information.

Handling / Directions for Use

NOTE: Protect product from light. Avoid moisture uptake as this can decrease the potency of the product.

1. Place G418 at room temperature (15 - 25°C) for 30 minutes before use.
2. Dissolve G418 in 1 mL of distilled water.

NOTE: If not used immediately, aliquot and store at -20°C. After thawing the aliquots, do not re-freeze. Sterilize solutions by filtration, not by autoclaving.

NOTE: The effective concentration of G418 needed to kill cells varies by cell type. Typical concentrations used for mammalian cell selection are 0.2 - 2 mg/mL; the optimal concentration must be determined for each cell type.

For more information, refer to the Technical Manual: Semi-Solid Cloning Testing Guidelines - ClonaCell™-CHO (Document #28096), available at www.stemcell.com or contact us to request a copy.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, and ClonaCell are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.