

Antibodies

Anti-Mouse CD11c Antibody, Clone N418, PerCP-Cy5.5



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #60002PS
#60002PS.1

Hamster (Armenian) monoclonal IgG
antibody against mouse CD11c, PerCP-
Cy5.5-conjugated

100 µg 0.2 mg/mL
25 µg 0.2 mg/mL

Product Description

The N418 antibody reacts with CD11c (α X integrin), a 150 kDa type 1 transmembrane glycoprotein that associates non-covalently with CD18 (β 2 integrin) to form a heterodimeric cell surface adhesion receptor. Through its interaction with ligands such as iC3b, fibrinogen, and CD54, the CD11c/CD18 receptor is involved in several immune response processes, including cell migration, stimulation of cytokine production by monocytes and macrophages, T cell proliferation, leukocyte recruitment, and phagocytosis. In mice, CD11c is expressed on dendritic cells, macrophages, monocytes, granulocytes, NK cells, and a subset of T cells.

| | |
|----------------------|--|
| Target Antigen Name: | CD11c |
| Alternative Names: | alphaX integrin, CR4, integrin alphaX chain, p150 |
| Gene ID: | 16411 |
| Species Reactivity: | Mouse |
| Host Species: | Hamster (Armenian) |
| Clonality: | Monoclonal |
| Clone: | N418 |
| Isotype: | IgG |
| Immunogen: | Mouse spleen dendritic cells |
| Conjugate: | PerCP-Cy5.5 (Peridinin chlorophyll protein complex-Cyanine5.5) |

Applications


| | |
|-----------------------|---|
| Verified: | FC |
| Reported: | FC |
| Special Applications: | This antibody clone has been verified for purity assessments of cells isolated with EasySep™ kits, including EasySep™ Mouse CD11c Positive Selection Kit II (Catalog #18780). |

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FACS: Fluorescence-activated cell sorting; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; RIA: Radioimmunoassay; WB: Western blotting

Properties

| | |
|------------------------|--|
| Formulation: | Phosphate-buffered saline, pH 7.2, containing 0.09% sodium azide and 0.1% gelatin |
| Purification: | The antibody was purified by affinity chromatography and conjugated with PerCP-Cy5.5 under optimal conditions. The solution is free of unconjugated PerCP-Cy5.5. |
| Stability and Storage: | Product stable at 2 - 8°C when stored undiluted. Do not freeze. Protect product from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com . |
| Directions for Use: | For flow cytometry, the suggested use of this antibody is ≤ 0.25 µg per 1×10^6 cells in 100 µL. It is recommended that the antibody be titrated for optimal performance for each application. |

Data

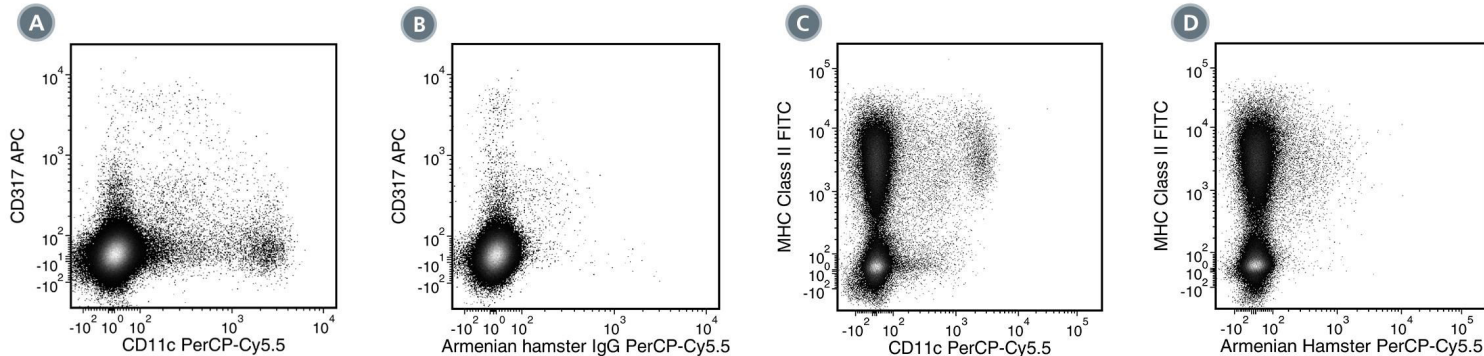

(A) Flow cytometry analysis of C57BL/6 mouse splenocytes labeled with Anti-Mouse CD11c Antibody, Clone N418, PerCP-Cy5.5 and an anti-mouse CD317 antibody, APC.

(B) Flow cytometry analysis of C57BL/6 mouse splenocytes labeled with an Armenian hamster IgG PerCP-Cy5.5 isotype control antibody and an anti-mouse CD317 antibody, APC.

(C) Flow cytometry analysis of C57BL/6 mouse splenocytes labeled with Anti-Mouse CD11c Antibody, Clone N418, PerCP-Cy5.5 and an anti-mouse MHC class II antibody, FITC.

(D) Flow cytometry analysis of C57BL/6 mouse splenocytes labeled with an Armenian hamster IgG PerCP-Cy5.5 isotype control antibody and an anti-mouse MHC class II antibody, FITC.

Related Products

For a complete list of antibodies, including other conjugates, sizes and clones, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/antibodies or contact us at techsupport@stemcell.com.

References

- Lai JD et al. (2018) Early cellular interactions and immune transcriptome profiles in human factor VIII-exposed hemophilia A mice. *J Thromb Haemost.* Epub ahead of print, DOI: 10.1111/jth.13936. (FC)
- Oderup C et al. (2013) Canonical and noncanonical Wnt proteins program dendritic cell responses for tolerance. *J Immunol* 190(12): 6126–34. (FC)
- Schneider D et al. (2012) Neonatal rhinovirus infection induces mucous metaplasia and airways hyperresponsiveness. *J Immunol* 188(6): 2894–904. (FC)
- Grewal JS et al. (2011) Salivary glands act as mucosal inductive sites via the formation of ectopic germinal centers after site-restricted MCMV infection. *FASEB J* 25(5): 1680–96. (IF, IHC)
- Bankoti J et al. (2010) Effects of TCDD on the fate of naive dendritic cells. *Toxicol Sci* 115(2): 422–34. (FC)
- Roland CL et al. (2009) Inhibition of vascular endothelial growth factor reduces angiogenesis and modulates immune cell infiltration of orthotopic breast cancer xenografts. *Mol Cancer Ther* 8(7): 1761–71. (FC, IHC)
- You Y et al. (2009) Cutting edge: Primary and secondary effects of CD19 deficiency on cells of the marginal zone. *J Immunol* 182(12): 7343–7. (IF, IHC)
- Cervantes-Barragan L et al. (2007) Control of coronavirus infection through plasmacytoid dendritic-cell-derived type I interferon. *Blood* 109(3): 1131–7. (FC)
- Turnquist HR et al. (2007) Rapamycin-conditioned dendritic cells are poor stimulators of allogeneic CD4+ T cells, but enrich for antigen-specific Foxp3+ T regulatory cells and promote organ transplant tolerance. *J Immunol* 178(11): 7018–31. (FC)
- Chin RK et al. (2006) Lymphotoxin pathway-directed, autoimmune regulator-independent central tolerance to arthritogenic collagen. *J Immunol* 177(1): 290–7. (IF, IHC)
- Kishimoto T et al. (Eds.). (1998) *Leucocyte Typing VI: White Cell Differentiation Antigens*. New York: Garland Publishing Inc.
- Barclay AN et al. (Eds.). (1997) *The Leucocyte Antigen Factsbook, Second Edition* (pp. 149–51). New York: Academic Press.
- Metlay JP et al. (1990) The distinct leukocyte integrins of mouse spleen dendritic cells as identified with new hamster monoclonal antibodies. *J Exp Med* 171(5): 1753–71. (FA/Blocking, FC, IHC, IP)

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2018 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, and EasySep are trademarks of STEMCELL Technologies Canada Inc. CyTOF is a registered trademark of Fluidigm Corporation. All other trademarks are the property of their respective holders. Alexa Fluor and Pacific Blue are trademarks of Life Technologies Corporation. Antibodies conjugated to Alexa Fluor® or Pacific Blue™ are licensed for internal research use only and sale is expressly conditioned on the buyer not using the antibody for manufacturing, performing a service or medical test, or otherwise generating revenue. For use other than research, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.