

## Anti-Mouse CD4 Antibody, Clone RM4-5, Biotin



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Antibodies

Rat monoclonal IgG2a antibody  
against mouse CD4, biotin-conjugated

Catalog #60017BT  
#60017BT.1

500 µg 0.5 mg/mL  
50 µg 0.5 mg/mL

FOR RESEARCH USE ONLY. NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES.

## Product Description

The RM4-5 antibody reacts with mouse CD4, an ~55 kDa single-chain type 1 transmembrane glycoprotein and member of the immunoglobulin (Ig) superfamily; CD4 contains four extracellular Ig-like domains. CD4 is expressed at relatively high levels by most thymocytes and a subpopulation of T cells (T-helper/inducer cells), and at lower levels on dendritic cells. In the mouse, CD4 is not expressed by monocytes/macrophages. CD4 binds to a non-polymorphic region of MHC II and acts as a co-receptor to the T cell receptor (TCR) in MHC II-restricted antigen recognition by enhancing the avidity of the association between the TCR and MHC II-antigen complex. CD4 also functions to amplify signals from the TCR to the cytoplasm through the interaction of its intracellular domain with cytoplasmic tyrosine kinases such as Lck. Binding of the RM4-5 antibody to CD4 inhibits ligand binding in vitro. Moreover, binding of the RM4-5 antibody can be blocked by the clone GK1.5 antibody.

Target Antigen Name:	CD4
Alternative Names:	L3T4, T4
Gene ID:	12504
Species Reactivity:	Mouse
Host Species:	Rat (DA)
Clonality:	Monoclonal
Clone:	RM4-5
Isotype:	IgG2a, kappa
Immunogen:	BALB/c mouse thymocytes
Conjugate:	Biotin

## Applications


Verified:	CellSep, FC
Reported:	FC, IHC
Special Applications:	This antibody clone has been verified for purity assessments of cells isolated with EasySep™ kits, including EasySep™ Mouse CD4+ T Cell Isolation Kit (Catalog #19852).

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; WB: Western blotting

## Properties

Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Purification:	The antibody was purified by affinity chromatography and conjugated with biotin under optimal conditions. The solution is free of unconjugated biotin.
Stability and Storage:	Product stable at 2 - 8°C when stored undiluted. Do not freeze. For product expiry date, please contact <a href="mailto:techsupport@stemcell.com">techsupport@stemcell.com</a> .
Directions for Use:	For flow cytometry the suggested use of this antibody is ≤ 0.25 µg per 1 × 10 <sup>6</sup> cells in 100 µL volume. It is recommended that the antibody be titrated for optimal performance for each application.

## Data

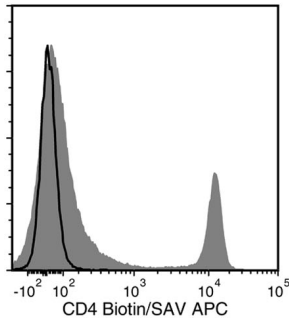

Flow cytometry analysis of C57BL/6 mouse splenocytes labeled with Anti-Mouse CD4 Antibody, Clone RM4-5, Biotin followed by streptavidin (SAV) APC (filled histogram) or a rat IgG2a, kappa biotin isotype control antibody followed by SAV APC (solid line histogram).

## Related Products

For a complete list of antibodies, including other conjugates, sizes and clones, as well as related products available from STEMCELL Technologies, please visit our website at [www.stemcell.com/antibodies](http://www.stemcell.com/antibodies) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

## References

1. Kruisbeek AM. In vivo depletion of CD4- and CD8-specific T cells. *Current Protocols in Immunology*. New York: John Wiley and Sons: 4.1.1-4.1.5, 1991 (Depletion, FA)
2. Nitta H, et al. Improved in situ immunodetection of leukocytes on paraffin embedded mouse spleen. *Cell Vision* 4(1): 73-80, 1997 (IHC)
3. Bosselut R, et al. Association of the adaptor molecule LAT with CD4 and CD8 coreceptors identifies a new coreceptor function in T cell receptor signal transduction. *J Exp Med* 190(10): 1517-26, 1999 (IP)
4. Muraille E, et al. Amastigote load and cell surface phenotype of infected cells from lesions and lymph nodes of susceptible and resistant mice infected with *Leishmania major*. *Infect Immun* 71(5): 2704-15, 2003 (IHC)
5. Irie J, et al. NOD.c3c4 congenic mice develop autoimmune biliary disease that serologically and pathogenetically models human primary biliary cirrhosis. *J Exp Med* 203(5): 1209-19, 2006 (IHC)
6. Leon-Ponte M, et al. Serotonin provides an accessory signal to enhance T-cell activation by signaling through the 5-HT7 receptor. *Blood* 109(8): 3139-46, 2007 (FC)
7. Bourdeau A, et al. TC-PTP-deficient bone marrow stromal cells fail to support normal B lymphopoiesis due to abnormal secretion of interferon- $\gamma$ . *Blood* 109(10): 4220-28, 2007 (FC)
8. Matsumoto M, et al. CD43 collaborates with P-selectin glycoprotein ligand-1 to mediate E-selectin-dependent T cell migration into inflamed skin. *J Immunol* 178(4): 2499-506, 2007
9. Shigeta A, et al. An L-selectin ligand distinct from P-selectin glycoprotein ligand-1 is expressed on endothelial cells and promotes neutrophil rolling in inflammation. *Blood* 112(13): 4915-23, 2008
10. Rodrigues-Manzanet R, et al. T and B cell hyperactivity and autoimmunity associated with niche-specific defects in apoptotic body clearance in TIM-4-deficient mice. *Proc Natl Acad Sci USA* 107(19): 8706-11, 2010

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485 MEDICAL DEVICE STANDARDS.

Copyright © 2014 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists and EasySep are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. Alexa Fluor® is a registered trademark of Life Technologies Corporation. This product is licensed for internal research use only and its sale is expressly conditioned on the buyer not using it for manufacturing, performing a service, or medical test, or otherwise generating revenue. For use other than research, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or [outlicensing@lifetech.com](mailto:outlicensing@lifetech.com). While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.