

## Anti-Mouse F4/80 Antibody, Clone BM8, PerCP-Cy5.5



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Antibodies

Rat monoclonal IgG2a antibody  
against human, mouse F4/80, PerCP-  
Cy5.5-conjugated

Catalog #60027PS	100 µg	0.2 mg/mL
#60027PS.1	25 µg	0.2 mg/mL

FOR RESEARCH USE ONLY. NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES.

## Product Description

The BM8 antibody reacts with the F4/80 antigen (also termed Ly-71 in mouse and EMR-1 in human), an ~160 kDa transmembrane glycoprotein belonging to the EGF-TM7 family of G-protein-coupled receptors. F4/80 is considered a marker of choice for the identification of mature tissue macrophages, being broadly but variably expressed by this cell type in the liver (Kupffer cells), skin (Langerhans cells), bone marrow stroma, pancreas, thymus, spleen (red pulp), lung, and other tissues. It is also expressed by circulating monocytes, eosinophils and a subset of dendritic cells. F4/80 expression levels increase following activation of macrophages. F4/80 is reportedly the only macrophage marker suitable for distinguishing destructive from non-destructive inflammatory processes in the pancreas. The protein is thought to play a role in the generation of CD8+ regulatory T cells involved in immune tolerance.

Target Antigen Name:	F4/80
Alternative Names:	EMR1, Ly71
Gene ID:	13733
Species Reactivity:	Human, Mouse
Host Species:	Rat
Clonality:	Monoclonal
Clone:	BM8
Isotype:	IgG2a, kappa
Immunogen:	Murine (BALB/c) macrophages obtained from 14-day-old bone marrow cell cultures
Conjugate:	PerCP-Cy5.5

## Applications


Verified:	FC
Reported:	FC
Special Applications:	This antibody clone has been verified for purity assessments of cells isolated with EasySep™ kits, including EasySep™ Mouse Monocyte Enrichment Kit (Catalog #19761) and EasySep™ Mouse CD11b Positive Selection Kit (Catalog #18770).

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; WB: Western blotting

## Properties

Formulation:	Phosphate-buffered saline, pH 7.2, containing 0.09% sodium azide and 0.1% gelatin
Purification:	The antibody was purified by affinity chromatography and conjugated with PerCP-Cy5.5 under optimal conditions. The solution is free of unconjugated PerCP-Cy5.5.
Stability and Storage:	Product stable at 2 - 8°C when stored undiluted. Do not freeze. Protect product from prolonged exposure to light. For product expiry date, please contact <a href="mailto:techsupport@stemcell.com">techsupport@stemcell.com</a> .
Directions for Use:	For flow cytometry the suggested use of this antibody is ≤ 1 µg per 1 × 10 <sup>6</sup> cells in 100 µL volume. It is recommended that the antibody be titrated for optimal performance for each application.

## Data

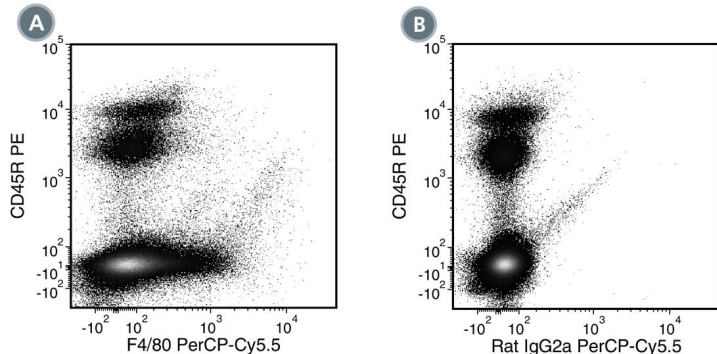

(A) Flow cytometry analysis of C57BL/6 mouse bone marrow cells labeled with Anti-Mouse F4/80 Antibody, Clone BM8, PerCP-Cy5.5 and Anti-Mouse CD45R Antibody, Clone RA3-6B2, PE (Catalog #60019PE).

(B) Flow cytometry analysis of C57BL/6 mouse bone marrow cells labeled with a rat IgG2a, kappa isotype control antibody, PerCP-Cy5.5 and Anti-Mouse CD45R Antibody, Clone RA3-6B2, PE.

## Related Products

For a complete list of antibodies, including other conjugates, sizes and clones, as well as related products available from STEMCELL Technologies, please visit our website at [www.stemcell.com/antibodies](http://www.stemcell.com/antibodies) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

## References

1. Leenen PJ, et al. Markers of mouse macrophage development detected by monoclonal antibodies. *J Immunol Methods* 174(1-2): 5-19, 1994 (FC)
2. Haidl ID, et al. The macrophage cell surface glycoprotein F4/80 is a highly glycosylated proteoglycan. *Eur J Immunol* 26(5): 1139-46, 1996
3. Murayama T, et al. Intraperitoneal administration of anti-c-fms monoclonal antibody prevents initial events of atherosclerosis but does not reduce the size of advanced lesions in apolipoprotein E-deficient mice. *Circulation* 99(13): 1740-46, 1999 (IHC)
4. Mackler AM, et al. Distribution and activation of uterine mononuclear phagocytes in peripartum endometrium and myometrium of the mouse. *Biol Reprod* 62(5): 1193-1200, 2000 (IHC)
5. Schaller E, et al. Inactivation of the F4/80 glycoprotein in the mouse germ line. *Mol Cell Biol* 22(22): 8035-43, 2002 (ICC, IHC)
6. Lin H-H, et al. The macrophage F4/80 receptor is required for the induction of antigen-specific efferent regulatory T cells in peripheral tolerance. *J Exp Med* 201(10): 1615-25, 2005 (IHC)
7. Poeckel D, et al. Dual 12/15- and 5-lipoxygenase deficiency in macrophages alters arachidonic acid metabolism and attenuates peritonitis and atherosclerosis in ApoE knock-out mice. *J Biol Chem* 284(42): 21077-89, 2009 (FC)
8. Koenen TB, et al. The inflammasome and caspase-1 activation: a new mechanism underlying increased inflammatory activity in human visceral adipose tissue. *Endocrinology* 152(10): 3769-78, 2011 (FC)
9. Meyer Zu Hörste M, et al. A novel mechanism involved in the pathogenesis of Graves ophthalmopathy (GO): clathrin is a possible targeting molecule for inhibiting local immune response in the orbit. *J Clin Endocrinol Metab* 96(11): E1727-36, 2011 (FC)
10. Balamurugan K, et al. FBXW7 $\alpha$  attenuates inflammatory signalling by downregulating C/EBP $\delta$  and its target gene Tlr4. *Nat Commun* 4: 1662, 2013 (FC)
11. Majumder M, et al. A practical and sensitive method of quantitating lymphangiogenesis in vivo. *Lab Invest* 93(7): 779-91, 2013 (IF, IHC)

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485 MEDICAL DEVICE STANDARDS.

Copyright © 2014 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists and EasySep are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. Alexa Fluor® is a registered trademark of Life Technologies Corporation. This product is licensed for internal research use only and its sale is expressly conditioned on the buyer not using it for manufacturing, performing a service, or medical test, or otherwise generating revenue. For use other than research, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or [outlicensing@lifetech.com](mailto:outlicensing@lifetech.com). While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.