

Small Molecules

(+)-Bay K8644

L-type calcium channel inhibitor

Catalog # 72362
72364

1 mg
10 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Bay K8644, originally described as a modulator of potential operated calcium channels, exists as two enantiomers that have opposite actions (Hess et al.; Nowycky et al.; O'Neill et al.; Ravens et al.; Yamamoto et al.). (+)-Bay K8644 is an L-type channel blocker that has negative inotropic and vasodilatory effects (Artigas et al.; Franckowiak et al.; Ravens et al.).

| | |
|--------------------|---|
| Molecular Name: | (+)-Bay K8644 |
| Alternative Names: | (+)-Bay-K 8644; (+)-Bay-R 4407; (R)-(+)-Bay K 8644; NI 105; R 4407 |
| CAS Number: | 98791-67-4 |
| Chemical Formula: | C ₁₆ H ₁₅ F ₃ N ₂ O ₄ |
| Molecular Weight: | 356.3 g/mol |
| Purity: | ≥ 98% |
| Chemical Name: | R-(+)-1,4-Dihydro-2,6-dimethyl-5-nitro-4-[2-(trifluoromethyl)phenyl]-3-pyridinecarboxylic acid methyl ester |
| Structure: | |

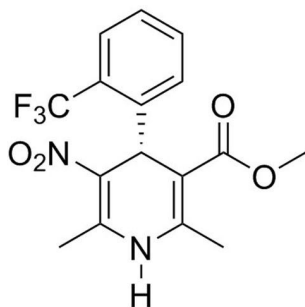

Properties

| | |
|----------------------|---|
| Physical Appearance: | A crystalline solid |
| Storage: | Product stable at -20°C as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt. |
| Solubility: | · DMSO ≤ 140 mM · Absolute ethanol ≤ 140 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 281 µL of fresh DMSO. |

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

REPROGRAMMING

· When combined with BIX01294 (Catalog #72042), enables reprogramming of mouse embryonic fibroblasts after transduction with OCT4 and KLF4 only (Shi et al.).

DIFFERENTIATION

· Increases neuronal differentiation from neural stem and progenitor cells (NSCs) derived from the brains of postnatal mice (D'Ascenzo et al.).

References

- Artigas P et al. (2003) Effects of the enantiomers of BayK 8644 on the charge movement of L-type Ca channels in guinea-pig ventricular myocytes. *J Membr Biol* 193(3): 215–27.
- D'Ascenzo M et al. (2006) Role of L-type Ca²⁺ channels in neural stem/progenitor cell differentiation. *Eur J Neurosci* 23(4): 935–44.
- Franckowiak G et al. (1985) The optical isomers of the 1,4-dihydropyridine BAY K 8644 show opposite effects on Ca channels. *Eur J Pharmacol* 114(2): 223–6.
- Hess P et al. (1984) Different modes of Ca channel gating behaviour favoured by dihydropyridine Ca agonists and antagonists. *Nature* 311(5986): 538–44.
- Nowycky MC et al. (1985) Long-opening mode of gating of neuronal calcium channels and its promotion by the dihydropyridine calcium agonist Bay K 8644. *Proc Natl Acad Sci USA* 82(7): 2178–82.
- O'Neill SK & Bolger GT. (1988) Enantiomer selectivity and the development of tolerance to the behavioral effects of the calcium channel activator BAY K 8644. *Brain Res Bull* 21(6): 865–72.
- Ravens U & Schöpper H-P. (1990) Opposite cardiac actions of the enantiomers of Bay K 8644 at different membrane potentials in guinea-pig papillary muscles. *Naunyn Schmiedeberg's Arch Pharmacol* 341(3).
- Shi Y et al. (2008) Induction of pluripotent stem cells from mouse embryonic fibroblasts by Oct4 and Klf4 with small-molecule compounds. *Cell Stem Cell* 3(5): 568–74.
- Yamamoto H et al. (1984) Bay K8644 differentiates between potential and receptor operated Ca²⁺ channels. *Eur J Pharmacol* 102(3-4): 555–7.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.