

## Small Molecules

SP600125

JNK pathway inhibitor; Inhibits JNK1, JNK2, and JNK3

Catalog # 72642

5 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

SP600125 is an inhibitor of c-Jun N-terminal kinase (JNK). The three isoforms of JNK are members of the MAP kinase superfamily that induce the expression of immediate-early genes in response to specific stress and inflammatory signals. Through these actions, the JNK enzymes modulate cell proliferation, apoptosis, differentiation, and autophagy. SP600125 is a potent and reversible inhibitor of JNK1-3 ( $IC_{50} = 0.11 \mu\text{M}$ ; Bennett et al.). It is cell-permeable and dose-dependently inhibits c-Jun phosphorylation in cells, blocking the expression of COX-2 and TNF- $\alpha$  in monocytes and IL-10, TNF- $\alpha$ , and IFN- $\gamma$  in T-cells (Bennett et al.).

Molecular Name:	SP600125
Alternative Names:	NSC 75890; Pyrazolanthrone; 1PMV
CAS Number:	129-56-6
Chemical Formula:	$C_{14}H_8N_2O$
Molecular Weight:	220.2 g/mol
Purity:	$\geq 98\%$
Chemical Name:	anthra[1,9-cd]pyrazol-6(2H)-one
Structure:	

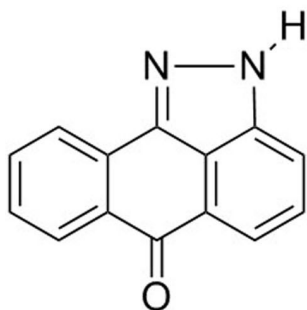

## Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at $-20^{\circ}\text{C}$ as supplied. Protect from prolonged exposure to light. For product expiry date, please contact <a href="mailto:techsupport@stemcell.com">techsupport@stemcell.com</a> .
Solubility:	<ul style="list-style-type: none"><li>· DMSO <math>\leq 90</math> mM</li><li>· Absolute ethanol <math>\leq 2.2</math> mM</li></ul> For example, to prepare a 10 mM stock solution in DMSO, resuspend 5 mg in 2.27 mL of fresh DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at  $-20^{\circ}\text{C}$ . Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### REPROGRAMMING

- Direct lineage reprogramming of fibroblasts to mature neurons, in combination with CHIR99021, RepSox, Forskolin, Gö6983, Valproic Acid and Y-27632 (Hu et al.).

### DIFFERENTIATION

- Inhibits BMP9-induced osteogenic differentiation in cultured mouse mesenchymal stem cells (MSCs) and in primary bone marrow stromal cells (Zhao et al.).
- Promotes adipogenic, but represses osteogenic differentiation of human MSCs (Bilkovski et al.; Liu et al.; Qiu et al.; Tominaga et al.).
- Causes cell death and inhibits neurogenesis when added during early stages of neuronal culture (Tiwari et al.).

## References

- Bennett BL et al. (2001) SP600125, an anthrapyrazolone inhibitor of Jun N-terminal kinase. *Proc Natl Acad Sci U S A* 98(24): 13681–6.
- Bilkovski R et al. (2010) Role of WNT-5a in the determination of human mesenchymal stem cells into preadipocytes. *J Biol Chem* 285(9): 6170–8.
- Hu W et al. (2015) Direct Conversion of Normal and Alzheimer’s Disease Human Fibroblasts into Neuronal Cells by Small Molecules. *Cell Stem Cell* 17(2): 204–212.
- Liu G et al. (2009) Canonical Wnts function as potent regulators of osteogenesis by human mesenchymal stem cells. *J Cell Biol* 185(1): 67–75.
- Qiu W et al. (2011) Activation of non-canonical Wnt/JNK pathway by Wnt3a is associated with differentiation fate determination of human bone marrow stromal (mesenchymal) stem cells. *Biochem Biophys Res Commun* 413(1): 98–104.
- Tiwari VK et al. (2012) A chromatin-modifying function of JNK during stem cell differentiation. *Nat Genet* 44(1): 94–100.
- Tominaga S et al. (2005) Negative regulation of adipogenesis from human mesenchymal stem cells by Jun N-terminal kinase. *Biochem Biophys Res Commun* 326(2): 499–504.
- Zhao Y et al. (2013) Activation of JNKs is essential for BMP9-induced osteogenic differentiation of mesenchymal stem cells. *BMB Rep* 46(8): 422–7.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.’S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2015 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.