

## Small Molecules

### Resveratrol

Antioxidant; COX-1 inhibitor

Catalog # 72862  
72864

100 mg  
500 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

Resveratrol is a potent phenolic antioxidant found in grapes and red wine that also has antiproliferative and anti-inflammatory activity (Rotondo et al.). Resveratrol is also a selective inhibitor of cyclooxygenase-1 (COX-1; Jang et al.). It inhibits COX and peroxidase activities of COX-1 with ED<sub>50</sub> values of 15 and 3.7 μM, respectively, with essentially no inhibition of the COX activity of COX-2.

Molecular Name:	Resveratrol
Alternative Names:	(E)-Resveratrol; trans-Resveratrol
CAS Number:	501-36-0
Chemical Formula:	C <sub>14</sub> H <sub>12</sub> O <sub>3</sub>
Molecular Weight:	228.2 g/mol
Purity:	≥ 98%
Chemical Name:	5-[(1E)-2-(4-hydroxyphenyl)ethenyl]-1,3-benzenediol
Structure:	

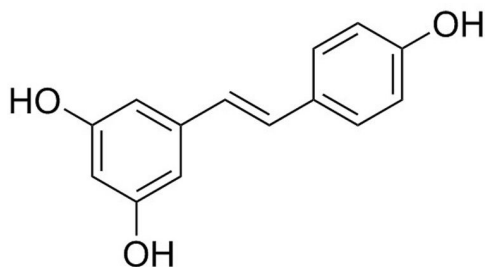

## Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt.
Solubility:	· DMSO ≤ 280 mM · Absolute ethanol ≤ 280 mM · PBS (pH 7.2) ≤ 430 μM For example, to prepare a 10 mM stock solution in DMSO, resuspend 100 mg in 43.8 mL of fresh DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### MAINTENANCE AND SELF-RENEWAL

- Inhibits differentiation of pig preadipocytes into adipocytes (Bai et al.).
- Promotes self-renewal of human mesenchymal stem cells (Peltz et al.).

### DIFFERENTIATION

- Induces osteogenic differentiation in human bone marrow-derived mesenchymal stem cell cultures (Boissy et al.; Dai et al.; Peltz et al.).
- Induces osteoblast differentiation while blocking adipocyte development in mouse mesenchymal stem cells and cultured primary rat bone marrow cells (Bäckesjö et al.).

## References

- Bäckesjö C-M et al. (2009) Activation of Sirt1 decreases adipocyte formation during osteoblast differentiation of mesenchymal stem cells. *Cells Tissues Organs* 189(1-4): 93–7.
- Bai L et al. (2008) Modulation of Sirt1 by resveratrol and nicotinamide alters proliferation and differentiation of pig preadipocytes. *Mol Cell Biochem* 307(1-2): 129–40.
- Boissy P et al. (2005) Resveratrol inhibits myeloma cell growth, prevents osteoclast formation, and promotes osteoblast differentiation. *Cancer Res* 65(21): 9943–52.
- Dai Z et al. (2007) Resveratrol enhances proliferation and osteoblastic differentiation in human mesenchymal stem cells via ER-dependent ERK1/2 activation. *Phytomedicine* 14(12): 806–14.
- Jang M et al. (1997) Cancer chemopreventive activity of resveratrol, a natural product derived from grapes. *Science* 275(5297): 218–20.
- Peltz L et al. (2012) Resveratrol exerts dosage and duration dependent effect on human mesenchymal stem cell development. *PLoS One* 7(5): e37162.
- Rotondo S et al. (1998) Effect of trans-resveratrol, a natural polyphenolic compound, on human polymorphonuclear leukocyte function. *Br J Pharmacol* 123(8): 1691–9.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

**This product is hazardous. Please refer to the Safety Data Sheet (SDS).**

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.