

## Small Molecules

### Sinomenine

Anti-inflammatory

Catalog # 72882  
72884

50 mg  
500 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

Sinomenine is a natural plant alkaloid commonly used to alleviate inflammation associated with rheumatoid arthritis (Wang and Li). It impairs signaling through Nuclear factor-kappa $\beta$  (NF- $\kappa$  $\beta$ ; Sun et al.; Wang and Li) and enhances the bioavailability of some compounds, at least in part through an inhibition of drug export by transporters like P-glycoprotein (Kesarwani et al.; Liu et al.). This product is supplied as the hydrochloride salt of the molecule.

Molecular Name:	Sinomenine (Hydrochloride)
Alternative Names:	Cocculine; Cucoline; NSC 76021
CAS Number:	6080-33-7
Chemical Formula:	C <sub>19</sub> H <sub>23</sub> NO <sub>4</sub> · HCl
Molecular Weight:	365.9 g/mol
Purity:	≥ 98%
Chemical Name:	9- $\alpha$ ,13- $\alpha$ ,14- $\alpha$ -Morphinan-6-one,7,8-didehydro-3,7-dimethoxy-4-hydroxy-17-methyl-,hydrochloride

Structure:

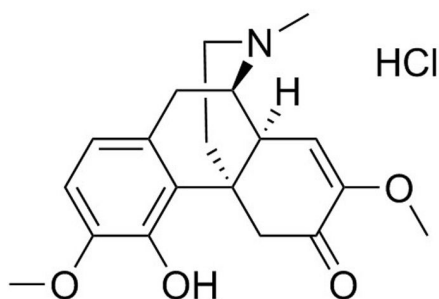

## Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com.
Solubility:	<ul style="list-style-type: none"><li>· PBS (pH 7.2) ≤ 13 mM</li><li>· DMSO ≤ 80 mM</li><li>· Absolute ethanol ≤ 13 mM</li></ul> For example, to prepare a 10 mM stock solution in DMSO, resuspend 50 mg in 13.7 mL of fresh DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### MAINTENANCE AND SELF-RENEWAL

- Promotes self-renewal in cultured human and mouse embryonic stem (ES) cells (Desbordes et al.).

## References

Desbordes SC et al. (2008) High-throughput screening assay for the identification of compounds regulating self-renewal and differentiation in human embryonic stem cells. *Cell Stem Cell* 2(6): 602–12.

Kesarwani K et al. (2013) Bioavailability enhancers of herbal origin: an overview. *Asian Pac J Trop Biomed* 3(4): 253–66.

Liu Z et al. (2014) Sinomenine sensitizes multidrug-resistant colon cancer cells (Caco-2) to doxorubicin by downregulation of MDR-1 expression. *PLoS One* 9(6): e98560.

Sun Y et al. (2014) A combination of sinomenine and methotrexate reduces joint damage of collagen induced arthritis in rats by modulating osteoclast-related cytokines. *Int Immunopharmacol* 18(1): 135–41.

Wang Q & Li X-K. (2011) Immunosuppressive and anti-inflammatory activities of sinomenine. *Int Immunopharmacol* 11(3): 373–6.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2015 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.