

Small Molecules

AG-490

Tyrosine kinase inhibitor; Inhibits EGFR, HER2, JAK2, JAK3, and STAT5a/b

Catalog # 72932
72934

10 mg
50 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

AG-490 is a member of the tyrphostin family of small molecule inhibitors of tyrosine kinases including EGFR (Gazit et al.), HER2 (Gazit et al.), JAK2 (Burger et al.), JAK3 (Brown et al.), and STAT5a/b (Wang et al.) with IC₅₀ values of 0.1, 13.5, 11, and 12 μM, respectively. AG-490 does not inhibit LCK, LYN, BTK, SYK, or SRC tyrosine kinases (Meydan et al.).

Molecular Name:	AG-490
Alternative Names:	Tyrphostin AG-490
CAS Number:	133550-30-8
Chemical Formula:	C ₁₇ H ₁₄ N ₂ O ₃
Molecular Weight:	294.3 g/mol
Purity:	≥ 98%
Chemical Name:	alpha-Cyano-(3,4-dihydroxy)-N-benzylcinnamide
Structure:	

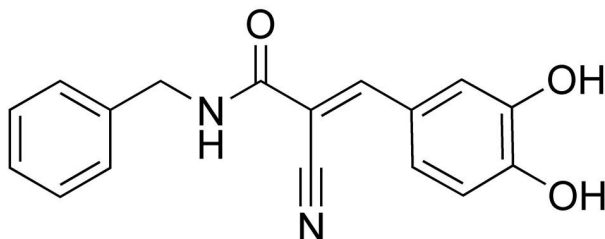

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt.
Solubility:	· DMSO ≤ 100 mM · Absolute ethanol ≤ 30 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 10 mg in 3.40 mL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

DIFFERENTIATION

- Induces differentiation of mouse (OG2) embryonic stem cells, with LIF (Catalog #78056) and BMP-4 (Chen et al.).

CANCER RESEARCH

- Blocks growth of acute lymphoblastic leukemia pre-B cells by inducing programmed cell death, via inhibition of JAK2 (Meydan et al.).
- Inhibits proliferation of cervical carcinoma cell lines (Soto-Cruz et al.).
- Inhibits STAT3 phosphorylation and induced translocation of beta-catenin to the cytoplasm in the colorectal cell line SW480 (Kawada et al.).
- Induces S phase arrest of GL15 glioblastoma cells via JAK2 inhibition (Sciacaluga et al.).
- Inhibits EGF-dependent proliferation of NIH 3T3 cell lines (Gazit et al.).
- Blocks IL-2 induced thymidine incorporation in T cell lines (Wang et al.).

References

- Brown GR et al. (2000) Naphthyl ketones: a new class of Janus kinase 3 inhibitors. *Bioorg Med Chem Lett* 10(6): 575–9.
- Burger R et al. (2009) Janus kinase inhibitor INCB20 has antiproliferative and apoptotic effects on human myeloma cells in vitro and in vivo. *Mol Cancer Ther* 8(1): 26–35.
- Chen S et al. (2006) Self-renewal of embryonic stem cells by a small molecule. *Proc Natl Acad Sci USA* 103(46): 17266–71.
- Gazit A et al. (1991) Tyrphostins. 2. Heterocyclic and alpha-substituted benzylidenemalononitrile tyrphostins as potent inhibitors of EGF receptor and ErbB2/neu tyrosine kinases. *J Med Chem* 34(6): 1896–907.
- Kawada M et al. (2006) Signal transducers and activators of transcription 3 activation is involved in nuclear accumulation of beta-catenin in colorectal cancer. *Cancer Res* 66(6): 2913–7.
- Meydan N et al. (1996) Inhibition of acute lymphoblastic leukaemia by a Jak-2 inhibitor. *Nature* 379(6566): 645–8.
- Sciacaluga M et al. (2007) Constitutive phosphorylation of Janus kinase 2 in the GL15 glioblastoma derived human cell line. *Oncol Rep* 17(1): 17–23.
- Soto-Cruz I et al. (2008) The Tyrphostin B42 Inhibits Cell Proliferation and HER-2 Autophosphorylation in Cervical Carcinoma Cell Lines. *Cancer Invest* 26(2): 136–44.
- Wang LH et al. (1999) JAK3, STAT, and MAPK signaling pathways as novel molecular targets for the tyrphostin AG-490 regulation of IL-2-mediated T cell response. *J Immunol* 162(7): 3897–904.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.