

Small Molecules

Rapamycin

Antibiotic; mTOR pathway inhibitor;
Inhibits FKBP-12

Catalog # 73362
73364

1 mg
10 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Rapamycin is a macrolide antibiotic and immunosuppressive compound that inhibits mammalian target of rapamycin (mTOR) signaling. It acts through formation of a complex with cytosolic FK-binding protein 12 (FKBP-12), which directly binds to mTOR complex 1 (mTORC1). Its immunosuppressive effects are mediated through inhibition of IL-2 signaling that is critical for T cell proliferation and activation (Gibbons et al.; Kay et al.). Rapamycin shows antifungal activity against *Candida albicans* and other fungi (Vézina et al.).

Molecular Name:	Rapamycin
Alternative Names:	AY 22989; Cypher; NSC 226080; Sirolimus
CAS Number:	53123-88-9
Chemical Formula:	C ₅₁ H ₇₉ NO ₁₃
Molecular Weight:	914.2 g/mol
Purity:	≥ 95%
Chemical Name:	(3S,6R,7E,9R,10R,12R,14S,15E,17E,19E,21S,23S,26R,27R,34aS)-9,10,12,13,14,21,22,23,24,25,26,27,32,33,34,34a-Hexadecahydro-9,27-dihydroxy-3-[(1R)-2-[(1S,3R,4R)-4-hydroxy-3-methoxycyclohexyl]-1-methylethyl]-10,21-dimethoxy-6,8,12,14,20,26-hexamethyl-23,27-ep

Structure:

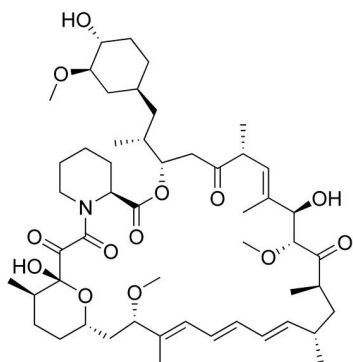

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt.
Solubility:	· DMSO ≤ 10 mM · Absolute ethanol ≤ 0.25 mM For example, to prepare a 5 mM stock solution in DMSO, resuspend 1 mg in 219 µL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

CANCER RESEARCH

- Inhibits growth of MDA-MB-468 human breast cancer cells in vitro, and inhibits tumor growth in a mouse xenograft model in vivo (Akcakanat et al.).
- Induces autophagy in malignant glioma cells (Takeuchi et al.).

References

- Akcakanat A et al. (2009) The rapamycin-regulated gene expression signature determines prognosis for breast cancer. *Mol Cancer* 8(1): 75.
- Gibbons et al. (2009) Mammalian target of rapamycin: Discovery of rapamycin reveals a signaling pathway important for normal and cancer cell growth. *Semin Oncol* 36(s3): s3-s17.
- Kay JE et al. (1991) Inhibition of T and B lymphocyte proliferation by rapamycin. *Immunology* 72(4): 544–9.
- Takeuchi H et al. (2005) Synergistic augmentation of rapamycin-induced autophagy in malignant glioma cells by phosphatidylinositol 3-kinase/protein kinase B inhibitors. *Cancer Res* 65(8): 3336–46.
- Vézina C et al. (1975) Rapamycin (AY-22,989), a new antifungal antibiotic. I. Taxonomy of the producing streptomycete and isolation of the active principle. *J Antibiot (Tokyo)* 28(10): 721–6.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.