

Antibodies

Mouse IgG3, kappa Isotype Control Antibody, Clone MG3-35, Alexa Fluor® 488

Mouse monoclonal IgG3, kappa isotype control antibody, Alexa Fluor® 488-conjugated

Catalog #60073AD
#60073AD.1

100 µg 0.5 mg/mL
25 µg 0.5 mg/mL



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

The MG3-35 antibody (IgG3, kappa) is suitable for use as an isotype-matched control antibody in several applications to estimate the degree of non-specific binding by an antigen-specific antibody. Ideally, the isotype control should have the same subclass of heavy chain (IgA, IgD, IgE, IgG, or IgM) and light chain (kappa or lambda) as the specific antibody being used. If a conjugated antibody is used, an isotype control conjugated to the same molecule (e.g. fluorochrome) should be chosen. The use of an appropriate isotype control helps confirm the specificity of the antigen-specific antibody and indicates non-specific binding that may result from binding to Fc receptors or other cell components. The MG3-35 antibody recognizes keyhole limpet hemocyanin (KLH) and has unknown binding specificity, having been screened on a variety of activated, resting, live, and fixed tissues from several species, including mouse, rat, human, and non-human primates.

| | |
|----------------------|----------------------|
| Target Antigen Name: | IgG3 Isotype Control |
| Alternative Names: | Not applicable |
| Gene ID: | Not applicable |
| Species Reactivity: | Not applicable |
| Host Species: | Mouse |
| Clonality: | Monoclonal |
| Clone: | MG3-35 |
| Isotype: | IgG3, kappa |
| Immunogen: | Trinitrophenol + KLH |
| Conjugate: | Alexa Fluor® 488 |

Applications


| | |
|-----------------------|---|
| Verified: | FC, ICC, IF |
| Reported: | FC, ICC, IF |
| Special Applications: | This antibody clone has been verified for use as an isotype control antibody for assessing non-specific binding to cells in flow cytometry and immunofluorescence microscopy applications (surface and intracellular staining). |

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FACS: Fluorescence-activated cell sorting; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; RIA: Radioimmunoassay; WB: Western blotting

Properties

| | |
|------------------------|---|
| Formulation: | Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide |
| Purification: | The antibody was purified by affinity chromatography and conjugated with Alexa Fluor® 488 under optimal conditions. |
| Stability and Storage: | Product stable at 2 - 8°C when stored undiluted. Do not freeze. Protect product from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com. |
| Directions for Use: | The suggested use of this antibody is at concentrations comparable to those of the specific antibody of interest. |

Data

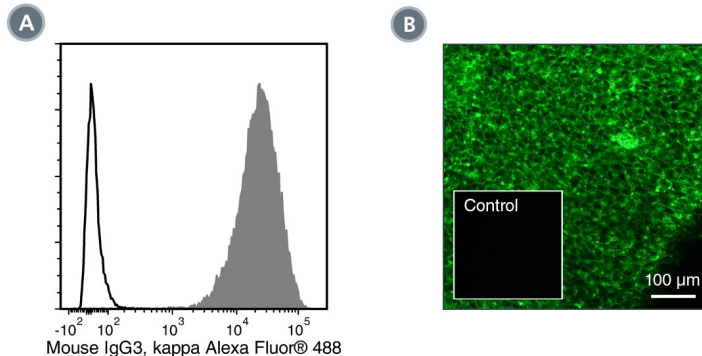

(A) Flow cytometry analysis of human induced pluripotent stem cells labeled with Mouse IgG3, kappa Isotype Control Antibody, Clone MG3-35, Alexa Fluor® 488 (solid line histogram), or a mouse IgG3, kappa positive control antibody (Anti-Human SSEA-4 Antibody, Clone MC-813-70, Alexa Fluor® 488; Catalog #60062AD) (filled histogram).

(B) Human embryonic stem cells were cultured in mTeSR™1 (Catalog #85850) on Corning® Matrigel®-coated glass slides, then fixed and labeled with Anti-Human SSEA-4 Antibody, Clone MC-813-70, Alexa Fluor® 488. Inset shows cells labeled with Mouse IgG3, kappa Isotype Control Antibody, Clone MG3-35, Alexa Fluor® 488.

Related Products

For a complete list of antibodies, including other conjugates, sizes and clones, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/antibodies or contact us at techsupport@stemcell.com.

References

1. Ameres S et al. (2014) CD8 T cell-evasive functions of human cytomegalovirus display pervasive MHC allele specificity, complementarity, and cooperativity. *J Immunol* 192(12): 5894–905. (FC)
2. Gao C et al. (2014) Carbohydrate sequence of the prostate cancer-associated antigen F77 assigned by a mucin O-glycome designer array. *J Biol Chem* 289(23): 16462–77. (FA)
3. De Riva A et al. (2013) Accelerated turnover of MHC class II molecules in nonobese diabetic mice is developmentally and environmentally regulated in vivo and dispensable for autoimmunity. *J Immunol* 190(12): 5961–71. (FC)
4. Fewou SN et al. (2012) Anti-ganglioside antibody internalization attenuates motor nerve terminal injury in a mouse model of acute motor axonal neuropathy. *J Clin Invest* 122(3): 1037–51. (FA)
5. Ribot JC et al. (2012) B7-CD28 costimulatory signals control the survival and proliferation of murine and human gamma delta T cells via IL-2 production. *J Immunol* 189(3): 1202–8. (FC)
6. Mitchell J & Sullam PM. (2009) Streptococcus mitis phage-encoded adhesins mediate attachment to α 2-8-linked sialic acid residues on platelet membrane gangliosides. *Infect Immun* 77(8): 3485–90. (FA)
7. Duan J et al. (2008) Microbial carbohydrate depolymerization by antigen-presenting cells: deamination prior to presentation by the MHC class II pathway. *Proc Natl Acad Sci USA* 105(13): 5183–8. (FA)

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, and EasySep are trademarks of STEMCELL Technologies Canada Inc. mTeSR is a trademark of WARF. Corning and Matrigel are registered trademarks of Corning Incorporated. All other trademarks are the property of their respective holders. Alexa Fluor® is a registered trademark of Life Technologies Corporation. Antibodies conjugated to Alexa Fluor® are licensed for internal research use only and sale is expressly conditioned on the buyer not using the antibody for manufacturing, performing a service or medical test, or otherwise generating revenue. For use other than research, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.