

## Cytokines

### Mouse Recombinant Flt3/Flk-2 Ligand

Fms-like tyrosine kinase 3/fetal liver kinase-2

Catalog #	78011.1	10 µg
	78011	100 µg
	78011.2	1000 µg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

Flt3/Flk-2 (Fms-like tyrosine kinase 3/fetal liver kinase-2) ligand is a hematopoietic cytokine that plays an important role as a co-stimulatory factor in the proliferation, differentiation, and survival of hematopoietic stem and progenitor cells and the development of the immune system (Lyman et al.; Rosnet et al.). Flt3/Flk-2 ligand, together with stem cell factor and thrombopoietin, is commonly used to promote expansion of primitive hematopoietic cells in culture. In combination with myeloid cytokines such as GM-CSF, G-CSF, or M-CSF, Flt3/Flk-2 ligand enhances the growth and numbers of clonogenic myeloid progenitor cells. In synergy with IL-3, IL-4, IL-7, IL-11, IL-12, IL-15, GM-CSF, and TNF- $\alpha$ , Flt3/Flk-2 ligand regulates the development of various lymphoid progenitor cells, including dendritic cell, B cell, T cell, and NK cell progenitors. In contrast, Flt3/Flk-2 ligand has no significant effect on erythropoiesis or megakaryopoiesis (Drexler & Quentmeier; Wodnar-Filipowicz).

Flt3/Flk-2 ligand exists as membrane-bound and soluble isoforms. Both isoforms are biologically active and signal through the class III tyrosine kinase receptor (Flt3/Flk-2, CD135; Rosnet et al.). Flt3/Flk-2 ligand is produced by a variety of cell types, including uncommitted and committed hematopoietic cells and stromal fibroblasts, whereas expression of the Flt3/Flk-2 receptor is restricted to hematopoietic stem and progenitor cells. Flt3/Flk-2 receptor is also expressed outside the hematopoietic system in the brain, placenta, and testis (Drexler & Quentmeier; Hannum et al.).

## Product Information

Alternative Names:	FL, FLT3L, Flt3-L, Flt3 Ligand, Fms-like tyrosine kinase 3 ligand
Accession Number:	P49772
Amino Acid Sequence:	MTPDCYFSHS PISSNFKVKF RELTDHLLKD YPVTAVNLQ DEKHCKALWS LFLAQRWIEQ LKTVAGSKMQ TLLEDVNTEI HFVTSCTFQP LPECLRFVQT NISHLLKDTC TQLLALKPCI GKACQNF SRC LEVQCQP DSS TLLPPRSPIA LEATELPEPR PRQ
Predicted Molecular Mass:	18.6 kDa
Species:	Mouse
Cross Reactivity:	Human, Rat
Formulation:	Lyophilized from a sterile-filtered solution containing sodium phosphate, pH 7.5.
Source:	E. coli


## Specifications

Activity:	The specific activity is $\geq 1 \times 10^5$ units/mg ( $EC_{50} \leq 10$ ng/mL) as determined by a cell proliferation assay using AML5 cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic Limulus amoebocyte lysate (LAL) analysis and is $\leq 1$ EU/ $\mu$ g protein.

## Preparation and Storage

Storage:	Store at $-20^{\circ}\text{C}$ to $-80^{\circ}\text{C}$ .
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. OPTIONAL: After reconstitution, if product will not be used immediately, dilute with concentrated bovine serum albumin (BSA) to a final BSA concentration of 0.1%. The effect of storage of stock solution on product performance should be tested for each application. As a general guide, do not store at $2 - 8^{\circ}\text{C}$ for more than 1 month or at $-80^{\circ}\text{C}$ for more than 3 months. Avoid repeated freeze-thaw cycles.

## Data

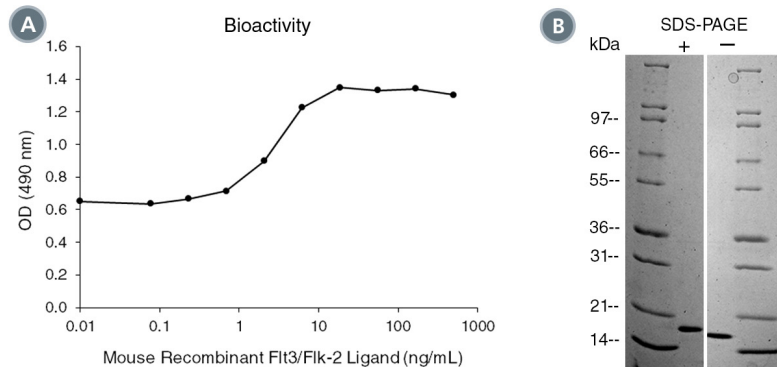

(A) The biological activity of Mouse Recombinant Flt3/Flk-2 Ligand was tested by its ability to promote the proliferation of AML5 cells. Cell proliferation was measured after 65 hours of culture using a fluorometric assay method. The EC<sub>50</sub> is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC<sub>50</sub> in the above example is 1.3 - 2.0 ng/mL.

(B) 1 µg of Mouse Recombinant Flt3/Flk-2 Ligand was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Mouse Recombinant Flt3/Flk-2 Ligand has a predicted molecular mass of 18.6 kDa.

## Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, visit [www.stemcell.com/cytokines](http://www.stemcell.com/cytokines) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

## References

- Drexler HG & Quentmeier H. (2004) Mini Review FLT3: receptor and ligand. *Growth Factors* 22(2): 71–3.
- Hannum C et al. (1994) Ligand for FLT3/FLK2 receptor tyrosine kinase regulates growth of haematopoietic stem cells and is encoded by variant RNAs. *Nature* 368(6472): 643–8.
- Lyman SD et al. (1993) Molecular cloning of a ligand for the flt3/flk-2 tyrosine kinase receptor: a proliferative factor for primitive hematopoietic cells. *Cell* 75(6): 1157–67.
- Rosnet O et al. (1991) Murine Flt3, a gene encoding a novel tyrosine kinase receptor of the PDGFR/CSF1R family. *Oncogene* 6(9): 1641–50.
- Wodnar-Filipowicz A. (2003) Flt3 ligand: role in control of hematopoietic and immune functions of the bone marrow. *News Physiol Sci* 18: 247–51.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.