

# Cytokines

## Human Recombinant TPO



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #	78070.1	10 µg
	78070	100 µg
	78070.2	1000 µg

## Product Description

Thrombopoietin (TPO) is a key regulator of megakaryocytopoiesis and thrombopoiesis in vitro and in vivo. TPO stimulates the proliferation and maturation of megakaryocytes and has an important role in regulating the level of circulating platelets in vivo (de Sauvage et al.; Bartley et al.; Foster et al.; Sohma et al.). TPO also promotes the survival, self-renewal, and expansion of hematopoietic stem cells and primitive multi-lineage progenitor cells and is commonly used with other cytokines such as stem cell factor and Flt3-Ligand to promote expansion of primitive hematopoietic cells in culture (Hitchcock & Kaushansky). TPO is a ligand for the cytokine receptor Mpl, which is expressed at all stages of megakaryopoiesis from hematopoietic stem and progenitor cells to mature platelets.

## Product Information

Alternative Names:	c-MPL ligand, Megakaryocyte colony-stimulating factor, MGDF
Accession Number:	P40225
Amino Acid Sequence:	MSPAPPACDL RVLKLLRDS HVLHSRLSQC PEVHPLPTPV LLPVDFSLG EWKTQMEETK AQDILGAVTL LLEGVMAARG QLGPTCLSSL LGQLSGQVRL LLAGLQSLLG TQLPPQGRIT AHKDPNAIFL SFQHLLRGKV RFLMLVGGST LCVRRAPPTT AVPSRTSLVL TLNEL
Predicted Molecular Mass:	18.8 kDa
Species:	Human
Cross Reactivity:	Mouse, Rat, Monkey
Formulation:	Lyophilized from a sterile filtered aqueous solution containing sodium phosphate, pH 7.5.
Source:	E. coli


## Specifications

Activity:	The specific activity is $\geq 2 \times 10^5$ units/mg ( $EC_{50} \leq 5$ ng/mL) as determined by a cell proliferation assay using MO7e cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic limulus amoebocyte lysate (LAL) analysis and is $\leq 1$ EU/µg protein.

## Preparation and Storage

Storage:	Store at -20°C to -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Reconstitution:	Centrifuge vial before opening. Resuspend the product in sterile water containing 0.1% bovine serum albumin (BSA) to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. Store at 2 - 8°C for up to 1 month or at -20°C to -80°C for up to 3 months. Avoid repeated freeze-thaw cycles. NOTE: If reconstituted product will be used immediately BSA is not required.

## Data

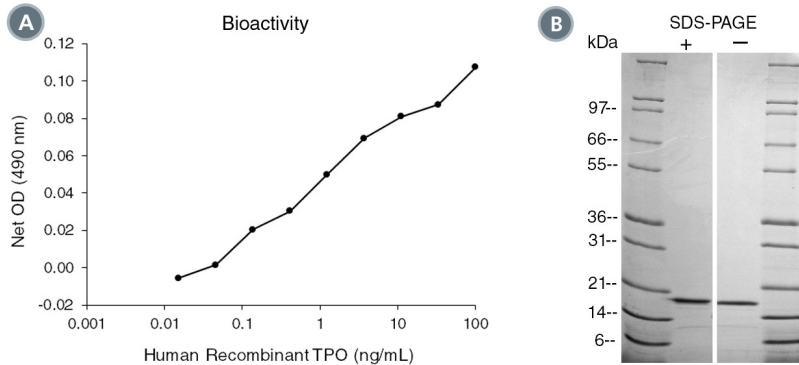

(A) The biological activity of Human Recombinant TPO was tested by its ability to promote the proliferation of MO7e cells growing in the presence of IL-3 and SCF. Cell proliferation was measured after 72 hours of culture using a fluorometric assay method. The EC<sub>50</sub> is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC<sub>50</sub> in the above example is 1.3 - 2 ng/mL.

(B) 1 µg of Human Recombinant TPO was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human Recombinant TPO has a predicted molecular mass of 18.8 kDa.

## Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, please visit our website at [www.stemcell.com/cytokines](http://www.stemcell.com/cytokines) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

## References

- Bartley TD et al. (1994) Identification and cloning of a megakaryocyte growth and development factor that is a ligand for the cytokine receptor Mpl. *Cell* 77(7): 1117-24.
- Foster DC et al. (1994) Human thrombopoietin: gene structure, cDNA sequence, expression, and chromosomal localization. *Proc Natl Acad Sci U S A* 91(26): 13023-7.
- Hitchcock IS & Kaushansky K. (2014) Thrombopoietin from beginning to end. *Br J Haematol* 165(2): 259-68.
- de Sauvage FJ et al. (1994) Stimulation of megakaryocytopoiesis and thrombopoiesis by the c-Mpl ligand. *Nature* 369(6481): 533-8.
- Sohma Y et al. (1994) Molecular cloning and chromosomal localization of the human thrombopoietin gene. *FEBS Lett* 353(1): 57-61.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.