

Small Molecules

8-Bromo-cAMP

cAMP pathway activator; Activates cAMP-dependent kinase

Catalog # 73602
73604

10 mg
100 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

8-Bromo-cAMP is a membrane-permeable cAMP derivative. It can activate cAMP-dependent protein kinase, with long-acting effects due to its resistance to cAMP phosphodiesterase (Schwede et al.). It can be used to study calcium-mediated pathways ($IC_{50} = 0.84 \text{ mM}$; Xaus et al.).

Molecular Name:	8-Bromo-cAMP
Alternative Names:	8-BrcAMP; 8-Bromoadenosine 3',5'-cyclic monophosphate; NSC 171719
CAS Number:	23583-48-4
Chemical Formula:	$C_{10}H_{11}BrN_5O_6P$
Molecular Weight:	408.1 g/mol
Purity:	$\geq 95\%$
Chemical Name:	8-bromo-adenosine cyclic 3',5'-(hydrogen phosphate)
Structure:	

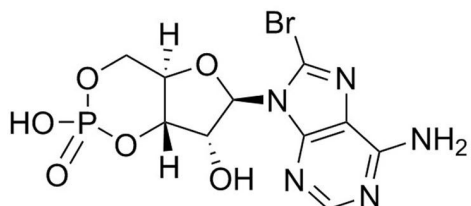

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect product from prolonged exposure to light. For long-term storage store with a desiccant. For product expiry date, please contact techsupport@stemcell.com .
Solubility:	· PBS (pH 7.2) $\leq 7.4 \text{ mM}$ For example, to prepare a 5 mM stock solution in PBS, resuspend 10 mg in 4.9 mL of PBS.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in PBS at -20°C . Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use.

Published Applications

REPROGRAMMING

- Improves the reprogramming efficiency of human neonatal foreskin fibroblast (HFF1) cells, in combination with valproic acid (Wang & Adjaye).

IMMUNOLOGY

- Inhibits M-CSF-dependent proliferation of macrophages (Xaus et al.).
- Protects neutrophils against TNF- α -induced apoptosis (Krakstad).

CANCER RESEARCH

- Induces a proliferative response in an IL-3-dependent leukemic cell line (Barge et al.).
- Induces membrane depolarization in pancreatic cancer cell lines (Sorio et al.).

References

Barge RM et al. (1997) 8-Bromo-cAMP induces a proliferative response in an IL-3 dependent leukemic cell line and activates Erk 1,2 via a Shc-independent pathway. *Biochim Biophys Acta* 1355(2): 141–6.

Krakstad C. (2004) cAMP protects neutrophils against TNF- α -induced apoptosis by activation of cAMP-dependent protein kinase, independently of exchange protein directly activated by cAMP (Epac). *J Leukoc Biol* 76(3): 641–647.

Schwede F et al. (2000) Cyclic nucleotide analogs as biochemical tools and prospective drugs. *Pharmacol Ther* 87(2-3): 199–226.

Sorio C et al. (2011) Defective CFTR expression and function are detectable in blood monocytes: development of a new blood test for cystic fibrosis. *PLoS One* 6(7): e22212.

Wang Y & Adjaye J. (2011) A cyclic AMP analog, 8-Br-cAMP, enhances the induction of pluripotency in human fibroblast cells. *Stem Cell Rev* 7(2): 331–41.

Xaus J et al. (1999) Adenosine inhibits macrophage colony-stimulating factor-dependent proliferation of macrophages through the induction of p27kip-1 expression. *J Immunol* 163(8): 4140–9.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.