

# Small Molecules

CAY10603

Epigenetic modifier; Inhibits histone deacetylase (HDAC) 6

Catalog # 73582

1 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

CAY10603 is a selective and potent inhibitor of histone deacetylase 6 (HDAC6;  $IC_{50} = 0.002$  nM; Kozikowski et al.). HDAC6 deacetylates lysine residues on the N-terminal part of histones. HDAC6 is also known to regulate heat shock protein 90 (hsp90) via hyperacetylation (Rao et al.).

Molecular Name:	CAY10603
Alternative Names:	HDAC6 inhibitor
CAS Number:	1045792-66-2
Chemical Formula:	$C_{22}H_{30}N_4O_6$
Molecular Weight:	446.5 g/mol
Purity:	≥ 95%
Chemical Name:	N-[4-[3-[[[7-(hydroxyamino)-7-oxoheptyl]amino]carbonyl]-5-isoxazolyl]phenyl]-1,1-dimethylethylester, carbamic acid

Structure:

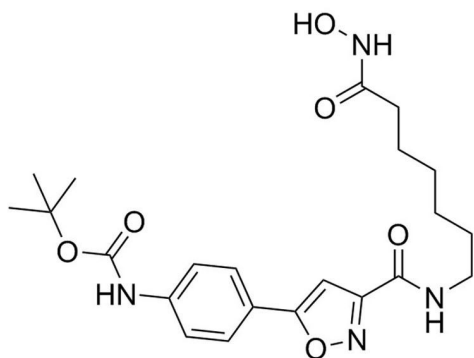

## Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at $-20^{\circ}C$ as supplied. Protect from prolonged exposure to light. For long-term storage store with a desiccant. Stable as supplied for 12 months from date of receipt.
Solubility:	· DMSO $\leq 30$ mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 224 $\mu$ L of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at  $-20^{\circ}C$ . Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### CANCER RESEARCH

· Inhibits growth of several pancreatic cancer cell lines ( $IC_{50} = 0.1 - 1 \mu M$ ) and SCCOHT, an ovarian cancer cell line (Kozikowski et al.; Wang et al.).

## References

Kozikowski AP et al. (2008) Use of the nitrile oxide cycloaddition (NOC) reaction for molecular probe generation: a new class of enzyme selective histone deacetylase inhibitors (HDACIs) showing picomolar activity at HDAC6. *J Med Chem* 51(15): 4370–3.

Rao R et al. (2008) HDAC6 inhibition enhances 17-AAG--mediated abrogation of hsp90 chaperone function in human leukemia cells. *Blood* 112(5): 1886–93.

Wang Y et al. (2015) Abstract 5381: Therapeutic potential of HDAC inhibitors in small cell carcinoma of the ovary, hypercalcemic type (SCCOHT). *Cancer Res* 75(15 Supplement): 5381.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.