

# Small Molecules

## RepSox

Activin/Nodal/TGF $\beta$  pathway inhibitor;  
Inhibits ALK5

Catalog # 73792  
73794

1 mg  
10 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

RepSox is a cell-permeable, selective inhibitor of the TGF- $\beta$  type 1 receptor (TGF $\beta$ RI) ALK5 ( $IC_{50}$  = 4, 18, and 23 nM for ALK5 autophosphorylation, TGF- $\beta$  cellular assay, and ALK5 binding in HepG2 cells, respectively; Gellibert et al.). This inhibitor demonstrates less potent activity ( $IC_{50}$  > 16  $\mu$ M) against 9 related kinases, including p38 MAPK and GSK3 (Gellibert et al.).

|                    |  |
|--------------------|--|
| Molecular Name:    | RepSox   |
| Alternative Names: | ALK5 Inhibitor II; E-616452; SJN 2511                          |
| CAS Number:        | 446859-33-2  |
| Chemical Formula:  | C <sub>17</sub> H <sub>13</sub> N <sub>5</sub>                 |
| Molecular Weight:  | 287.3 g/mol  |
| Purity:            | ≥ 98%  |
| Chemical Name:     | 2-[5-(6-methylpyridin-2-yl)-1H-pyrazol-4-yl]-1,5-naphthyridine |
| Structure:         |  |

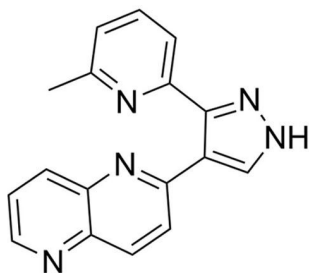

## Properties

|                      |  |
|----------------------|--|
| Physical Appearance: | A crystalline solid  |
| Storage:             | Product stable at -20°C as supplied. Protect product from prolonged exposure to light. For long-term storage store with a desiccant.<br>Stable as supplied for 12 months from date of receipt. |
| Solubility:          | · DMSO ≤ 30 mM<br>· Absolute ethanol ≤ 3.5 mM<br>For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 348 $\mu$ L of DMSO.  |

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### REPROGRAMMING

- Enhances reprogramming of mouse embryonic fibroblasts that have been transduced with OCT4, KLF4, and c-MYC (Ichida et al.; Subramanyam et al.).
- Direct lineage reprogramming of fibroblasts to mature neurons, in combination with CHIR99021 (Catalog #72052), Valproic Acid (Catalog #72292), Forskolin (Catalog #72112), SP600125 (Catalog #72642), Gö6983 (Catalog #72462), and Y-27632 (Catalog #72302) (Hu et al.).

### DIFFERENTIATION

- Alone or in combination with Forskolin, Dexamethasone (Catalog #72092), and Nicotinamide (Catalog #07154), induces differentiation of human pancreatic progenitor cells into insulin-producing cells (Kunisada et al.; Rezania et al.).

## References

- Gellibert F et al. (2004) Identification of 1,5-naphthyridine derivatives as a novel series of potent and selective TGF-beta type I receptor inhibitors. *J Med Chem* 47(18): 4494–506.
- Hu W et al. (2015) Direct conversion of normal and Alzheimer's disease human fibroblasts into neuronal cells by small molecules. *Cell Stem Cell* 17(2): 204–12.
- Ichida JK et al. (2009) A small-molecule inhibitor of TGF-beta signaling replaces Sox2 in reprogramming by inducing Nanog. *Cell Stem Cell* 5(5): 491–503.
- Kunisada Y et al. (2012) Small molecules induce efficient differentiation into insulin-producing cells from human induced pluripotent stem cells. *Stem Cell Res* 8(2): 274–84.
- Rezania A et al. (2011) Production of functional glucagon-secreting  $\alpha$ -cells from human embryonic stem cells. *Diabetes* 60(1): 239–47.
- Subramanyam D et al. (2011) Multiple targets of miR-302 and miR-372 promote reprogramming of human fibroblasts to induced pluripotent stem cells. *Nat Biotechnol* 29(5): 443–8.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

**This product is hazardous. Please refer to the Safety Data Sheet (SDS).**

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.