

Small Molecules

Sirtinol

Epigenetic modifier; Inhibitor of Sir2p, SIRT1, and SIRT2

Catalog # 73824

10 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Sirtinol is a cell-permeable inhibitor of sirtuin NAD⁺-dependent histone deacetylases (HDACs), including yeast Sir2p (IC₅₀ = 68 μM; Grozinger et al.), human SIRT1 (IC₅₀ = 131 μM; Schuetz et al.), and SIRT2 (IC₅₀ = 38 μM; Grozinger et al.). It does not inhibit human HDAC1 activity (Grozinger et al.). Sirtinol also activates the auxin signaling pathway in Arabidopsis (Dai et al.).

Molecular Name:	Sirtinol
Alternative Names:	Sir Two Inhibitor Naphthol
CAS Number:	410536-97-9
Chemical Formula:	C ₂₆ H ₂₂ N ₂ O ₂
Molecular Weight:	394.5 g/mol
Purity:	≥ 98%
Chemical Name:	2-[[[(Z)-(2-oxonaphthalen-1-ylidene)methyl]amino]-N-(1-phenylethyl)benzamide
Structure:	

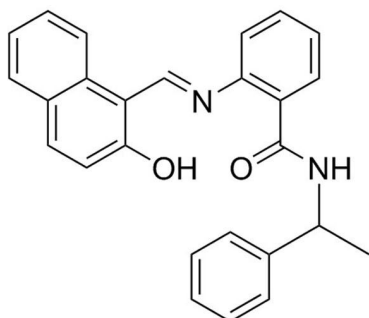

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect product from prolonged exposure to light. For long-term storage store with a desiccant. Stable as supplied for 12 months from date of receipt.
Solubility:	· DMSO ≤ 25 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 253 μL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

IMMUNOLOGY

- Suppresses inflammatory signaling in human dermal microvascular endothelial cells (Orecchia et al.).

CANCER RESEARCH

- Enhances apoptosis and cell death in leukemic cells, but not in healthy leukocytes or hematopoietic progenitor cells (Cea et al.).
- Induces senescence-like growth arrest in human cancer cell lines MCF-7 and H1299 (Ota et al.; Wang et al.).

References

- Cea M et al. (2011) Synergistic interactions between HDAC and sirtuin inhibitors in human leukemia cells. *PLoS One* 6(7): e22739.
- Dai X et al. (2005) Genetic and chemical analyses of the action mechanisms of sirtinol in Arabidopsis. *Proc Natl Acad Sci USA* 102(8): 3129–34.
- Grozinger CM et al. (2001) Identification of a class of small molecule inhibitors of the sirtuin family of NAD-dependent deacetylases by phenotypic screening. *J Biol Chem* 276(42): 38837–43.
- Orecchia A et al. (2011) Sirtinol treatment reduces inflammation in human dermal microvascular endothelial cells. *PLoS One* 6(9): e24307.
- Ota H et al. (2006) Sirt1 inhibitor, Sirtinol, induces senescence-like growth arrest with attenuated Ras-MAPK signaling in human cancer cells. *Oncogene* 25(2): 176–85.
- Schuetz A et al. (2007) Structural basis of inhibition of the human NAD⁺-dependent deacetylase SIRT5 by suramin. *Structure* 15(3): 377–89.
- Wang J et al. (2012) Sirtinol, a class III HDAC inhibitor, induces apoptotic and autophagic cell death in MCF-7 human breast cancer cells. *Int J Oncol* 41(3): 1101–9.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.