

Cytokines

Human Recombinant FGF-10 (KGF-2)

Fibroblast growth factor 10



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #	78037	10 µg
	78037.1	50 µg
	78037.2	500 µg
	78037.3	1000 µg

Product Description

Fibroblast growth factor 10 (FGF-10) is a member of the fibroblast growth factor (FGF) family predominantly expressed by mesenchymal fibroblasts during embryonic development (Emoto et al.; Igarashi et al.). It binds with high affinity to fibroblast growth factor receptor 2-IIIb (FGFR2-IIIb), and also has a weaker affinity for FGFR1-IIIb (Beer et al.). FGF-10 and FGF-7 have similar receptor binding properties and target cell specificities but are differentially regulated by components of the extracellular matrix (Emoto et al.; Igarashi et al.). FGF-10 has been shown to mediate epithelial-mesenchymal interactions, which are essential to lung development (Sekine et al.; Ware & Matthay). FGF-10 also has a role in mobilization and proliferation of lung-resident mesenchymal stem cells (MSCs) and protection and repair against acute lung injury (Tong et al.; Ware & Matthay), as well as endodermal differentiation of human pluripotent stem cells to insulin-producing pancreatic-like cells (Takeuchi et al.).

Product Information

Alternative Names:	FGFA, Fibroblast growth factor-10, Keratinocyte growth factor-2
Accession Number:	O15520
Amino Acid Sequence:	LGQDMVSPEA TNSSSSFSS PSSAGRHVRS YNHLQGDVRW RKLFSFTKYF LKIEKNGKVS GTKKENCPYS ILEITSVEIG VVAVKAINSN YYLAMNKKKGK LYGSKEFNND CKLKERIEEN GYNTYASFNW QHNGRQMYVA LNGKGAPRRG QKTRRKNTSA HFLPMVVHS
Predicted Molecular Mass:	19.3 kDa
Species:	Human
Cross Reactivity:	Mouse, Rat
Formulation:	Lyophilized after dialysis against phosphate-buffered saline.
Source:	E. coli


Specifications

Activity:	The specific activity is $\geq 5 \times 10^4$ units/mg ($EC_{50} \leq 20$ ng/mL) as determined by a cell proliferation assay using 4MBr-5 cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic Limulus amoebocyte lysate (LAL) analysis and is ≤ 0.2 EU/µg protein.

Preparation and Storage

Storage:	Store at -80°C .
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. As a general guide, do not store at $2 - 8^{\circ}\text{C}$ for more than 1 week or at -20°C for more than 3 months. Avoid repeated freeze-thaw cycles.

Data

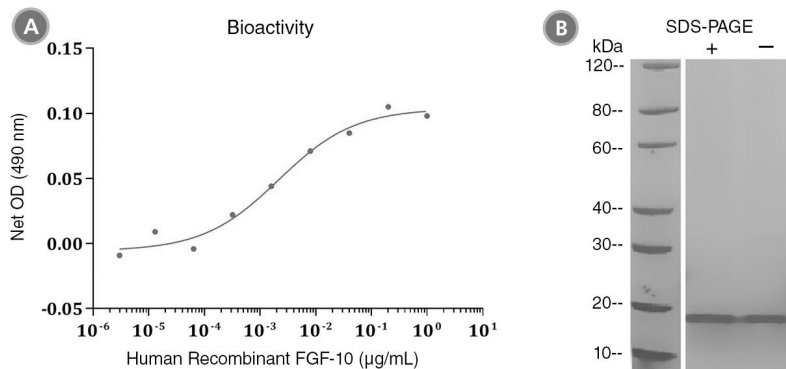

(A) The biological activity of Human Recombinant FGF-10 (KGF-2) was tested by its ability to promote the proliferation of 4MBr-5 cells. Cell proliferation was measured using a fluorometric assay method. The EC₅₀ is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC₅₀ in the above example is less than 20 ng/mL.

(B) 2 µg of Human Recombinant FGF-10 (KGF-2) was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human Recombinant FGF-10 (KGF-2) has a predicted molecular mass of 19.3 kDa.

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, visit www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Beer HD et al. (2000) Fibroblast growth factor (FGF) receptor 1-IIIb is a naturally occurring functional receptor for FGFs that is preferentially expressed in the skin and the brain. *J Biol Chem* 275(21): 16091–7.
- Emoto H et al. (1997) Structure and expression of human fibroblast growth factor-10. *J Biol Chem* 272(37): 23191–4.
- Igarashi M et al. (1998) Characterization of recombinant human fibroblast growth factor (FGF)-10 reveals functional similarities with keratinocyte growth factor (FGF-7). *J Biol Chem* 273(21): 13230–5.
- Sekine K et al. (1999) Fgf10 is essential for limb and lung formation. *Nat Genet* 21(1): 138–41.
- Takeuchi H et al. (2014) Endodermal differentiation of human pluripotent stem cells to insulin-producing cells in 3D culture. *Sci Rep* 4: 4488.
- Tong L et al. (2016) Fibroblast growth factor-10 (FGF-10) mobilizes lung-resident mesenchymal stem cells and protects against acute lung injury. *Sci Rep* 6: 21642.
- Ware LB & Matthay MA. (2002) Keratinocyte and hepatocyte growth factors in the lung: roles in lung development, inflammation, and repair. *Am J Physiol Lung Cell Mol Physiol* 282(5): L924–40.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2019 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.