

## Small Molecules

### SAHA

Epigenetic modifier; Inhibits HDAC1 and HDAC3

Catalog # 73902  
73904

100 mg  
500 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

SAHA (suberoylanilide hydroxamic acid) is a potent inhibitor of class I and class II histone deacetylases (HDACs). This compound inhibits the activities of HDAC1 and HDAC3 ( $IC_{50}$  = 10 nM and 20 nM, respectively; Richon et al.). SAHA exhibits antiproliferative properties and can cause cell cycle arrest at G1 (Marks & Breslow).

|                    |   |
|--------------------|---|
| Molecular Name:    | SAHA  |
| Alternative Names: | Suberoylanilide Hydroxamic Acid; Vorinostat |
| CAS Number:        | 149647-78-9                                 |
| Chemical Formula:  | $C_{14}H_{20}N$                             |
| Molecular Weight:  | 264.3 g/mol                                 |
| Purity:            | ≥ 98%                                       |
| Chemical Name:     | N-Hydroxy-N'-phenyloctanediamide            |
| Structure:         |   |

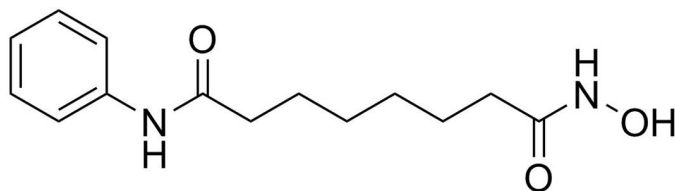

## Properties

|                      |  |
|----------------------|--|
| Physical Appearance: | A crystalline solid  |
| Storage:             | Product stable at -20°C as supplied. Protect product from prolonged exposure to light. For long-term storage store with a desiccant.<br>Stable as supplied for 12 months from date of receipt. |
| Solubility:          | · DMSO ≤ 75 mM<br>· Absolute ethanol ≤ 0.9 mM<br>For example, to prepare a 10 mM stock solution in DMSO, resuspend 10 mg in 3.78 mL of DMSO.   |

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### REPROGRAMMING

- Improves reprogramming efficiency from human and mouse somatic to induced pluripotent stem cells (Huangfu et al.).

### CANCER RESEARCH

- Downregulates homologous recombination and leads to apoptosis in ovarian cancer cells when used in combination with olaparib (Konstantinopoulos et al.).
- Downregulates Nanog expression, which leads to inhibition of proliferation in head and neck cancer stem cells (Kumar et al.).

## References

- Huangfu D et al. (2008) Induction of pluripotent stem cells by defined factors is greatly improved by small-molecule compounds. *Nat Biotechnol* 26(7): 795–7.
- Konstantinopoulos PA et al. (2014) Suberoylanilide hydroxamic acid (SAHA) enhances olaparib activity by targeting homologous recombination DNA repair in ovarian cancer. *Gynecol Oncol* 133(3): 599–606.
- Kumar B et al. (2015) Suberoylanilide hydroxamic acid (SAHA) reverses chemoresistance in head and neck cancer cells by targeting cancer stem cells via the downregulation of nanog. *Genes Cancer* 6(3–4): 169–81.
- Marks PA & Breslow R. (2007) Dimethyl sulfoxide to vorinostat: development of this histone deacetylase inhibitor as an anticancer drug. *Nat Biotechnol* 25(1): 84–90.
- Richon VM et al. (1998) A class of hybrid polar inducers of transformed cell differentiation inhibits histone deacetylases. *Proc Natl Acad Sci USA* 95(6): 3003–7.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

**This product is hazardous. Please refer to the Safety Data Sheet (SDS).**

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.