

Cytokines

Human Recombinant BMP-4



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Bone morphogenetic protein 4

Catalog #	78211	20 µg
	78211.1	100 µg
	78211.2	500 µg
	78211.3	1000 µg

Product Description

Bone morphogenetic protein 4 (BMP-4) is a member of the highly conserved transforming growth factor β (TGF- β) superfamily. Mature BMP-4 is a disulfide-linked homodimeric protein consisting of two 116-amino-acid residue subunits, and is generated by the proteolytic removal of the signal peptide and propeptides (Xiao et al.). BMP-4 binds to type I and type II receptors on cells. This binding results in the phosphorylation of receptor 1, which in turn results in the phosphorylation of Smad proteins, which then go on to act as transcription factors (Zhang et al.). BMPs have been shown to be key regulators of embryogenesis and are known to play a role in the growth and differentiation of various cell types, including embryonic stem (ES) cells, induced pluripotent stem (iPS) cells, mesenchymal cells, epithelial cells, hematopoietic cells, and neuronal cells (Chadwick et al.; Graham et al.; Jones et al.; Lengerke et al.; Zhang et al.).

Product Information

Alternative Names:	BMP-2B, BMP2B1, MCOPS6, OFC11, ZYME
Accession Number:	P12644
Amino Acid Sequence:	SPKHHSQRAR KKNKNCRRHS LYVDFSDVGW NDWIVAPPGY QAFYCHGDCP FPLADHLNST NHAIVQTLVN SVNSSIPKAC CVPTLSAIS MLYLDEYDKV VLKQYQEMVV ECGGCR
Predicted Molecular Mass:	13.1 kDa monomer; 26.2 kDa dimer
Species:	Human
Cross Reactivity:	Mouse, Rat
Formulation:	Lyophilized from a sterile-filtered solution containing 0.1% trifluoroacetic acid.
Source:	CHO

Specifications



Activity:	The specific activity is $\geq 6.7 \times 10^4$ units/mg ($EC_{50} \leq 15$ ng/mL) as determined by alkaline phosphatase activity induced in ATDC-5 cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic Limulus amoebocyte lysate (LAL) analysis and is ≤ 1 EU/ μ g protein.

Preparation and Storage

Storage:	Store at -20°C to -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex.

OPTIONAL: After reconstitution, if product will not be used immediately, dilute with concentrated bovine serum albumin (BSA) to a final BSA concentration of 0.1%. The effect of storage of stock solution on product performance should be tested for each application. As a general guide, do not store at 2 - 8°C for more than 1 month or at -20°C to -80°C for more than 3 months. Avoid repeated freeze-thaw cycles.

Data

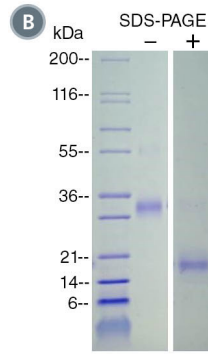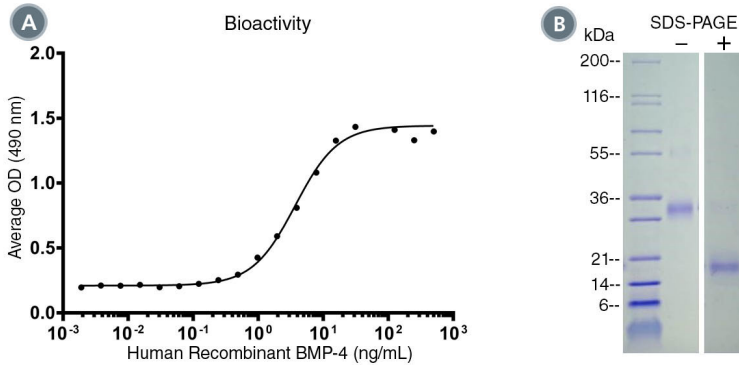

(A) The biological activity of Human Recombinant BMP-4 was tested by its ability to induce alkaline phosphatase production in ATDC-5 cells. Alkaline phosphatase production was measured using a fluorometric assay method. The EC₅₀ is defined as the effective concentration of the growth factor at which alkaline phosphatase activity is at 50% of maximum. The EC₅₀ in the above example is 3.99 ng/mL.

(B) Human Recombinant BMP-4 was resolved with SDS-PAGE under reducing (+) conditions and non-reducing (-) conditions and visualized by Coomassie Blue staining. The predicted molecular mass is 13.1 kDa (monomer) or 26.2 kDa (dimer), but these migrate to an apparent molecular mass of 20 kDa or 34 kDa, respectively, due to variable glycosylation.

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, visit www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Chadwick K et al. (2003) Cytokines and BMP-4 promote hematopoietic differentiation of human embryonic stem cells. *Blood* 102(3): 906–15.
- Graham A et al. (1994) The signalling molecule BMP4 mediates apoptosis in the rhombencephalic neural crest. *Nature* 372(6507): 684–6.
- Jones C et al. (1991) Involvement of bone morphogenetic protein-4 (BMP-4) and Vgr-1 in morphogenesis and neurogenesis in the mouse. *Development* 111(2): 531.
- Lengerke C et al. (2009) Hematopoietic development from human induced pluripotent stem cells. *Ann N Y Acad Sci* 1176(1): 219–27.
- Xiao YT et al. (2007) Bone morphogenetic protein. *Biochem Biophys Res Commun* 362(3): 550–3.
- Zhang P et al. (2008) Short-term BMP-4 treatment initiates mesoderm induction in human embryonic stem cells. *Blood* 111(4): 1933–41.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2019 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.